

MCA4A

4 Input Multichannel Analyzer
2 Input Multichannel Scaler

User Manual

© Copyright FAST ComTec GmbH
Grünwalder Weg 28a, D-82041 Oberhaching
Germany



Version 2.21, May 24, 2024

Software Warranty

FAST ComTec warrants proper operation of this software only when used with software and hardware supplied by FAST ComTec. FAST ComTec assumes no responsibility for modifications made to this software by third parties, or for the use or reliability of this software if used with hardware or software not supplied by FAST ComTec. FAST ComTec makes no other warranty, expressed or implied, as to the merchantability or fitness for an intended purpose of this software.

Software License

You have purchased the license to use this software, not the software itself. Since title to this software remains with FAST ComTec, you may not sell or transfer this software. This license allows you to use this software on only one compatible computer at a time. You must get FAST ComTec's written permission for any exception to this license.

Backup Copy

This software is protected by German Copyright Law and by International Copyright Treaties. You have FAST ComTec's express permission to make one archival copy of this software for backup protection. You may not otherwise copy this software or any part of it for any other purpose.

**Copyright © 2015-2024 FAST ComTec GmbH
D 82041 Oberhaching, Germany
All rights reserved**

This manual contains proprietary information; no part of it may be reproduced by any means without prior written permission of FAST ComTec, Grünwalder Weg 28a, D 82041 Oberhaching, Germany. Tel: +49 89 6651800, FAX: +49 89 66518040, Email: info@fastcomtec.com.

The information in this manual describes the hardware and the software as accurately as possible, but is subject to change without notice.

Table of Contents

1 Introduction.....	8
2 Installation procedure.....	9
2.1 Hard- and Software Requirements.....	9
2.2 Driver Installation.....	10
2.3 Software Installation.....	12
2.4 Hardware Installation.....	12
2.5 Getting Started.....	13
2.5.1 Initial offset calibration.....	13
2.5.2 Basic Scope measurement.....	15
2.5.3 Simple Pulse Height Analysis.....	18
2.5.4 Measuring a typical Pre-Amplifier Signal.....	22
2.5.5 Multiscaling measurement.....	26
2.5.6 Basic Usage of the RealTimeClock Option.....	28
2.6 Installing more than one MCA4A Module.....	29
3 Hardware Description.....	31
3.1 Overview.....	31
3.2 Front Panel Connectors.....	31
3.2.1 Analog ADC Inputs.....	31
3.2.2 Changing the Input range.....	32
3.2.3 GATE / MCS Inputs.....	33
3.2.4 FEATURE I/O Connector.....	34
3.3 MCA section, Analog-to-Digital Converter.....	35
3.3.1 Scope mode.....	35
3.3.2 Pulse height analysis.....	36
3.3.3 Sampled voltage analysis.....	38
3.3.4 Single channel analyzer (SCA).....	38
3.3.5 Real time, live time and system time counter.....	38
3.3.6 List mode data acquisition.....	38
3.3.7 Event tagging, spectrum multiscaling.....	38
3.4 Multichannel scaler.....	39
3.4.1 Introduction.....	39
3.4.2 Dwell time.....	39
3.4.3 Sweep counter.....	39
3.4.4 Spectrum sequencing.....	39
3.5 Octal 125 MHz Counter / Timer.....	40
3.6 Time Base / Reference Clock.....	40
4 Software Description.....	41
4.1 Server functions.....	41
4.1.1 Initialization files.....	42
4.1.2 Action menus.....	42
4.1.3 File menu.....	42
4.1.4 Settings Overview.....	44
4.1.5 MCS Settings.....	45
4.1.6 Presets.....	46
4.1.7 ADC Settings.....	47
4.1.8 Pulse-width histogramming.....	51
4.1.9 Scaler Setting.....	52
4.1.10 System Definition.....	53
4.1.11 Calculated Spectra dialog.....	54
4.1.12 Conditions.....	58
4.2 File formats.....	60
4.3 Control Language.....	62
4.4 Controlling the MCA4 Windows Server via DDE.....	68
4.4.1 Open Conversation.....	68
4.4.2 DDE Execute.....	68
4.4.3 Close Conversation.....	70
4.5 Controlling the MCA4 Windows Server via DLL.....	72
5 MPANT Software.....	73

5.1 File Menu.....	74
5.2 Window Menu.....	77
5.3 Region Menu.....	78
5.4 Options Menu.....	83
5.5 Action Menus.....	93
6 Programming and Software Options.....	94
7 Appendix.....	95
7.1 Absolute Maximum Ratings.....	95
7.2 Recommended Operating Conditions.....	95
7.3 Specifications.....	95
7.3.1 ADC inputs.....	95
7.3.2 GATE / MCS inputs.....	95
7.3.3 SCA outputs.....	96
7.3.4 Digital I/O ports.....	96
7.3.5 Reference clock I/O.....	96
7.3.6 Octal Counter inputs.....	96
7.3.7 RTC RESET input.....	96
7.4 Performance.....	97
7.4.1 General.....	97
7.4.2 Analog-to-Digital Converter.....	97
7.4.3 Counters.....	99
7.4.4 Multichannel scaler.....	99
7.5 Physical.....	101
7.6 Accessories.....	101
Optional:	101
7.7 Trouble Shooting.....	103
7.8 Personal Notes.....	104

Table of Figures

Fig. 2.1: Device Manager.....	10
Fig. 2.2: Driver Installation.....	10
Fig. 2.3: Choosing the right folder.....	11
Fig. 2.4: Device Manager after successful driver installation.....	11
Fig. 2.5: MCA4A rear panel.....	12
Fig. 2.6: Software startup screen.....	13
Fig. 2.7: MCA4A Settings window.....	13
Fig. 2.8: ADC Settings dialog (2x).....	14
Fig. 2.9: ADC settings for the first measurement.....	15
Fig. 2.10: Scale parameters for the scope display.....	16
Fig. 2.11: Scale parameters for the spectrum display.....	16
Fig. 2.12: 1st scope measurement.....	17
Fig. 2.13: Single display options dialog.....	18
Fig. 2.14: Setup a FFT analysis.....	18
Fig. 2.15: FFT analysis display.....	18
Fig. 2.16: Setup a Gaussian Fit PHA measurement.....	19
Fig. 2.17: Scope picture of a 1MHz sine wave.....	20
Fig. 2.18: Switch to firmware based PHA.....	20
Fig. 2.19: Setup the internal Test Pulser.....	22
Fig. 2.20: Setup the internal Test Pulser.....	23
Fig. 2.21: Scope Picture of the original Test Pulser Signals.....	23
Fig. 2.22: Firmware based shaping.....	24
Fig. 2.23: Scope Picture of the inverted and positive going Test Pulser Signals.....	24
Fig. 2.24: Enable the Trapezoidal Filter.....	25
Fig. 2.25: Effects of different Pole Zero Settings.....	25
Fig. 2.26: Scope Picture of the trapezoidal shaped Test Pulser Signals.....	25
Fig. 2.27: Modified shaper parameters.....	26
Fig. 2.28: MCA4 Settings: select MCS Mode.....	26
Fig. 2.29: MCS Settings for a simple measurement.....	27
Fig. 2.30: Simple MCS spectrum.....	27
Fig. 2.31: Basic RTC Experiment.....	28
Fig. 2.32: Triggering the RTC with the TRG signal.....	28
Fig. 2.33: Basic Time Spectra.....	29
Fig. 2.34: RTC Spectrum Definition.....	29
Fig. 2.35: Editing MCA4.INI.....	30
Fig. 3.1: MCA4 front panel.....	31
Fig. 3.2: Left: jumper setting for different input ranges, right: jumper positions.....	32
Fig. 3.3: Simplified analog input and ADC circuitry.....	32
Fig. 3.4: Circuit diagram of each GATE / MCS input.....	33
Fig. 3.5: Circuit diagram of the FEATURE I/O port.....	34
Fig. 3.6: Digital I/O port details.....	35
Fig. 3.7: GO-LINE details.....	35
Fig. 3.8: Pulswidth check - discard "false" pulses.....	36
Fig. 3.9: Gaussian shaped pulse fitting.....	37
Fig. 3.10: Flattop pulse fitting modes.....	37
Fig. 3.11: Pile-up detection.....	38
Fig. 3.12: Octal 125 MHz Counter Connector.....	40
Fig. 4.1: MCA4 Server Window.....	41
Fig. 4.2: MCA4A.ini File.....	42
Fig. 4.3: Data Operations dialog for MPA data (left) and selected spectra (right).....	43
Fig. 4.4: Replay Settings dialog.....	43
Fig. 4.5: About MCA4 dialog box and Installed items.....	44
Fig. 4.6: Settings overview, left: PHA mode, right: MCS mode, more than one MCA4 module.....	44
Fig. 4.7: MCS Settings for one (left) or two stop inputs (right).....	45
Fig. 4.8: Thresholds Dialog.....	46
Fig. 4.9: ADC Presets Dialog.....	46
Fig. 4.10: ADC Settings dialog, right: Scope mode.....	47
Fig. 4.11: Left: Setting Tag bits, right: advanced Settings.....	49

Fig. 4.12: Setting a ROI in the waveform for the pulse-height analysis of flattop pulses.....	50
Fig. 4.13: How to switch off Auto Scaling	50
Fig. 4.14: Setting parameters for flattop PHA.....	50
Fig. 4.15: Pulse Width Histogramming.....	51
Fig. 4.16 Scaler Settings dialog box.....	52
Fig. 4.17 Left: MCA4-Scaler Status utility, right: as ratemeter.....	52
Fig. 4.18: System Definition dialog box.....	53
Fig. 4.19: Dual parameter and Calculated spectra dialog box.....	54
Fig. 4.20: Multi Display Setting.....	54
Fig. 4.21: Multi Time Display Setting, left with installed RTC option and right without RTC option.....	55
Fig. 4.22: Calculated Spectrum Setting, with (left) and without (right) RTC and Coincidence options.....	56
Fig. 4.23: Sum of Counts dialog.....	56
Fig. 4.24: Sum of counts spectrum before and after calibration.....	57
Fig. 4.25: Conditions.....	58
Fig. 4.26: ROI Conditions dialog.....	58
Fig. 4.27: Combine Conditions dialog.....	58
Fig. 4.28: Conditions dialog with tag bits enabled.....	59
Fig. 4.29: Tag Conditions dialog.....	59
Fig. 4.30: Opening the DDE conversation with the MCA4 server in LabVIEW.....	68
Fig. 4.31: Executing a MCA4 command from a LabVIEW application.....	69
Fig. 4.32: Getting the total number of spectrum data with LabVIEW.....	69
Fig. 4.33: Getting the data with LabVIEW.....	70
Fig. 4.34: Closing the DDE communication in LabVIEW.....	70
Fig. 4.35: Control Panel of the demo VI for LabVIEW.....	71
Fig. 5.1: MPANT main window.....	73
Fig. 5.2: MPANT Map and Isometric display.....	74
Fig. 5.3: File New Display dialog.....	75
Fig. 5.4: Compare dialog.....	76
Fig. 5.5: Compare two spectra.....	76
Fig. 5.6: Print dialog.....	77
Fig. 5.7: Slice and rectangular ROI Editing dialog.....	80
Fig. 5.8: Polygonal ROI Editing dialog.....	80
Fig. 5.9: Single Gaussian Peak Fit.....	81
Fig. 5.10: Evaluation of a non-Gaussian peak.....	82
Fig. 5.11: Log file Options.....	82
Fig. 5.12: Colors dialog.....	83
Fig. 5.13: Color Palette dialog.....	83
Fig. 5.14: Single Display Options dialog	84
Fig. 5.15: Map Display Options dialog.....	85
Fig. 5.16: Isometric Display Options dialog.....	85
Fig. 5.17: Axis Parameters dialog.....	86
Fig. 5.18: Scale Parameters dialog.....	86
Fig. 5.19: Calibration dialog.....	87
Fig. 5.20: Comments dialog.....	88
Fig. 5.21: Settings dialog.....	88
Fig. 5.22: Data Operations dialog.....	89
Fig. 5.23: System Definition dialog.....	89
Fig. 5.24: Spectra dialog.....	90
Fig. 5.25: Slice dialog.....	90
Fig. 5.26: Replay dialog.....	91
Fig. 5.27: Tool Bar dialog.....	91
Fig. 5.28: Function keys dialog.....	92
Fig. 5.29: MPANT with four systems enabled.....	93
Fig. 7.1: ADC noise spectrum (generated with a ramp generator and Gaussian fit).....	98
Fig. 7.2: ADC integral linearity plot for 32k ADC range.....	98
Fig. 7.3: ADC differential non-linearity plot for 32k ADC range.....	98
Fig. 7.4: ADC DNL distribution.....	98
Fig. 7.5: Typical peak widths of Gaussian fitted pulses (32k range, 4 μ s shaping).....	99
Fig. 7.6: MCS noise spectrum (32k range, 50ns dwelltime).....	100
Fig. 7.7: MCS differential non-linearity plot (32k range, 50ns dwelltime).....	100
Fig. 7.8: MCS DNL distribution.....	101
Fig. 7.9: Counter Cable BNC or Lemo.....	102

1 Introduction

The MCA4A series is a versatile family of advanced USB3.0-connected, software controlled Multichannel Analyzers and Multiscalers.

The design is capable of converting incoming analog pulses (PHA mode) on up to 4 channels at virtually no conversion induced dead-time. Pulse Height Analysis can be performed using different fitting algorithms optimized for different input pulse forms and shaping times. A digital Pole Zero Filter and Trapezoidal Pulse Shaper enables to even process typical Pre-Amplifier signals.

PHA data is transferred via the USB 3.0 interface in list-mode at a max. rate of > 1 M events / sec on each input. Four ADCs take 16bit data samples continuously every 8 ns, providing spectra with 32k resolution using pulse-height analysis programmed in the FPGA firmware.

The availability of the simultaneous usage of the oscilloscope (Scope) mode together with a Fast Fourier Transformation (FFT) and software pulse height analysis allows for an easy setup and evaluation of the experiment.

The PHA data can be tagged with the 8 digital I/O bits and also with the 8ns Real Time Clock.

A set of eight 125MHz Counter/Scalers is also provided.

In Multiscaling (MCS) mode the internal memory is used to accumulate spectra of up to 16M bins. An automatic sequential mode allows to acquire up to 65536 sequential spectra – each with a pre-settable number of sweeps – with no dead-time between each cycle or synchronized with an external signal. The maximum count rate is 400 MHz, dwelltime is ≥ 30 ns for one STOP input or ≥ 60 ns for two STOP inputs.

A versatile 8 bit digital I/O port provides controlling external devices or to react on additional external signals.

The “GO-LINE” compatibility allows to synchronize begin and end of the experiment across many measurement devices.

The 32 bit or 64 bit MPANT Windows Software is able to handle 6 MCA4A providing up to 24 MCA channels. In-field upgrade options enable upgrading of MCS mode, more MCA channels, 8 ns RTC time tagging, and multiparameter coincidences for 2 or 4 ADCs including Replay of list files.

2 Installation procedure

2.1 Hard- and Software Requirements

For operating the MCA4A you need a standard PC with a high-speed USB 2.0 or USB 3.0 port and Microsoft Windows XP or higher. USB 1.0 ports do not work. The maximum throughput rate is approximately 32 Mbyte/sec on a USB 2.0 port and 120 Mbyte/sec on a USB 3.0 port. This throughput rate is somewhat dependent on the computer hardware and operating system.

We do not expect any problems with host compatibility as we do not use something else than any USB hard disk or memory stick. Of course, it is necessary that the PC has a true USB 2.0 or USB 3.0 / USB 3.1 port.

Older PCs that do not have USB 3.0 on the motherboard can use a PCI / PCIe card providing USB 3.0 / USB 3.1 Gen1 ports.

First check you have all shipped equipment available:

- Transport case
- MCA4A module
- Power supply module 100...240 V AC / 12 V DC out
- Power line cord
- USB 3.0 (A/A) cable, 3 m
- User manual
- CDR or USB stick with operating software

2.2 Driver Installation

To install the driver, please start the device manager. You can do it by right-clicking the computer shortcut on the desktop and selecting "Manage" from the drop-down-menu, then "Device Manager". Another way is to select in the Start Menu Settings..Control Panel, then System..Device Manager.

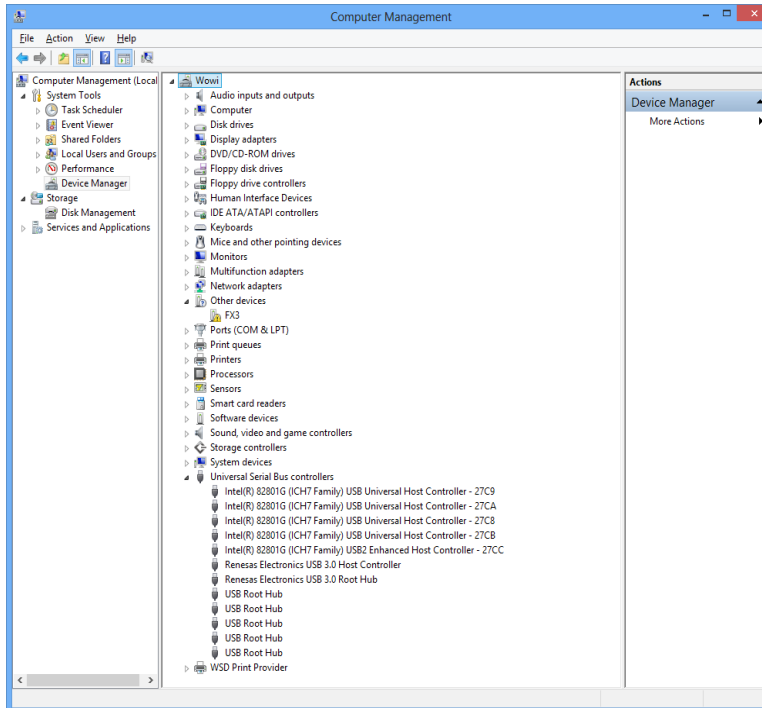


Fig. 2.1: Device Manager

The new device is displayed as "FX3" in a class "Other devices". Please right-click it and select "Update driver" from the drop-down menu.

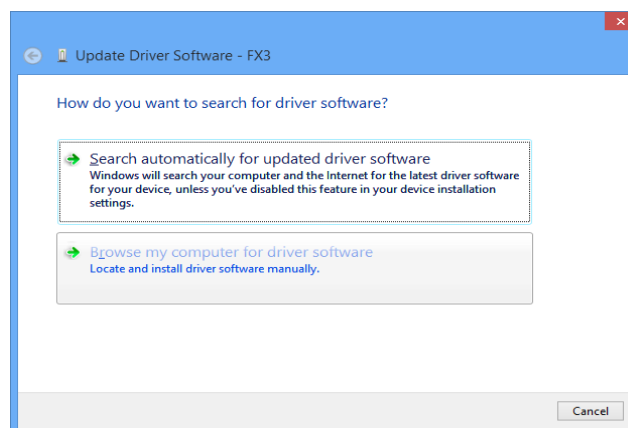


Fig. 2.2: Driver Installation

Windows offers to search automatically for updated driver software or browse the computer for the driver software. Don't try the automatic installation, it will probably not work. Select to browse manually for the driver and select for example for Windows 7 or 8 or 10 the driver\fastusb3 on the CD or Software stick. If you have an MCA4 that you want to operate with the same software, please install the driver in the driver\fastusb2 folder.

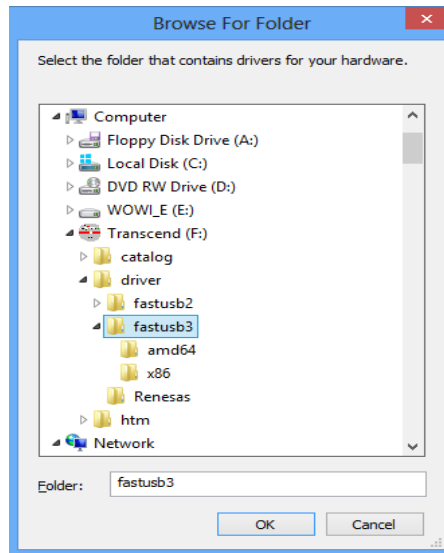


Fig. 2.3: Choosing the right folder

Please don't go deeper into the folder structure, the processor type (x86, amd64, etc.) is found automatically. After pressing "OK" the driver installation will proceed. The device manager will then show a "FAST ComTec USB3 Device" in a class "libusb (WINUSB) devices" as shown in Fig. 2.4.

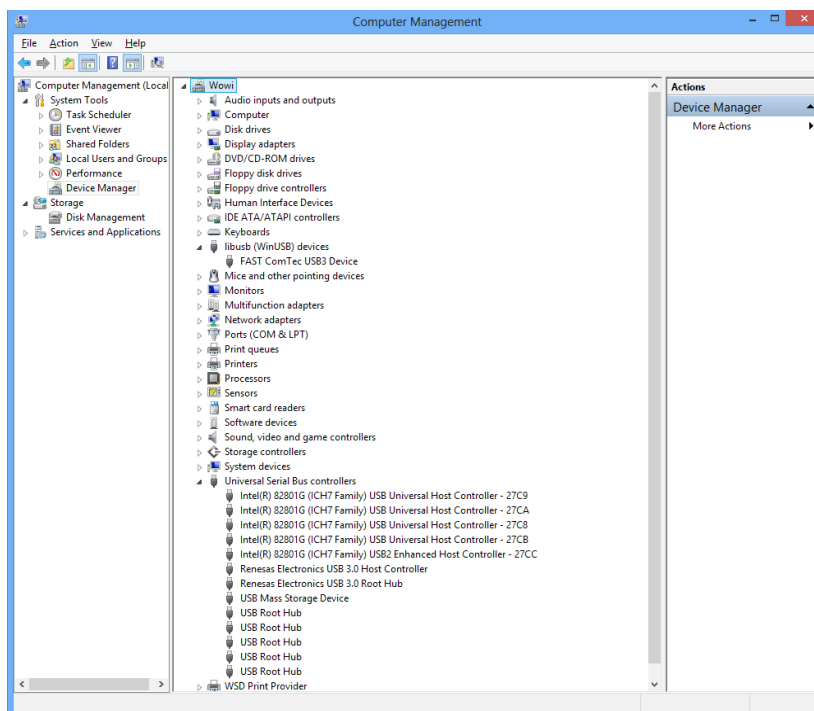


Fig. 2.4: Device Manager after successful driver installation

If there is an exclamation mark at the driver and there is an error message that the device could not start, please remove the power from the MCA4A device and reconnect it. If you don't see the device at all, it is possible that the USB3 driver must be updated. For a Renesas USB3 controller there is an update program on the software disk in the driver\Renesas folder.

Please check the Power Saving Options in the PC Settings. If the PC goes to sleep during a running acquisition and shuts the USB port down, the software will hang up. So it is recommended to disable automatic going to sleep in the energy options.

2.3 Software Installation

To install the MCA4A software on your hard disk insert the installation medium (CD or USB Stick) and start the installation program setup.exe for the 32-bit software or setupx64.exe for the native 64-bit software.

A directory called C:\MCA4A for the 32-bit software or C:\mca4ax64 for the native 64-bit software is created on the hard disk and all MCA4A and MPANT files are transferred to this directory. Drive C: is taken as default drive and \mca4ax64 as default directory. It is not mandatory that the MCA4A operating software is located in this directory. You may specify a different directory during the installation or may copy the files later to any other directory.

The Setup program will install a shortcut on the desktop. The icon directly starts MCA4A.EXE. The server program will automatically call the MPANT.EXE program when it is executed. The MCA4A Server program controls the MCA4A module but provides no graphics display capability by itself. By using the MPANT program, the user has complete control of the MCA4A along with the MPANT display capabilities.

To run the MCA4A software, simply double click on the “MCA4A Server Program” icon. To close it, close the MCA4A server in the Taskbar.

2.4 Hardware Installation



Fig. 2.5: MCA4A rear panel

Installation of the MCA4A is as easy as connecting a cable. As the USB 3.0 interface is compatible to USB 2.0 it is just a question of the maximum transfer speed which cables or host ports are used. Connect the power supply to the “POWER” connector and an USB cable to the “USB” port. At the host computer plug in the USB cable into a high-speed capable USB port. Since USB is a hot-pluggable interface the sequence of applying power and connecting the USB port is not important.

When the power is applied the “POWER ON” LED on the rear of the MCA4A should be lit and a second later “POWER GOOD” as well.

When the software and driver are already installed the computer will detect the MCA4A automatically as soon as it is connected.

2.5 Getting Started

2.5.1 Initial offset calibration

First step is to start the MCA4A software by double clicking the corresponding icon. This will automatically start the MPANT program as well. On startup the MCA4A Server is iconized and one does not have to worry about it since all hardware settings are accessible from the MPANT program which actually is the graphical user interface and which will appear now on your screen.

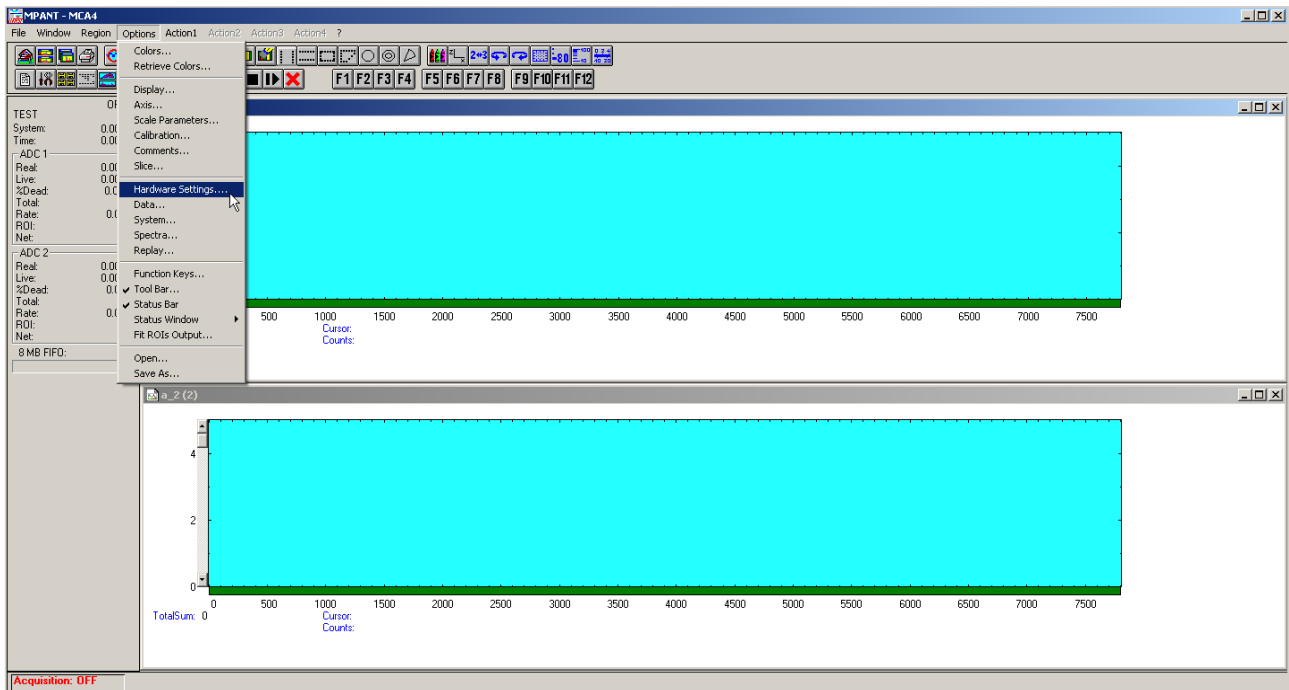



Fig. 2.6: Software startup screen

At the very first start of the software and the MCA4A a calibration of the internal hardware based offsets should be performed and the resulting data stored. Click on *Options – Hardware Settings...* or  (ref. Fig. 2.6) to find the MCA4A Settings window pop up (Fig. 2.7).

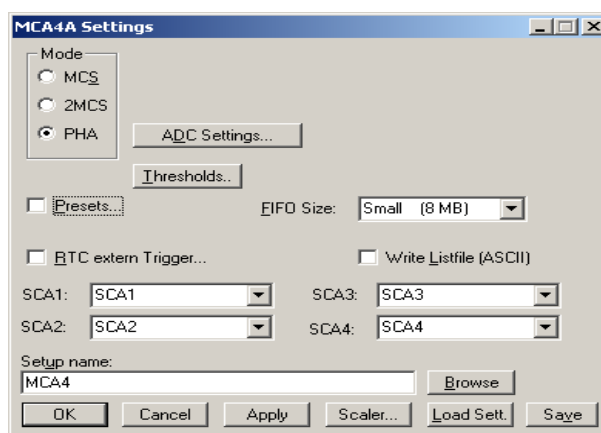


Fig. 2.7: MCA4A Settings window

Here select PHA Mode and then click on *ADC Settings...* to get the ADC Settings dialog window (ref. Fig. 2.8 - left).

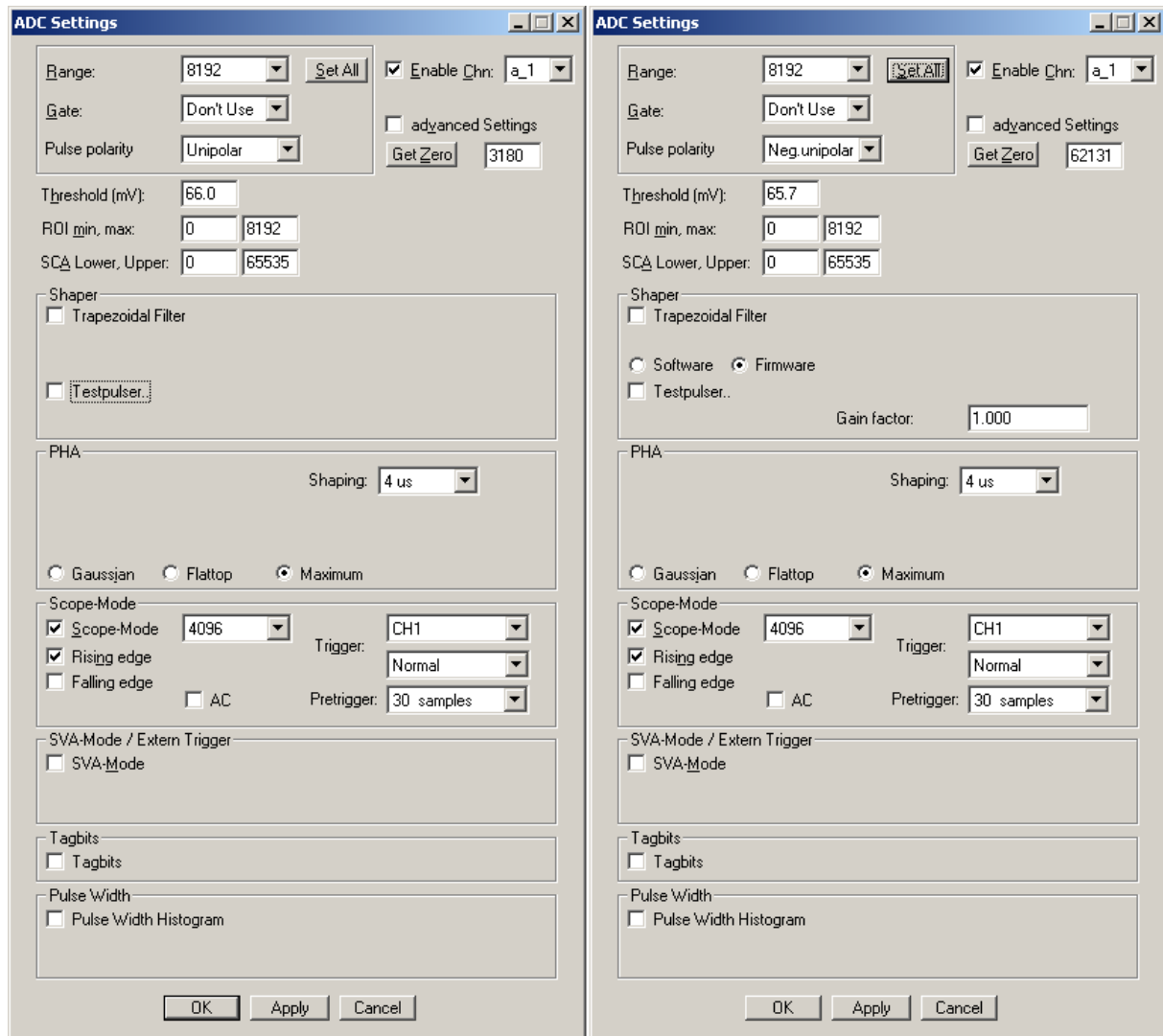


Fig. 2.8: ADC Settings dialog (2x)

Make sure you have no signals connected to the MCA4A inputs! Click *Get Zero* and then *Apply* to perform a measurement of the input offsets. The zero value is expected to be somewhere around 3,200.

Since the MCA4A provides unipolar, bipolar, negative unipolar and negative bipolar input ranges now switch to the *Neg. unipolar* range by selecting *advanced settings – Neg. unipolar Pulse polarity*. If you have more than one ADC input you may set the same polarity for all inputs simultaneously by clicking the *Set All* button. You will be asked if you want new offsets for all channels, then new zero offsets are measured automatically after Pulse polarity switching for all channels. The zero value (offset) should be in the magnitude of some 62,000 now (ref. Fig. 2.8 - right hand side). For the bipolar input range the zero value should be around 32,000.

Switch back to the positive unipolar input range, click *Apply* or *Set All* and *OK* and in the MCA4A Settings (ref. Fig. 2.7) dialog on *Save* to store the measured offset data.


NOTE:

A calibration of the internal offsets (zero values) for all input ranges (pos. unipolar, bipolar and neg. unipolar Pulse polarity) should be performed at the first startup
Please let the MCA4A warm up for some 30 minutes before calibration is performed

2.5.2 Basic Scope measurement

This might be the best thing to start with. In case you don't know the details of your input signals you have the option to watch and analyze them like with any oscilloscope. You will be able to find problems with the setup like ringing or reflections. Even defective cables can be identified. It is also easy to adjust your amplifiers' offset and pole zero settings without the need of an additional oscilloscope.

We use a standard waveform generator (e.g. 10kHz square wave) connected to the "ADC 1" input (ref. Fig. 3.1).

In the MPANT window click on *Options – Hardware Settings...* or  (ref. Fig. 2.6) to find the MCA4A Settings window pop up (Fig. 2.7. Here select PHA Mode and then click on *ADC Settings...* to get the ADC Settings and Presets dialog window.

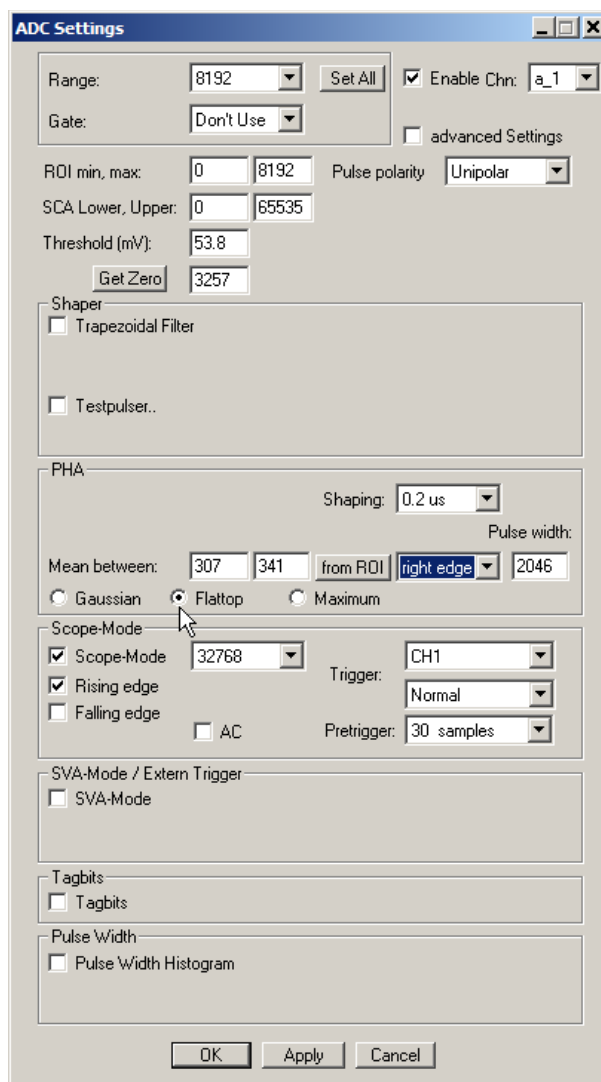



Fig. 2.9: ADC settings for the first measurement

Referring Fig. 2.9 do *Enable* ADC input channel 1 (*Chn: a_1*). Choose *PHA – Flattop* to allow pulse height analysis of the square wave signal. Set *Scope-Mode - Scope-Mode, Rising edge, Trigger: CH1* and the length of the scope window (here the maximum range of 32768). The trigger *Threshold* is set to 53.8mV.

Before we will start the measurement we first set the display windows scalings.

Click on the scope display window (A1 (3) in Fig. 2.12) to make it the active window in MPANT. Select *Options – Scale parameters...* or click  and set the *Counts Range* to 0 ... 65535 and deselct *Auto scale* and *Minimum auto scale* (ref. Fig. 2.10). Press OK. This makes sure we will see the full 16 bit ADC conversion range.

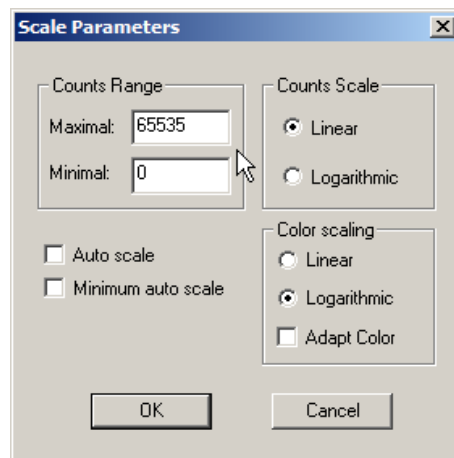


Fig. 2.10: Scale parameters for the scope display

After activation of the spectrum display (a_1 (2) in Fig. 2.12) again select *Options – Scale parameters...* and enable *Autoscale* (ref. Fig. 2.11). Press OK.

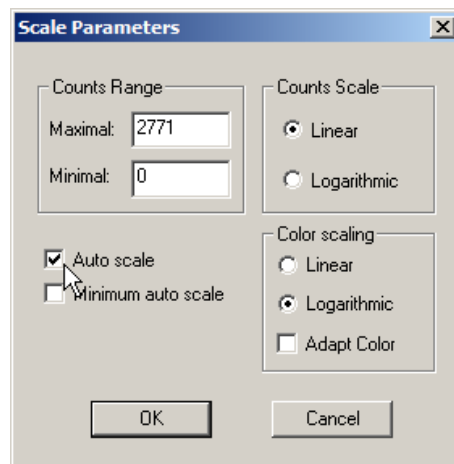




Fig. 2.11: Scale parameters for the spectrum display

Now select *Action1 – Start* or click  and start our first measurement.  halts the measurement.

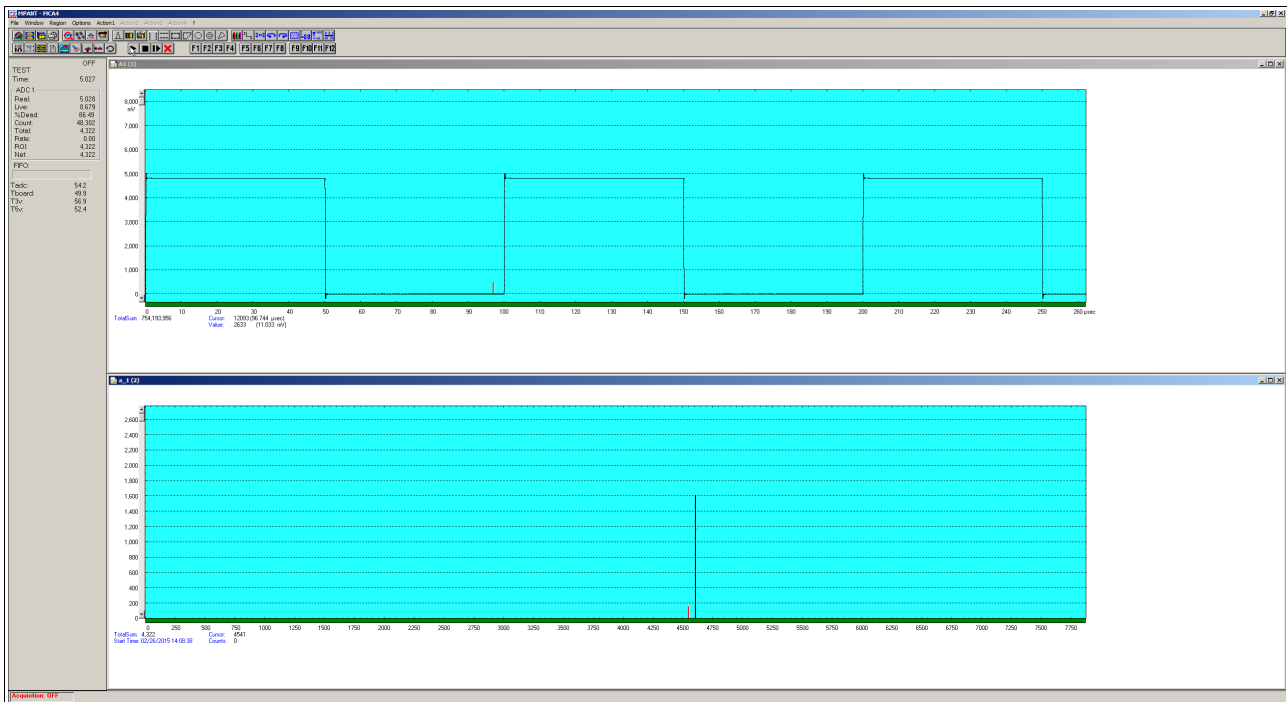


Fig. 2.12: 1st scope measurement


In the upper client window (A1) you see the scope picture of our square wave while in the lower (a_1) the spectrum of the analysed pulse height is displayed.

Holding the left mouse button in a display window does set and move the cursor which values are shown just below the x-scale. Hold the right mouse button to spread a range selection which can be used to zoom into the window.

If the trigger threshold was not chosen properly and no sampling begins then select *Trigger: auto*. Measure the sampled signal by holding the left mouse button and use the cursor values below the scope window x-scale to find a good value for the trigger threshold and revert to normal CH1 trigger.

Looking at Fig. 2.12 one can easily see the square wave had an amplitude of some 4.8V with little over- and undershoot at the edges. The amplitude was stable enough which can be seen from the narrow peak (needle) in the lower pulse height analysis spectrum.

An additional feature is the provision of a Fast Fourier Transformation. This shows the input signal in the frequency domain. It helps to find noise pickup or parasitic signal sources that would downgrade experimental results.

Activate the scope window (A1) and click on *Options – Display...* or  to get the Single Display Options dialog (Fig. 2.13).

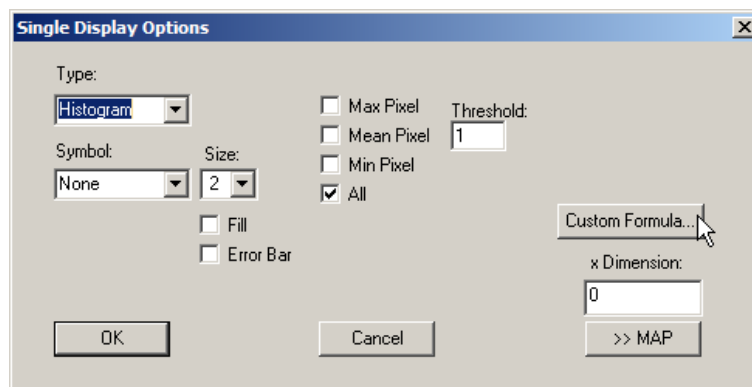


Fig. 2.13: Single display options dialog

Here select *Custom Formula...* to get the *Custom-transformed spectra* setup. Choose *Formula #: 6* (ref. Fig. 2.14, left side) which is the FFT and click *Create*.

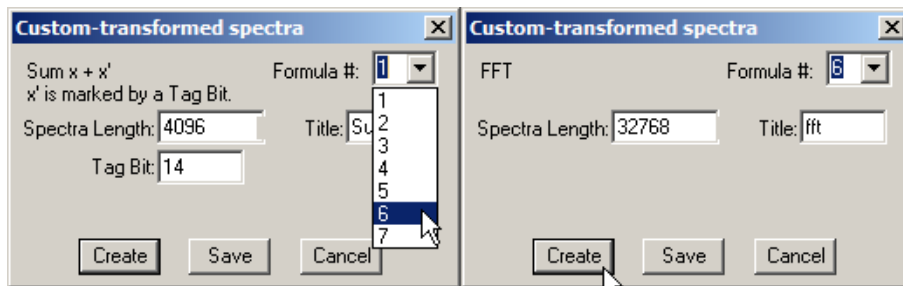



Fig. 2.14: Setup a FFT analysis

This creates a third display window containing the Fast Fourier Transformation data (ref. Fig. 2.15). Most probably you might want to switch the scaling from linear to logarithmic or vice versa (use the  button).

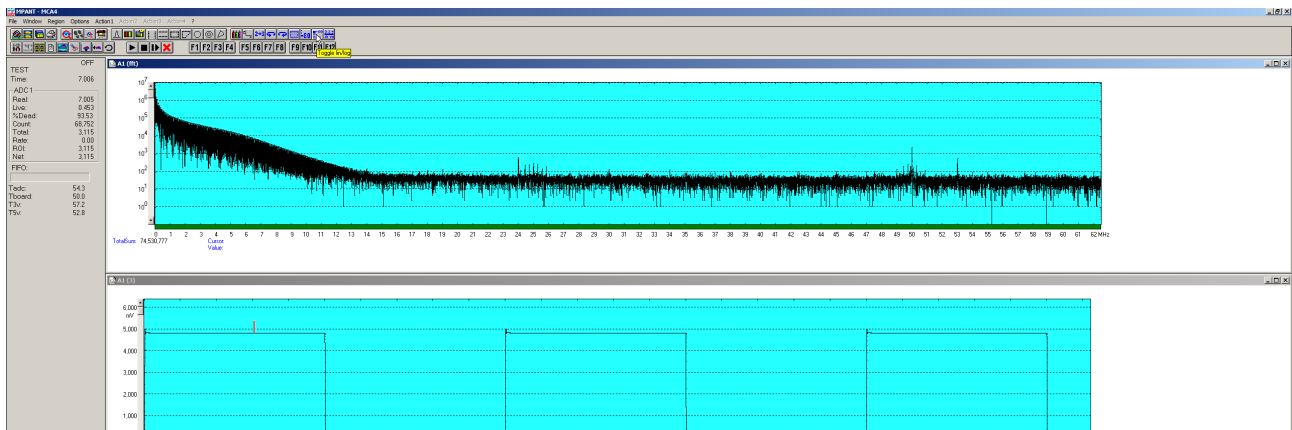



Fig. 2.15: FFT analysis display

2.5.3 Simple Pulse Height Analysis

We already did analyse the pulse height of a square wave signal in the previous chapter. Since we had the scope mode active we did the PHA in software. Now we will setup an experiment using the internal firmware based pulse height analysis. This provides far higher conversion rates (> 1MHz) without loss.

We use a general sine wave generator. The upper half of a sine wave actually does not look so different from a Gaussian shaped amplifier signal. We set the frequency to 1MHz, the offset such that the signal is mostly in the positive voltage range. E.g. set the amplitude to 6V_{PP} and the Offset to 3V. Connect the signal to the ADC 1 input (ref. Fig. 3.1).

First we want to verify what the input signal actually looks like. In the MPANT window click on *Options – Hardware Settings...* or  and select *Scope-Mode* with a range of 4096 and *Trigger = CH1*. As PHA fit algorithm choose *Gaussian*. We also need to select the shaping time. After the first measurement you can get help selecting the proper shaping time with the *guess* button (ref. Fig. 2.16). For a 1MHz sine wave 0.5 μs is a good choice. Click *OK*.

Installation procedure

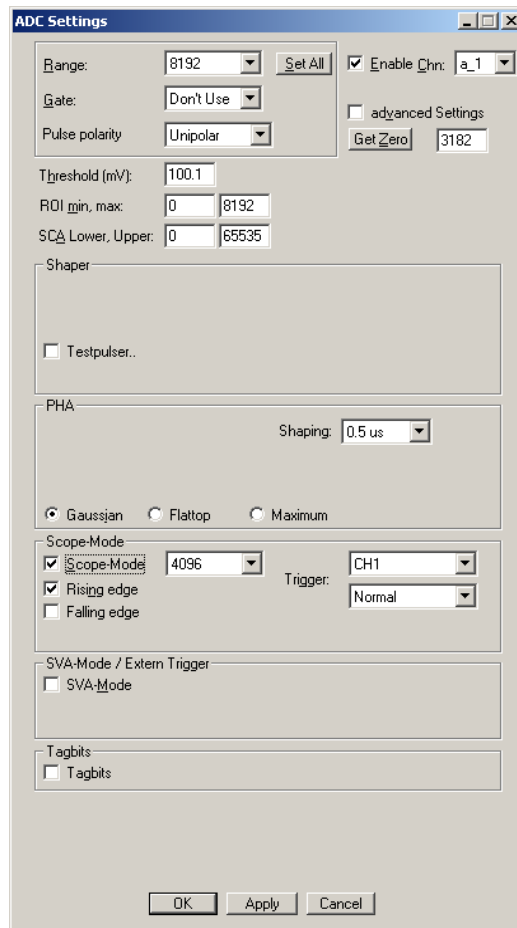




Fig. 2.16: Setup a Gaussian Fit PHA measurement

Select *Action1 – Start* or click  and start the scope measurement -  halts the measurement. The resulting picture is in Fig. 2.17. You clearly see the sine wave with 1MHz and 6V_{PP} amplitude. The PHA spectrum below shows a rising peak near channel 5559.

As can be seen near the mouse pointer the converted signal rate is only some 2kHz. To be able to process higher pulse rates the firmware based PHA fitting must be enabled and scope mode switched off.

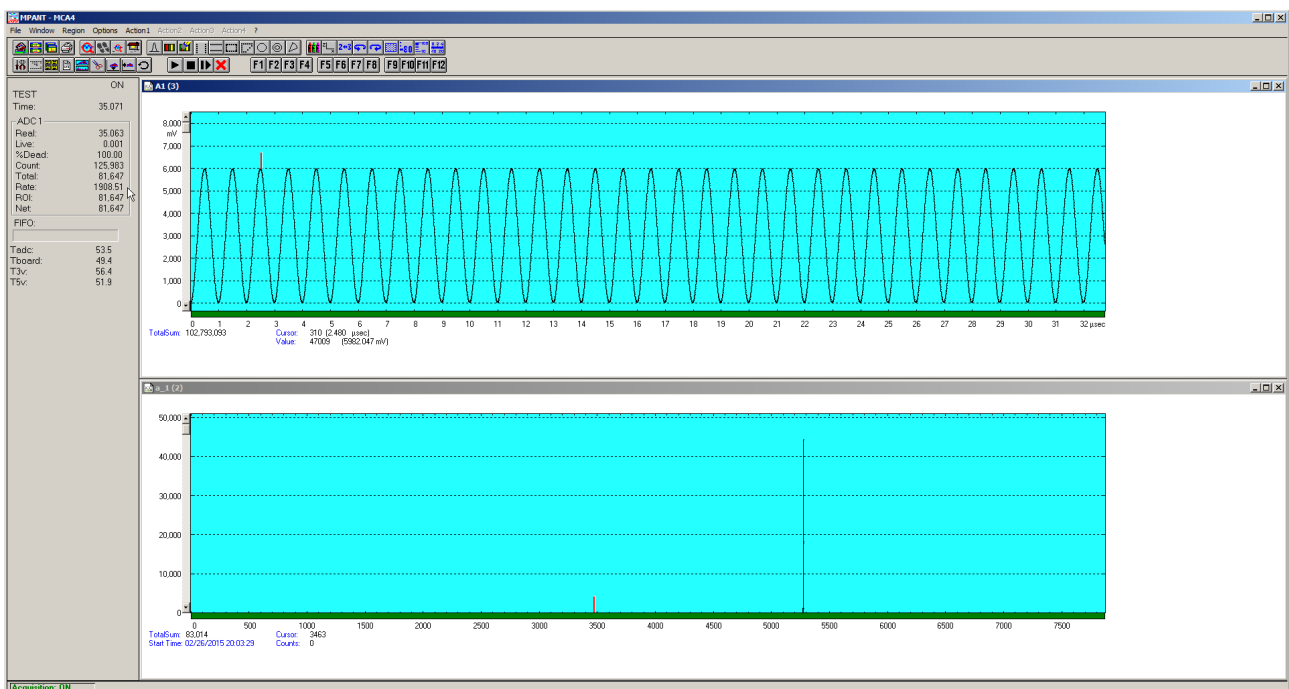


Fig. 2.17: Scope picture of a 1MHz sine wave

Again go to the *ADC Settings and Presets* dialog and deselect *Scope-Mode* (ref. Fig. 2.18).

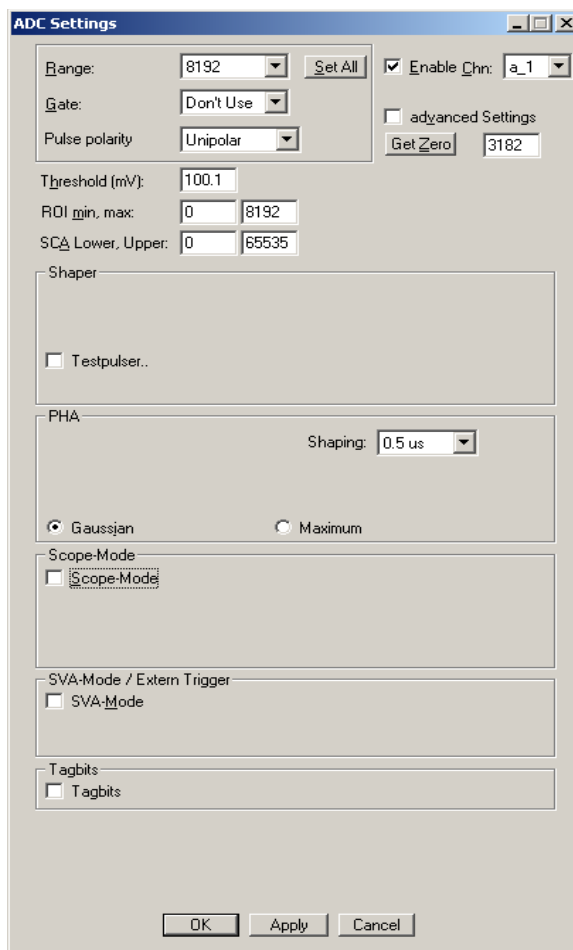


Fig. 2.18: Switch to firmware based PHA

Now only one display (the PHA spectrum) remains in the MPANT window. After starting the measurement again a peak begins to rise approximately at channel 5559 and you can see a data rate of almost 1MHz. This rate depends a bit on the computer and the USB hubs used.

2.5.4 Measuring a typical Pre-Amplifier Signal

We will now use the MCA4A's internal test pulser to produce pre-amplifier like input signals and setup the pole zero filter and the trapezoidal shaper.

Since most pre-amplifier signals are typically negative we switch to the negative unipolar pulse polarity and use the scope mode to see how the input signals look like. In the ADC Settings dialog select *Neg. unipolar*, set *PHA* to *Maximum* and *Scope-Mode* with *Falling edge* trigger. Then setup the Shaper with *Software* (ref. Fig. 2.19).

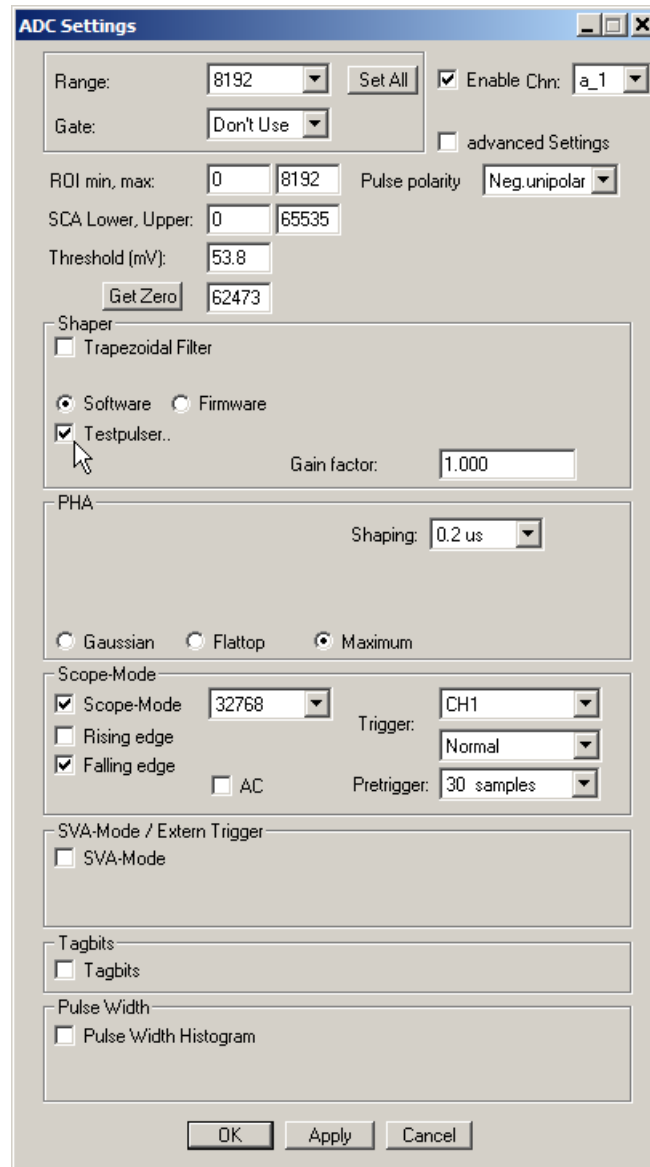


Fig. 2.19: Setup the internal Test Pulser

When clicking into the *Testpulser..* checkbox a Test pulser setup window opens. Setup the pulser with a fast rise time of "2", a long fall time of "200", an amplitude setting of "50" and a period of "5" x 32.8µs as in Fig. 2.20.

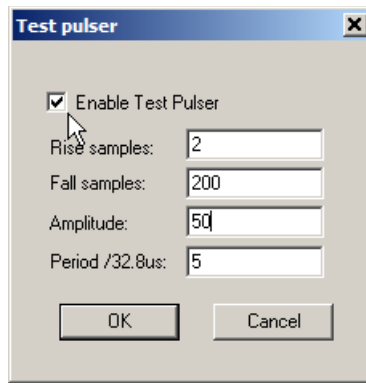


Fig. 2.20: Setup the internal Test Pulser

Starting the measurement we see the negative going input pulses as in Fig. 2.21.

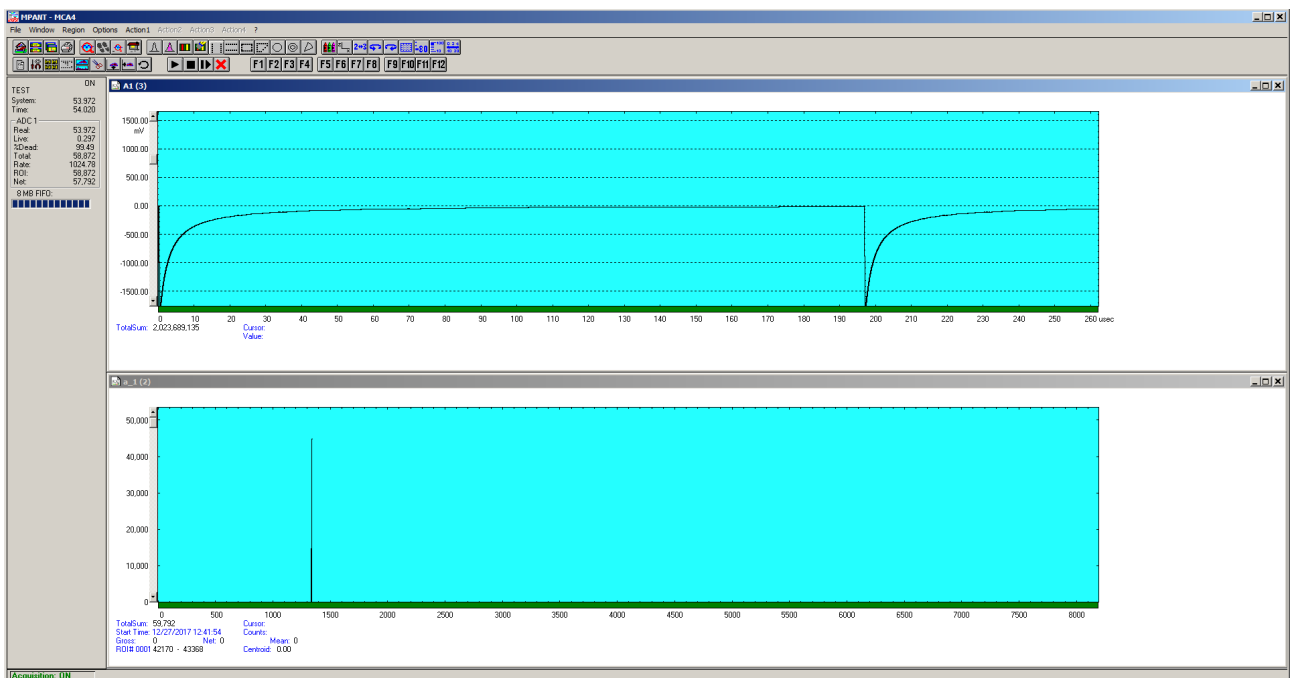


Fig. 2.21: Scope Picture of the original Test Pulser Signals

Now switch to *Firmware* based shaping (ref. Fig. 2.22). This inverts the input signals which is necessary since the firmware based shaping and pulse height analysis procedures are applicable for positive going pulses only.

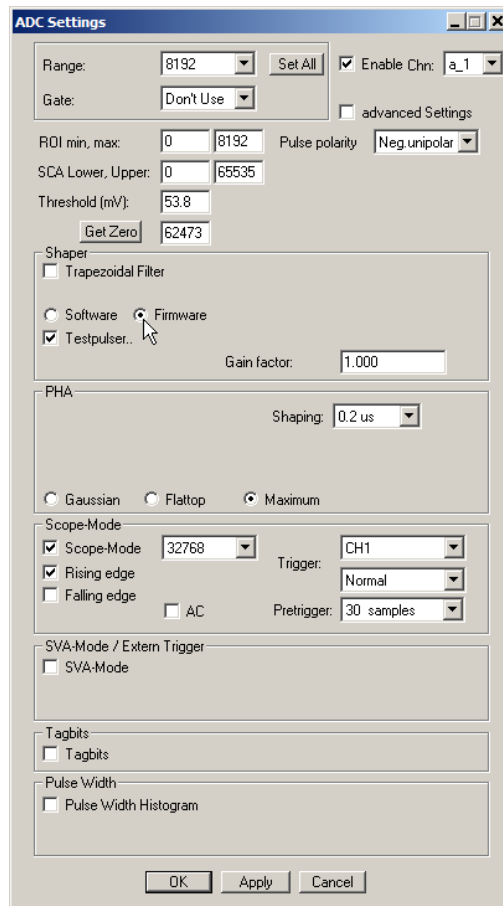


Fig. 2.22: Firmware based shaping

Now we can see positive going, inverted signals that are applicable to the firmware shaping and analysis procedures.

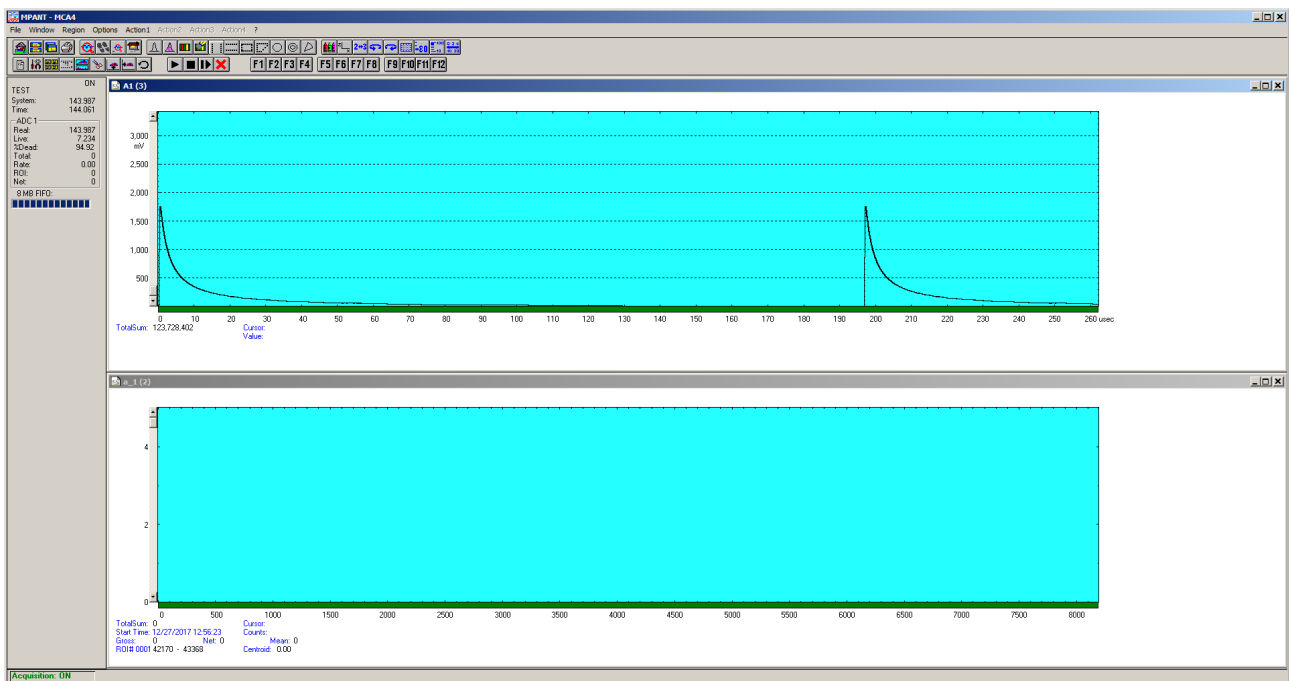


Fig. 2.23: Scope Picture of the inverted and positive going Test Pulsar Signals

Installation procedure

In the ADC Settings dialog enable the Trapezoidal Filter (ref. Fig. 2.24).

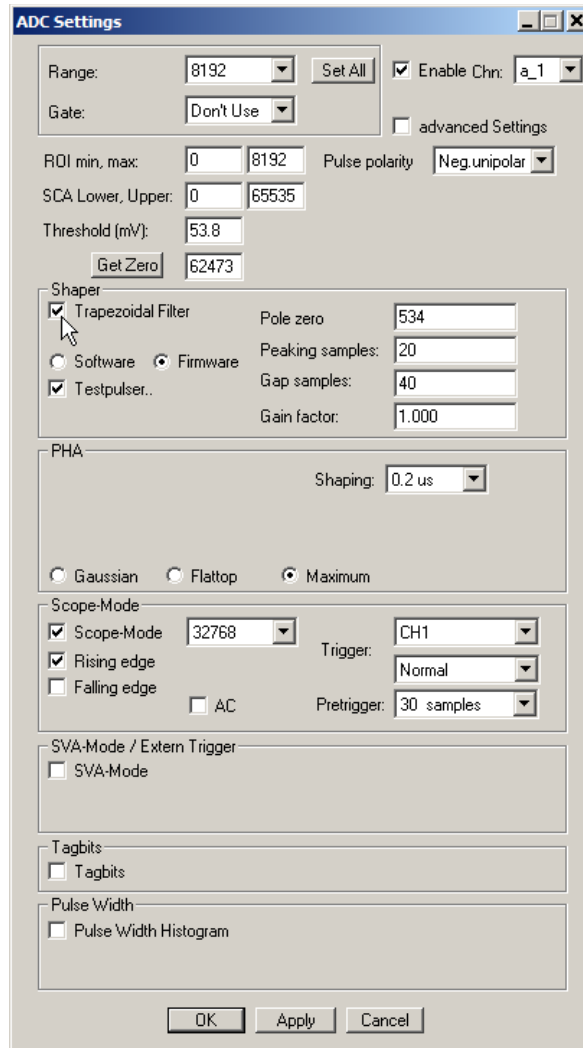


Fig. 2.24: Enable the Trapezoidal Filter

The resulting scope picture shows the trapezoidal shaped pre-amplifier pulses as in Fig. 2.26. With the now open Pole Zero Adjust window you may change the settings of the pole zero filter and optimize the resulting shaped pulse form.

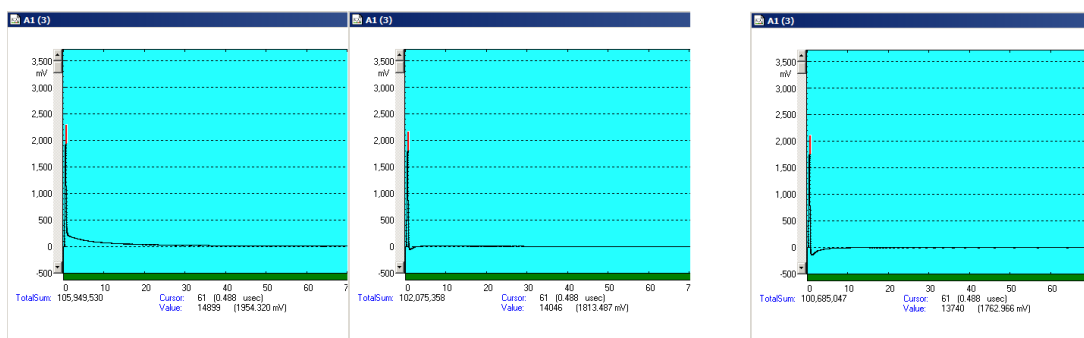


Fig. 2.25: Effects of different Pole Zero Settings

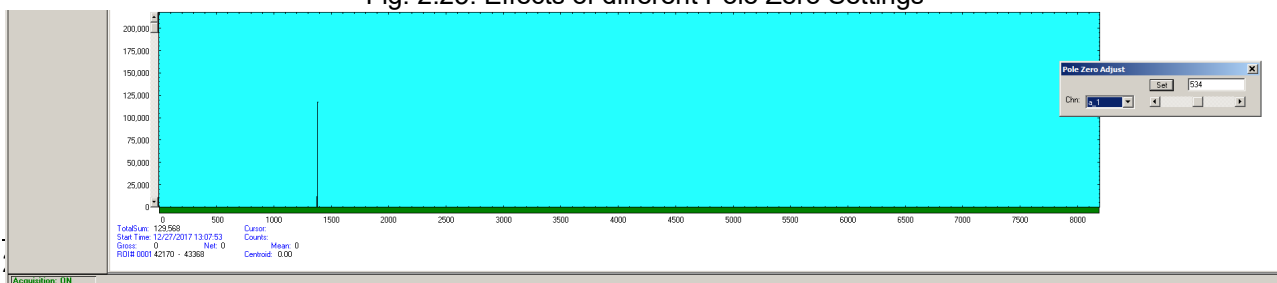


Fig. 2.26: Scope Picture of the trapezoidal shaped Test Pulsar Signals

Also try to optimize the Peaking and Gap samples. In Fig. 2.27 you see the resulting waveform with Peaking samples = 40 and Gap samples = 80.

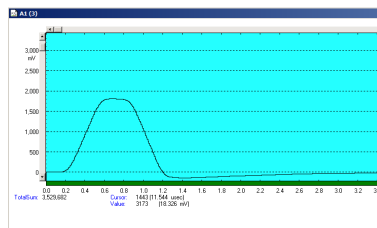


Fig. 2.27: Modified shaper parameters

2.5.5 Multiscaling measurement

To setup a simple multiscaling measurement we need a square wave generator that provides 2 synchronous signals. This can be achieved with e.g. a function generator with additional SYNC output. Or you use 1 signal that is connected via a power splitter to the MCA4.

In our example we have a function generator that is set to square wave with 0...3.3V amplitude like a low voltage TTL signal and a frequency of 100kHz. We connect the signal output to STOP 1 and the SYNC output to START.

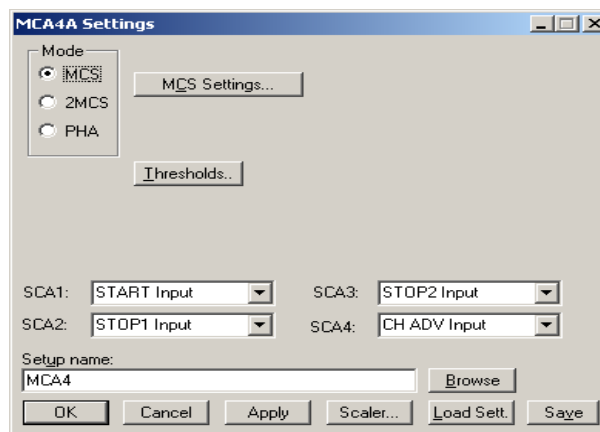



Fig. 2.28: MCA4 Settings: select MCS Mode

Now click on *Options – Hardware Settings...* or  to get the MCA4 Settings dialog pop up. Select *Mode – MCS* and click on *MCS Settings...* (ref. Fig. 2.28) to get to the MCS Settings window (Fig. 2.29).

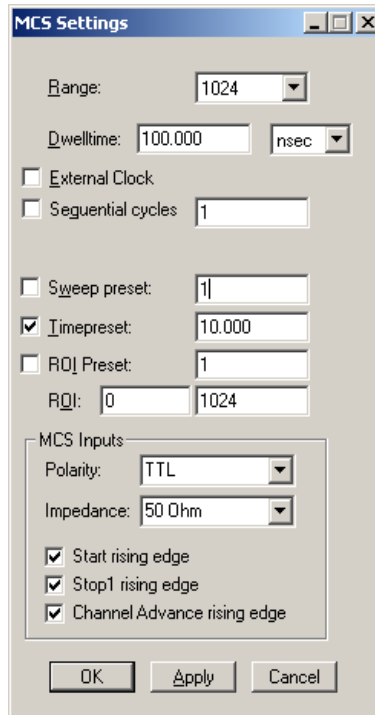


Fig. 2.29: MCS Settings for a simple measurement

Set Range = 1024, Dwelltime = 100ns, MCS Inputs Polarity to TTL, Impedance = 50 Ohm and click OK.

After starting the measurement a garden fence like spectrum begins to grow (ref. Fig. 2.30). The distance of the peaks is 100 channels corresponding to $100 \times 100\text{ns} = 10\mu\text{s}$ or 100kHz.

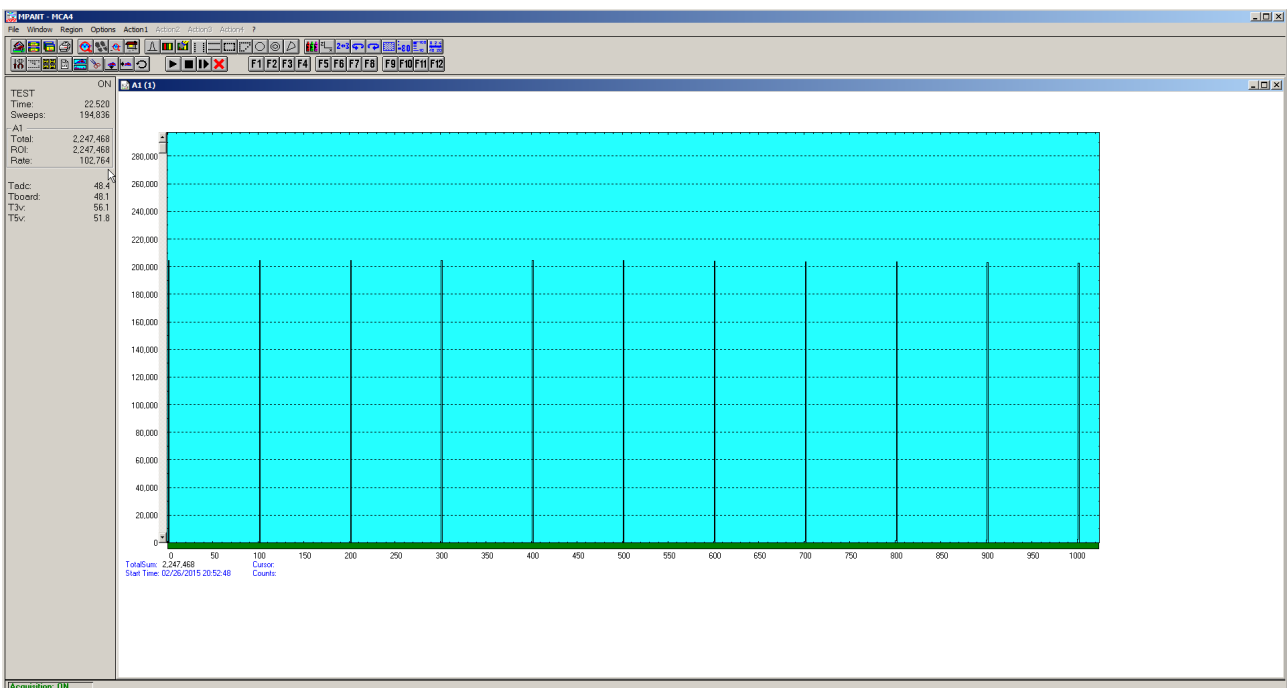


Fig. 2.30: Simple MCS spectrum

2.5.6 Basic Usage of the RealTimeClock Option

To familiarize with the usage of the RealTimeClock (RTC) option a simple experiment is setup. The intention is to measure the arrival time of single ADC events relatively to a start (trigger) signal like it might be done in Time-of-Flight or similar experiments. To do so a variable delay is used to shift analog output pulses relatively to the TTL trigger pulse that on the other hand resets (reloads) the RTC via the RTC RESET (TRG) input. Thus, in terms of TOF, the trigger signal acts as start and the ADC signal as stop input. The delay time is then measured with a resolution of 8 ns (or 1 msec if the 8 ns RTC option is not enabled) and a time spectrum of the very ADC is accumulated.

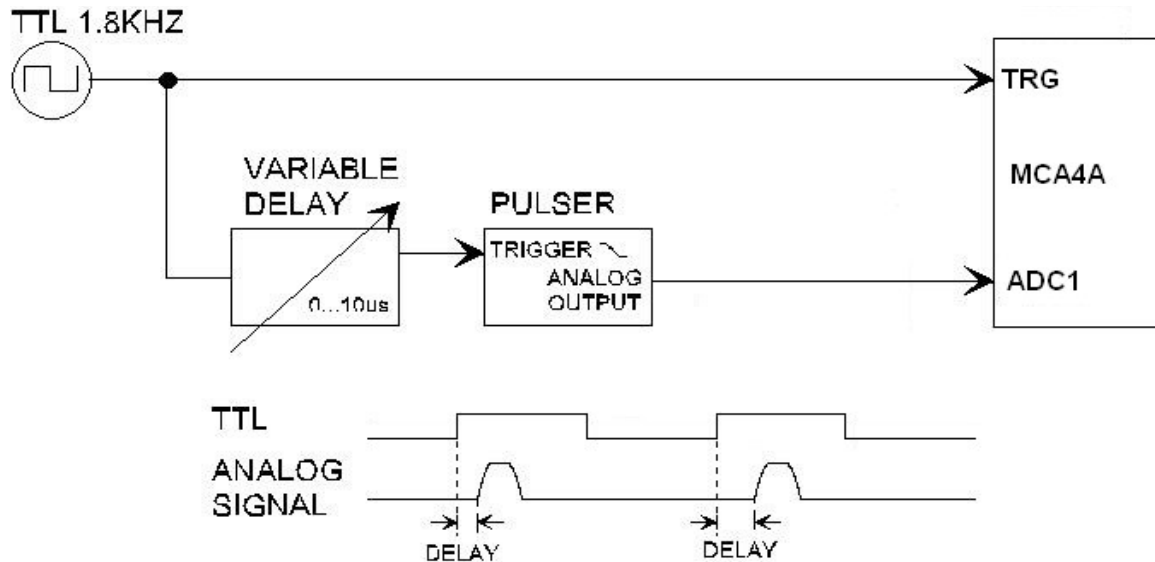


Fig. 2.31: Basic RTC Experiment

The RTC is set to reload (restart) with a RTC RESET signal at the TRG connector (ref. Fig. 2.32). The timer value capture command is derived from ADC 1 signal crossing the threshold. On-board Histogramming must be disabled so the event data contain a time stamp. If you have the RTC option for that channel, the time stamp is in units of 8 ns, otherwise in units of 1 msec.

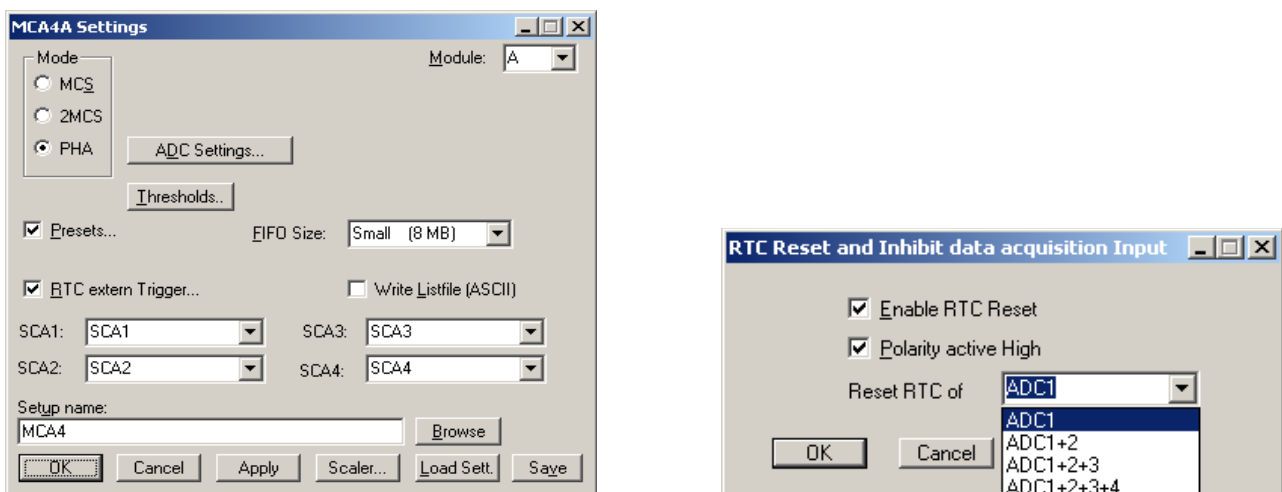


Fig. 2.32: Triggering the RTC with the TRG signal

The TRG is a 50 Ohm terminated TTL input.

To visualize the timing relationship a spectrum is defined that shows the arrival time on the x-axis, the ADC pulse height on the y-axis and the count rate in z-direction (ref. Fig. 2.34). Also a one dimensional spectrum is defined (set ADC range to 1) to show just a time spectrum (projection onto the x-axis).

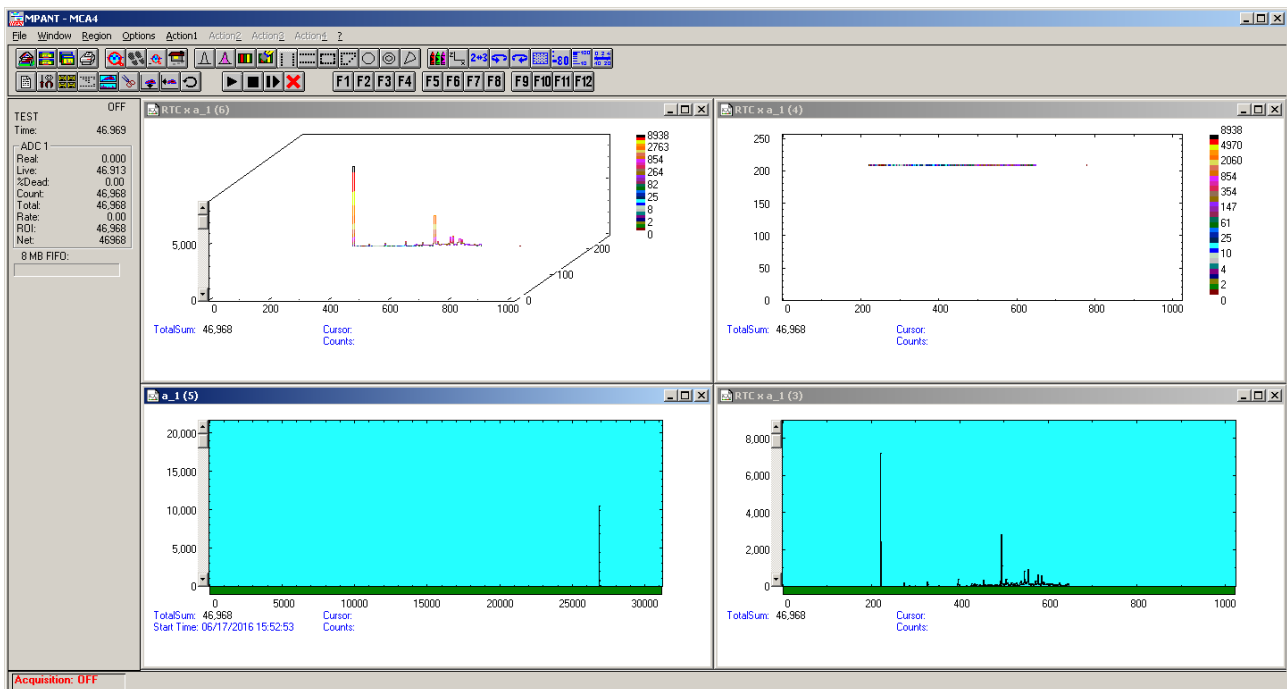


Fig. 2.33: Basic Time Spectra

In Fig. 2.33 the resulting spectra can be seen. Watching the map display (window (4)) amplitude variations show up as vertical lines whereas delay time variations result in horizontal lines. Window (6) is the corresponding three-dimensional view, window (5) the pulse height analysis and (3) the time spectrum. The acquisition will not start before a trigger signal is sent to the TRG connector at the rear panel. The Real timer is reset by each new signal there, so the real time cannot longer be used for normalizing. The system time, computer run time and live time may be used instead of it. The trigger signal should be a short pulse. For high polarity, the RTC is stopped and hold at zero as long as the signal level is high, and the acquisition is inhibited. After the pulse the acquisition is running for the time the signal level remains low.

2.6 Installing more than one MCA4A Module

The software is able to operate up to 6 MCA4A modules with a single computer. Just connect them at an USB port, switch on the power, and install the driver when you are prompted. When starting the software, you will be asked to edit the MCA4.INI file as more than one MCA4A modules are found, and the notepad editor will be automatically started with MCA4.INI loaded. Edit the line "devices=1" accordingly and save the file.

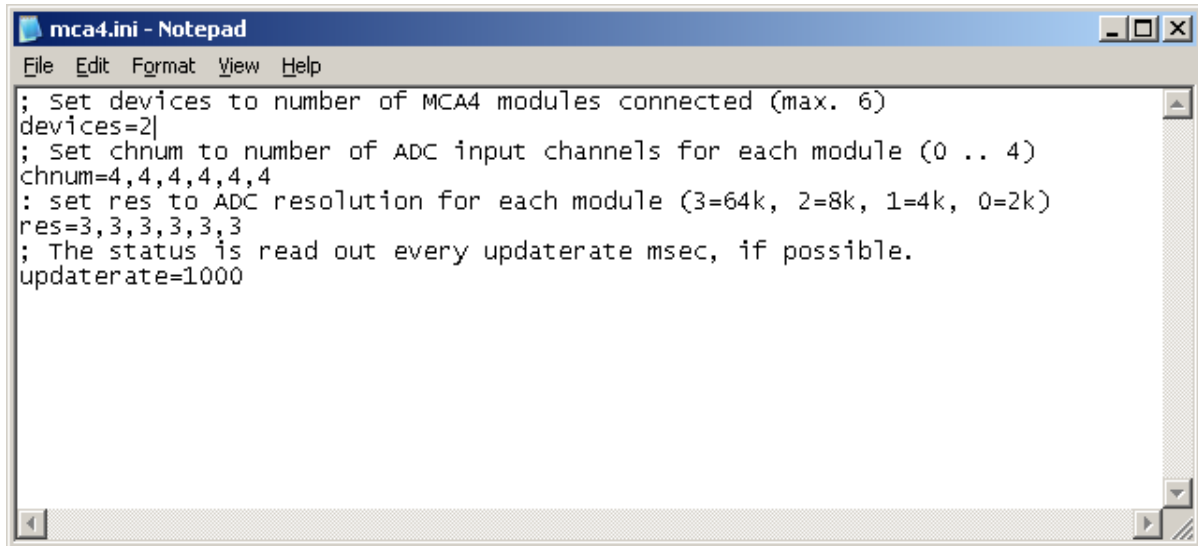


Fig. 2.35: Editing MCA4.INI

That is all you have to do, next time when you start the software all modules will be found and can be operated.

3 Hardware Description

3.1 Overview

The MCA4A series is a family of USB 3.0 (compatible with USB 2.0) connected multichannel analyzers. All settings are software selectable. No jumper (except for the optional ADC full scale ranges), switch etc. configurations are necessary.

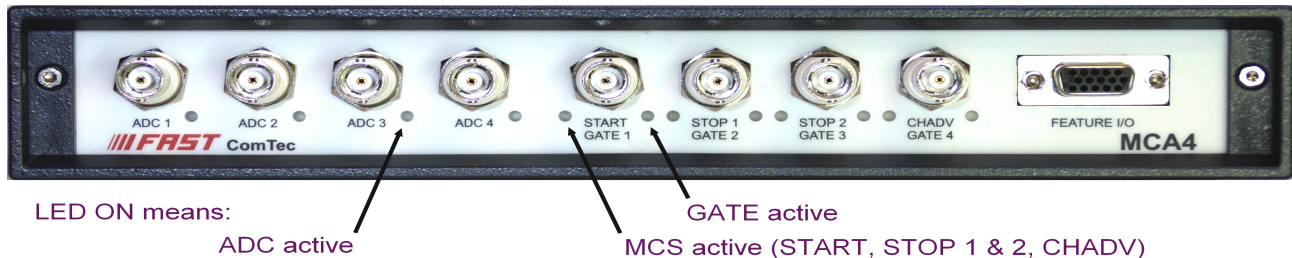


Fig. 3.1: MCA4 front panel

A 1Gbit internal DDR2 RAM is configurable as large FIFO buffer for listmode data transfer or as histogramming memory for onboard spectrum acquisition.

Four onboard discriminators at the four rightmost BNC connectors allow for either positive or negative GATE / MCS input signals. The input voltage range is from -1.2 V to + 5.0V. The corresponding threshold voltages are software selectable in steps of 1mV. The input impedance can be switched between 50 Ohm to GND and 4.7kOhm with either pull-down or pull-up to +3.3V.

Four outputs at the feature connector provide several selectable signals (e.g. SCA, BUSY, etc.) for external usage.

The legendary 'GO'-line (compatible with other FAST ComTec products) allows for easy synchronization or triggering of other measurement equipment.

Furthermore a versatile 8 bit digital I/O port allows for a whole bunch of experimental control, monitoring or whatsoever other tasks. Moreover, the 8 bit digital I/O port can be used for event and spectra tagging to allow for multi-detector configurations, sequential data acquisition etc.

3.2 Front Panel Connectors

3.2.1 Analog ADC Inputs

The analog ADC inputs are 4 BNC type connectors on the front panel of the MCA4A. Each is paired with a green LED that indicates that the input is active. They accept negativ or positiv uni-polar or bi-polar signals.

The input signal is continuously sampled by a 125MHz / 16bit (11, 12, 13 bit options are available) analog-to-digital converter. The converted data is buffered in a 32k x 16bit FIFO for subsequent evaluation. This buffer FIFO can be segmented in 1, 2, 4 or 8 equal partitions. Each segment holding the complete data of a detected input pulse for evaluation. This allows for virtually dead-time-less conversion as incoming pulses can be stored while a previous pulse is still processed by the PHA fitting algorithms. In absolute maximum mode (no fit) actually no conversion induced dead time exists.

3.2.2 Changing the Input range

The full scale voltage range is jumper selectable between 2.0V, 2.5V, 8.0V and 10.0V (default).

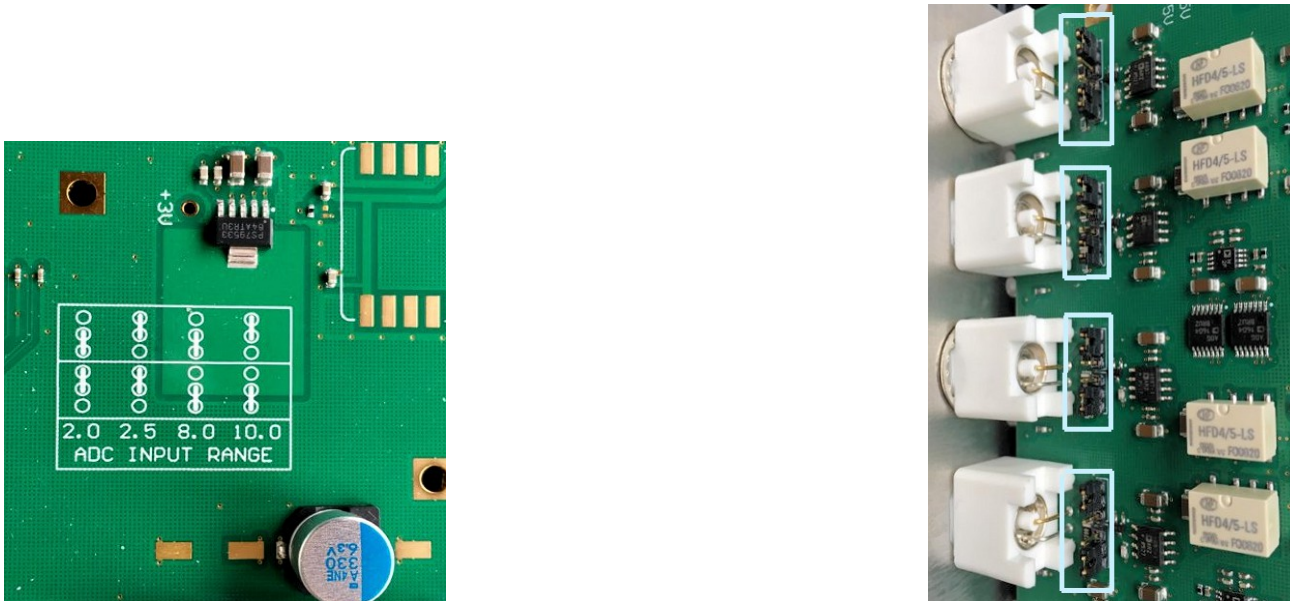


Fig. 3.2: Left: jumper setting for different input ranges, right: jumper positions

The default setting is 10 V corresponding to a VPP value of 11019 (mV) for the full 16 bit range 0..65535. This value can be edited in the ADC settings, if the “advanced Settings” checkbox is crossed, see Fig. 4.11. Standard values for the 8 V, 2.5 V and 2 V ranges are 8875, 2640 and 2120. Enter the new VPP value if you have changed the input range and save it. You can do a precise calibration of VPP if you connect a DC voltage U to the input that you measure by a good digital voltmeter and get the corresponding value V(U) by clicking the “Get Zero” button. The VPP value is then

$$VPP = U * 65535 / (V(U) - V(0)).$$

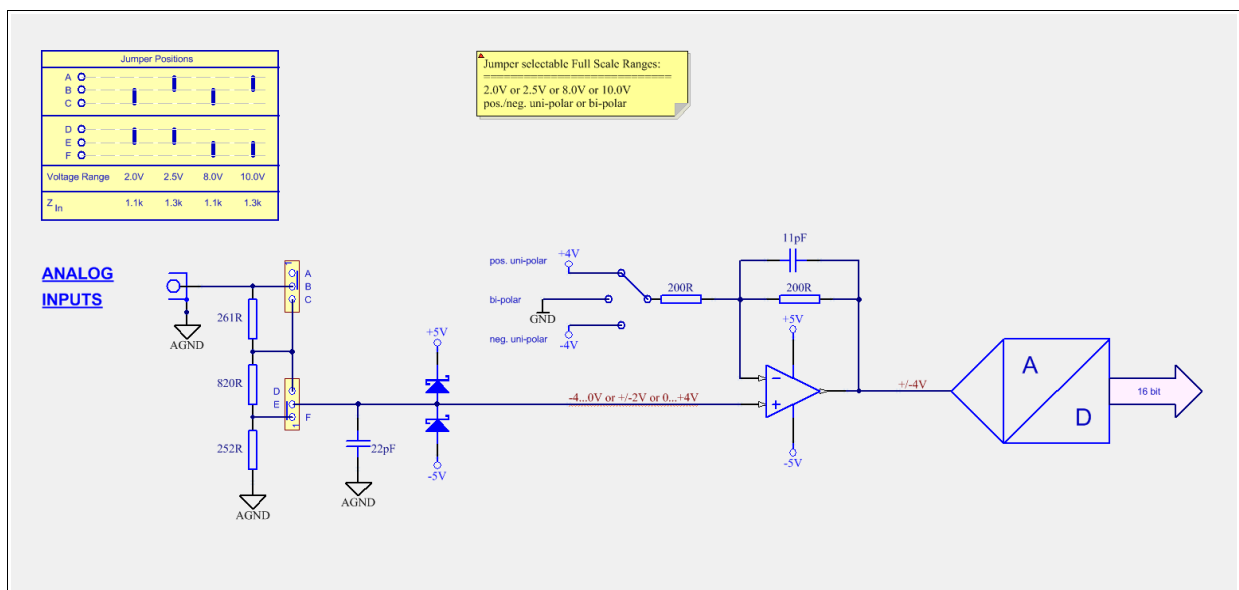


Fig. 3.3: Simplified analog input and ADC circuitry

3.2.3 GATE / MCS Inputs

The GATE / MCS inputs provide a switchable 4.7kOhm pullup to +3.3V, 4.7kOhm pulldown to GND or 50Ohm to GND input impedance. Each input is equipped with a fast discriminator. The threshold level can be chosen from +1.3V for TTL signals or -300mV for fast NIM signals.

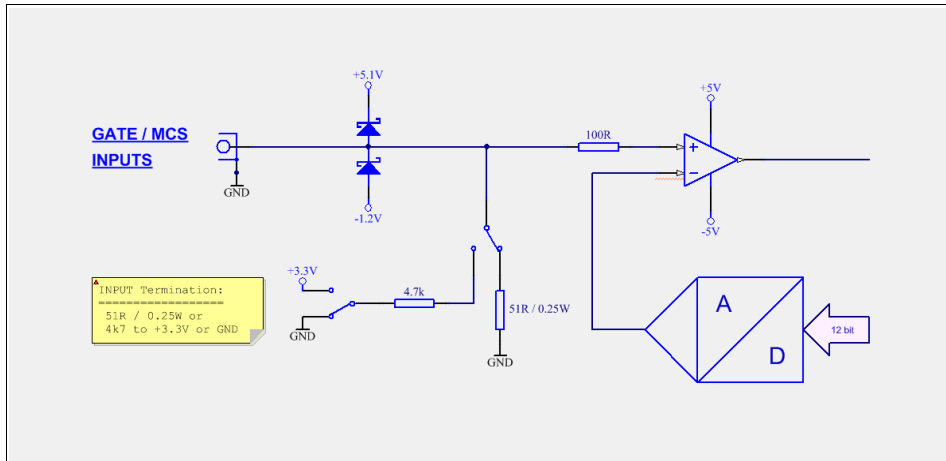


Fig. 3.4: Circuit diagram of each GATE / MCS input

In PHA mode these inputs can be used as GATE. PHA event gating is evaluated at the threshold crossing. COINC and ANTI coincidence is available.

In SVA mode the selected edge of the GATE signal samples the corresponding ADC input signal.

In MCS mode these inputs are used as START, STOP 1 & 2 and external channel advance (CHADV).

NOTE:

For clean switching and oscillation free operation the input slew rate should be higher than $25V/\mu s$.

3.2.4 FEATURE I/O Connector

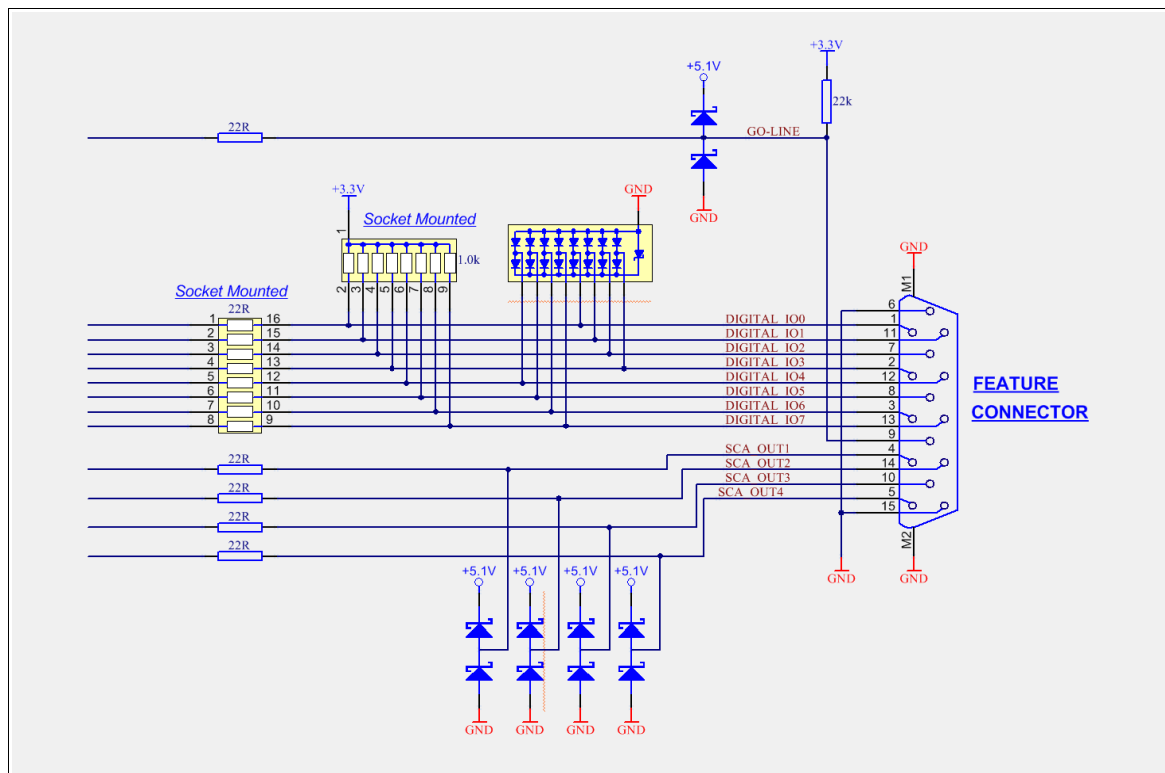


Fig. 3.5: Circuit diagram of the FEATURE I/O port

Single-channel-analyzer (SCA) Outputs

When a maximum¹ of the sampled analog ADC data is detected the corresponding maximum value is evaluated for lower- and upper-level thresholds (LLD and ULD). If the maximum falls in between LLD and ULD the pulse is accepted as valid and a SCA output pulse of 120 ns is initiated.

Other signals can be selected at the nominal SCA ports:

1. SCA (120ns pulse for each valid SCA event)
2. ADC BUSY (over-the-threshold OR input FIFO full → dead time)
3. GATE input signal (as discriminated)
4. MCS event (2.5ns pulse for each detected count edge)
5. SVA event (as internally validated)
6. ...

Digital I/O ports

A very versatile 8 bit digital I/O port is implemented on the FEATURE I/O connector. Since the corresponding resistors are socket mounted they can be easily user configured in a most flexible way.

This I/O port is fully software controllable and each single (1-bit) port is individually configurable. It might be used for external alerts, sample changer control, status input / outputs, spectrum / event tagging, etc.

As can be seen from Fig. 3.6 each individual bit might be configured as input only (tri-stated output), pull-up (open drain) or driver output (small R-I/O) with readback capability. Wired OR / AND connections are feasible.

¹ A peak maximum is detected when the actual input data is 6% lower than the previous maximum value since the threshold crossing

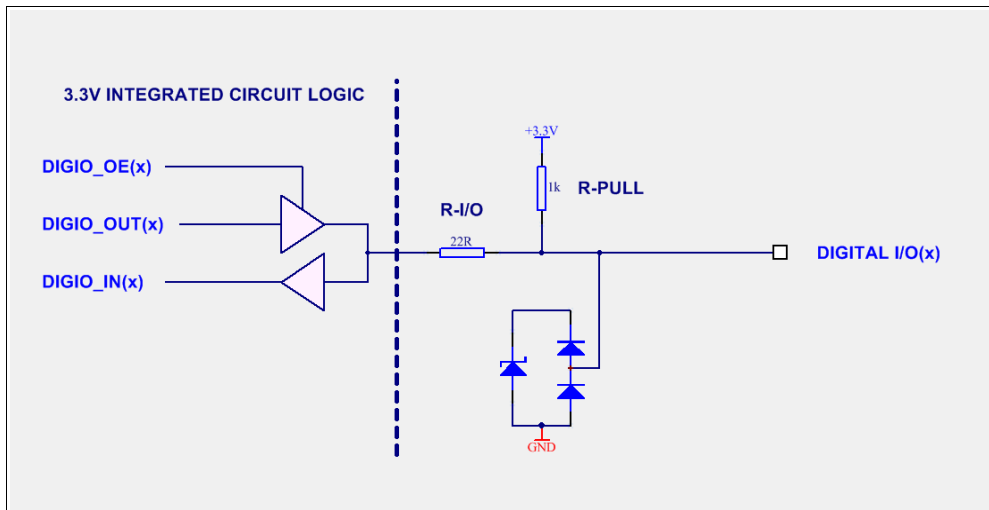


Fig. 3.6: Digital I/O port details

GO-Line

The system-wide open-drain wired-AND 'GO-Line' enables any connected device to start and stop all participating measurement equipment simultaneously. The 'GO' line may be enabled, disabled, set and reset under software control.

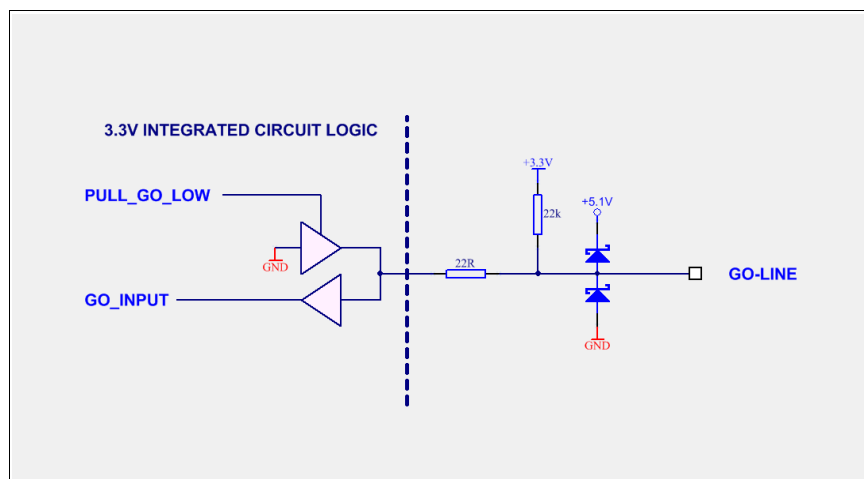


Fig. 3.7: GO-LINE details

3.3 MCA section, Analog-to-Digital Converter

The 16bit analog to digital converter samples the input signals at a continuous 125MHz rate.

The data can be evaluated in several ways. Each input channel might be configured independently.

- Scope mode
- Pulse height analysis (PHA)
- Sampled voltage analysis (SVA)

3.3.1 Scope mode

After a trigger event scope pictures of 4k, 8k, 16k or 32k data samples are acquired.

The trigger source can be chosen from channel 1, 2, 3, 4 or external (START input). The trigger edge (rising or falling) is also selectable.

The available trigger modes are SINGLE, NORMAL and AUTO.

In SINGLE mode only the first trigger event is displayed. In NORMAL mode subsequent trigger events start a new display. In AUTO mode the measurement is triggered by real trigger events or - in absence of those - automatically.

Note:

In scope mode a small data FIFO size should be chosen (ref. Fig. 4.6) for better display response

For further evaluation FFT and software PHA are available (ref. chapter 5 MPANT Software).

Be aware that due to the large number of data being transferred the effective event rate is low compared to a firmware based (onboard) pulse height analysis. Thus, the scope mode is mostly for setting up the experiment.

3.3.2 Pulse height analysis

Onboard (firmware based) signal evaluation and pulse height analysis allows for high event rates and only little conversion induced dead time.

The incoming ADC data stream is evaluated for a rising edge threshold crossing which is taken as the start of a pulse. The input signal is sampled at 8ns intervals and the data is stored in the 32k input buffer FIFO. This 32k buffer FIFO is segmentable in 1, 2, 4 or 8 equal sized segments. This provides for virtually no analysis induced dead-time as a new pulse can be stored while a previous is still processed.

After the input signal fell below the threshold again the pulse sampling is finished and processing begins. Optionally a hysteresis for the threshold crossings can be chosen. This might improve the spectrum for low amplitudes, slowly falling input pulses and high noise levels.

Each PHA event is tagged with the RTC (1ms or optionally 8ns) and – if selected – the TAG inputs at the pulse start. The tagging is executed at the leading edge threshold crossing.

After detecting a completed input pulse the stored data is checked for the minimum and maximum number of samples. This enables to discard false (too short and too long) pulses (ref. Fig. 3.8). Then the data is evaluated using the selected fitting algorithm.

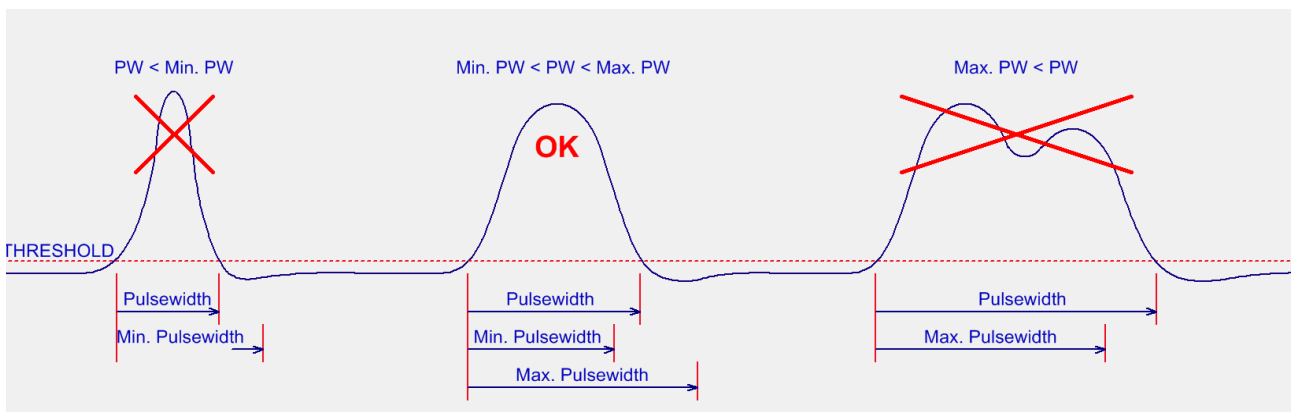


Fig. 3.8: Pulsewidth check - discard "false" pulses

Absolute maximum (none):

No calculation is necessary and the absolute maximum value is taken.

Fit for Gaussian shaped pulses:

A curve is fitted through the selected number of samples ($2 \times i1 + 1$) around the absolute maximum. This provides for a reduction of noise and thus, better peak widths in the acquired spectrum.

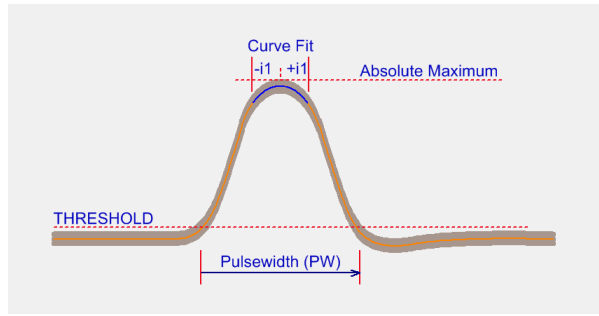


Fig. 3.9: Gaussian shaped pulse fitting

The conversion time (ref. chapter 7.4.2) depends on the number of samples used for the fit.

NOTE:

The optimum number of samples used for the fit depends on the shaping time of the input pulses.

Flat top pulse fitting:

For flattop pulses a chosen number of samples is averaged. The area (start and end) of the averaged data is software selectable.

Three different modes are selectable for defining the averaging area (ref. Fig. 3.10).

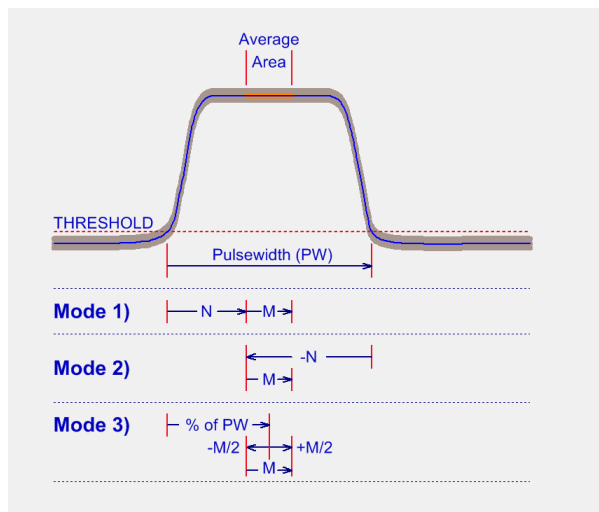


Fig. 3.10: Flattop pulse fitting modes

Mode 1 (left edge) is good for pulse forms that settle near the beginning

Mode 2 (right edge) is good for pulse forms that settle towards the end.

Mode 3 (CFT) is particularly useful when the pulse width is largely dependent on the pulse height.

Pile-up rejection

A pile-up rejection feature is implemented. The incoming data is checked for a possible second peak before the signal fell below the threshold again (ref. Fig. 3.11). If such a second peak is detected the pulse data is marked with a "pile-up detected" flag. It is also possible to discard the event without analyzing it further, i.e. no fitting is performed, when this feature is enabled (ref. chapter 4.1.7 ADC Settings and Fig. 4.11).

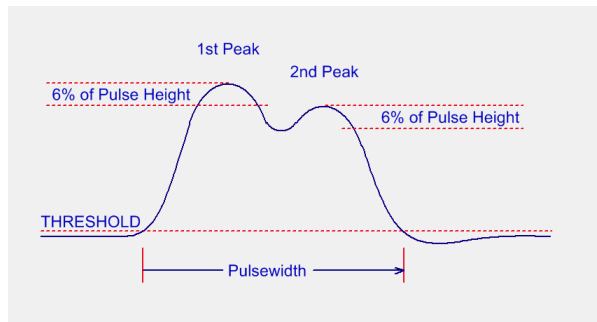


Fig. 3.11: Pile-up detection

A second option to discard pile-up events is to use the GATE input as REJECT signal if the external electronics (e.g. pre-amplifier) provides such a signal. In this mode an active REJECT signal while the ADC input signal is above the threshold (i.e. a pulse is active) discards the very event.

3.3.3 Sampled voltage analysis

In SVA mode the corresponding GATE input is used to initiate the sampling of the current analog input voltage. The rising or falling edge of the GATE signal is selectable.

Simultaneously each SVA event is tagged with the RTC (1ms or optionally 8ns) and – if selected – the TAG inputs.

3.3.4 Single channel analyzer (SCA)

When a maximum of the input signal is detected the corresponding absolute maximum value is checked if it falls in between the lower and upper threshold levels. If so, a 120ns long SCA pulse is generated. This is available at the FEATURE I/O connector (ref. Fig. 3.5).

3.3.5 Real time, live time and system time counter

Each of the 4 ADC channels is equipped with 2 independent and presettable 48 bit counters. Also a 48 bit and presettable system timer is provided.

The timer counters provide a resolution of 1ms (8ns with RTC option).

Each pulse event is tagged with the corresponding system time value when the leading edge threshold occurs.

3.3.6 List mode data acquisition

The analyzed data (pulse height) is patched with additional information like event time, TAG bits or pulse width and transferred to the computer in 64 bit words per event.

3.3.7 Event tagging, spectrum multiscaling

The 8 digital I/O ports can be used for event tagging as well. Each bit is individually selectable. The digital I/O line reading is continuously updated every 2.5 μ s.

In PHA mode of operation the digital I/O port reading is sampled at the leading edge threshold crossing.

3.4 Multichannel scaler

3.4.1 Introduction

In MCS mode of operation the MCA4 counts STOP events relative to the START (trigger) signal. The time resolution is equal to the selected dwell time. The total measurement time is the spectrum length in time bins multiplied by the dwell time. The maximum STOP count rate is 400MHz.

The maximum number of time bins (spectrum length) is 16M. It can be selected in steps of 1. This time granularity allows to optimize the sweep length for the specific experiment.

Histogramming is performed in the onboard DDR2 memory providing the low dwell times (30 ns / 60 ns) and high count rates (400 MHz).

3.4.2 Dwell time

The internal dwell timer provides automatic time bin advancement in the range of min. 30ns (1 channel MCS) or min. 60ns (dual channel MCS) up to $((2^{48} - 1) \times 10\text{ns}) = 781\text{h}$ in steps of 10ns.

External channel advance is supported via the corresponding CHADV input. Please note that the external channel advance signal is internally synchronized to a free running 100MHz clock. An integrated control circuitry assures the minimum cycle time of 100ns, thus guaranteeing no triggering on too fast, ringing or spiky signals.

3.4.3 Sweep counter

The presettable sweep counter is 48 bit wide and increments on each valid START (trigger) event.

3.4.4 Spectrum sequencing

Sequences of spectra consisting of a preset number of sweeps may be automatically generated.

The number of sequences is in the range of 2 to 65536.

When the preset number of sweeps is accumulated the sequence counter is incremented and a new set of sweeps is started.

3.5 Octal 125 MHz Counter / Timer

The MCA4A includes eight 48bit synchronous 125 MHz up-/down-counters. 2 counters are presettable.

All inputs are TTL like and are terminated with 50Ω to ground. For each counter input the polarity is software selectable.

The counters actually either count the software selected rising or falling COUNT input edges (counter mode) or work as 125 MHz timers (RTC option is required).

Counting can be individually gated or switched between up- and down-counting by the individual GATE inputs. Also a common gate is provided. The usage can be individually selected for each channel.

Also a common LOAD is individually selectable. This actually clears the chosen counters.

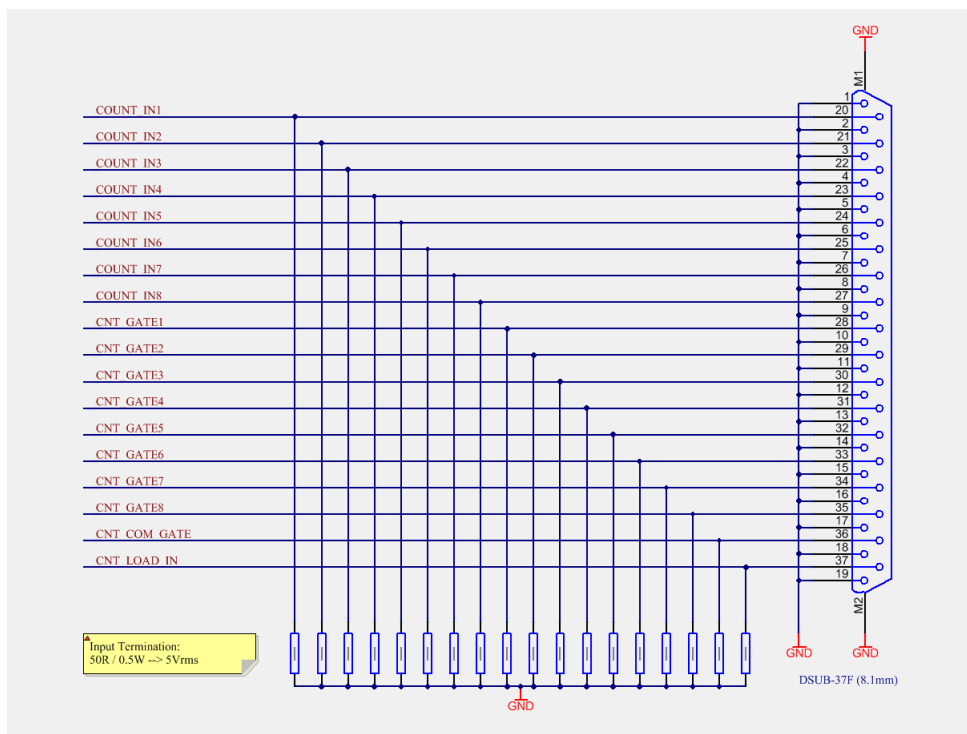


Fig. 3.12: Octal 125 MHz Counter Connector

The system wide GO-Line may also control the counting action. This can be set individually for each counter.

3.6 Time Base / Reference Clock

The internal reference is a 10 MHz crystal oscillator. This can be output at the CLK rear side BNC connector. Or an external clock signal may be used. This external reference frequency might be in the range from 5 to 100 MHz in steps of 1 MHz.

4 Software Description

The software is a layer of several parts. A driver (WINUSB.SYS) provides the communication via the USB port. A server program (MCA4A.EXE) controls the hardware and the data. A graphical user interface (MPANT.EXE) controls the server via a DLL (DMCA4.DLL) and provides graphical displays. This chapter describes now the server program.

The window of the MCA4 server program is shown here. It enables the full control of the MCA4A to perform measurements and save data. This program has no own spectra display, but it provides - via a DLL ("dynamic link library") - access to all functions, parameters and data. The server can be completely controlled from the MPANT software that provides all necessary graphic displays. The view of the server window depends on the number of active input channels, here all 4 ADC channels are active.

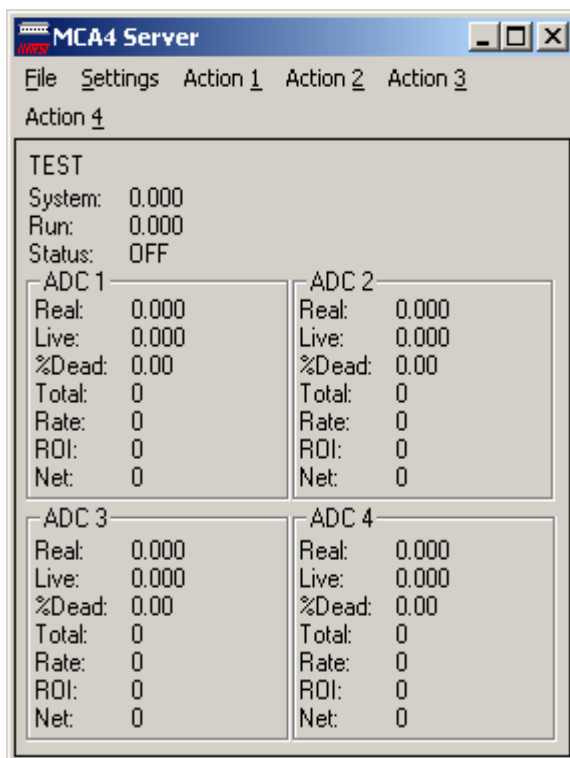


Fig. 4.1: MCA4 Server Window

4.1 Server functions

To start the software, just double click a shortcut icon linking to the server program. The server program performs a test whether the hardware works well on this computer, then starts MPANT and gets iconized. Usually you will control everything from MPANT, but it is possible to work with the server alone and independently from MPANT.


Fig. 4.1 shows the opened server window. The upper part of the window contains parameters that belong to all ADCs. The filename (TEST), "System" is the runtime of the system as counted by hardware, "Run" means the run time of the acquisition in seconds as measured by the PC clock, "Status" is the acquisition status (ON / OFF / READ OUT). The lower part contains in boxes labeled with the ADC name (ADC1, ADC2, ..) counters belonging to each ADC. "Real" is the real time in seconds, measured by the hardware. "Live" is the active time counted from the start of the acquisition, i.e. the time of pulse processing subtracted. "%Dead" means the percentage of the dead time (time of pulse processing) relative to the total time. "Total" is the number of counts in the spectrum as counted by the software. "Rate" is the number of counts per second within the ROI (region of interest) that is defined in the settings. The ROI is usually spanned over the total spectrum and rate is then the total count rate. "ROI" is the sum of counts within the ROI, default value is the total sum. "Net" is the net sum of the ROI, a background is subtracted as calculated by linear interpolation from the spectrum points at the ROI boundaries.

4.1.1 Initialization files

At program start the configuration files MCA4A.INI and MCA4A.SET are loaded.

Parameters that can be set by editing the MCA4A.INI file are the update rate in msec for the refresh of the status, and the voltage calibration factors. The number of ADC channels and resolution is set automatically.

The file MCA4A.SET contains the default settings. It is not necessary to edit this file since it is saved automatically. Instead of this .SET file any other setup file can be used if its name without the appendix 'A.SET' is used as command line parameter (e.g. MCA4 TEST to load TESTA.SET). The command line parameter can be used either by calling the program from a "DOS box" or by changing the properties of an icon on the desktop that points to the program.

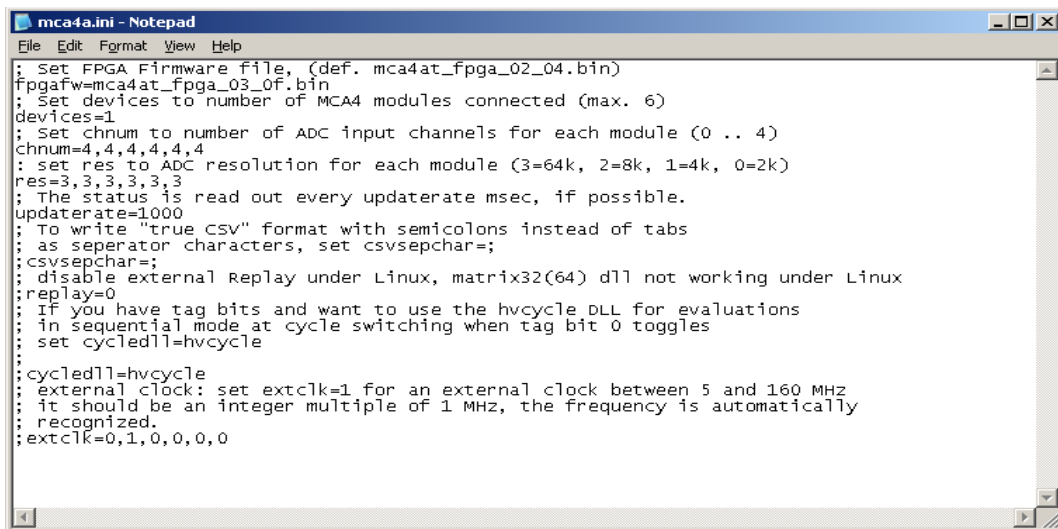

```
mca4a.ini - Notepad
File Edit Format View Help
; Set FPGA Firmware file, (def. mca4at_fpga_02_04.bin)
fpgafw=mca4at_fpga_03_0f.bin
; Set devices to number of MCA4 modules connected (max. 6)
devices=1
; Set chnum to number of ADC input channels for each module (0 .. 4)
chnum=4,4,4,4,4,4
; set res to ADC resolution for each module (3=64k, 2=8k, 1=4k, 0=2k)
res=3,3,3,3,3,3
; The status is read out every updatarate msec, if possible.
updatarate=1000
; To write "true csv" format with semicolons instead of tabs
; as separator characters, set csvsepchar=;
csvsepchar=;
; disable external Replay under Linux, matrix32(64) dll not working under Linux
replay=0
; If you have tag bits and want to use the hvcycle DLL for evaluations
; in sequential mode at cycle switching when tag bit 0 toggles
; set cycledll=hvcycle
;
;cycledll=hvcycle
; external clock: set extclk=1 for an external clock between 5 and 160 MHz
; it should be an integer multiple of 1 MHz, the frequency is automatically
; recognized.
extclk=0,1,0,0,0,0
```

Fig. 4.2: MCA4A.ini File

4.1.2 Action menus

The server program normally is shown as an icon in the taskbar. After clicking the icon it is opened to show the status window. Using the "Start" menu item from the Action1 menu a measurement can be started. In the status window every second the acquired events, the counting rate and the time are shown. Clicking the "Halt" menu item the measurement is stopped and via "Continue" proceeded. As it is possible in the System Settings dialog to create up to four systems each containing different MCA's combined for simultaneous starting and stopping, there are four Action menus Action1..Action4.

4.1.3 File menu

Clicking in the File menu on the Data... item opens the Data Operations dialog box.

This dialog allows to edit the data format settings and perform operations like Save, Load, Add, Subtract, Smooth and Erase.

The Radio Buttons **MPA**, **Selected Spectra** and **New Spectra** provide a choice between handling of the complete data set (MPA) or selected spectra, or to load new selected spectra for compare.

Mark the checkbox **Save at Halt** to write a MPA file containing the configuration and all spectra at the stop of a measurement. The file name can be entered. If the checkbox **auto incr.** is crossed, a 3-digit number is appended to the file name that is automatically incremented with each saving. The format of the data can be ASCII (extension for separated spectra .ASC), binary (.DAT), GANAAS (.SPE), EMSA (.EMSA) and CSV (.CSV).

If **Separate Header** is not checked, the Header and data is saved together in a file with extension .MP, otherwise the file with extension .MP contains only the header and the data is written separately into a file with appropriate extension.

The buttons **Save**, **Load**, and **Erase** perform the respective operation. With **Add** and **Sub** spectra can be added or subtracted from the present data. The checkbox **calibr.** can be checked to use a calibration and to shift the data then according to the calibration. The **Smooth** button performs an n-point smoothing of selected single spectra. The number of points to average can be set with the **Pts** edit field between 2 and 21. Check the **Write Listfile** checkbox to write a listfile during a run. **No Histogramming** prevents calculating any spectra to save computing time and concentrate the system on writing the list file.

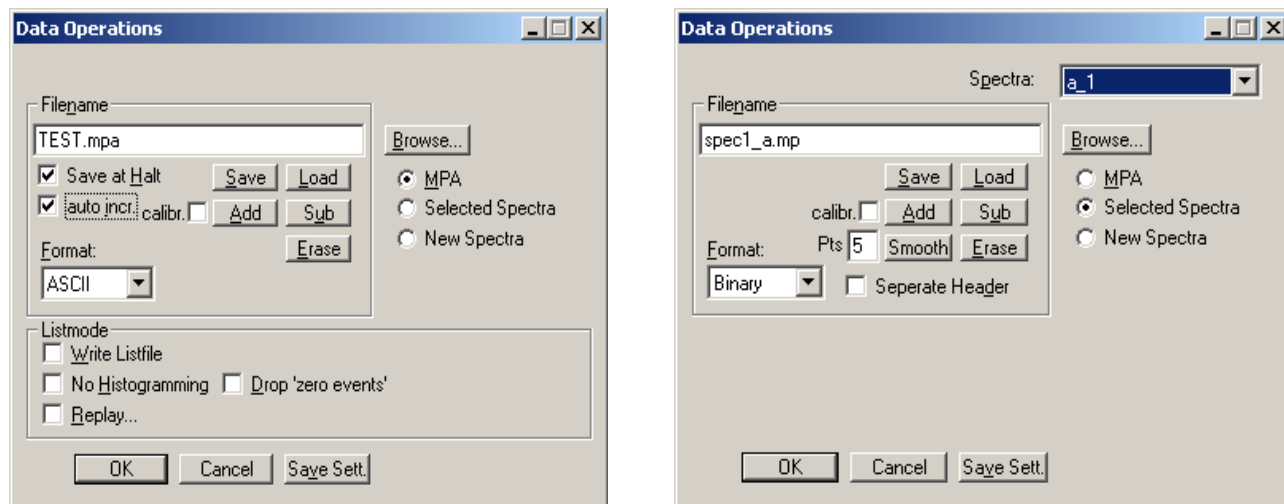


Fig. 4.3: Data Operations dialog for MPA data (left) and selected spectra (right)

The menu item **File – R**eplay... opens the Replay dialog. It enables evaluation of list data. The Replay software licence is usually programmed into the MCA4 module or by an optional USB dongle.

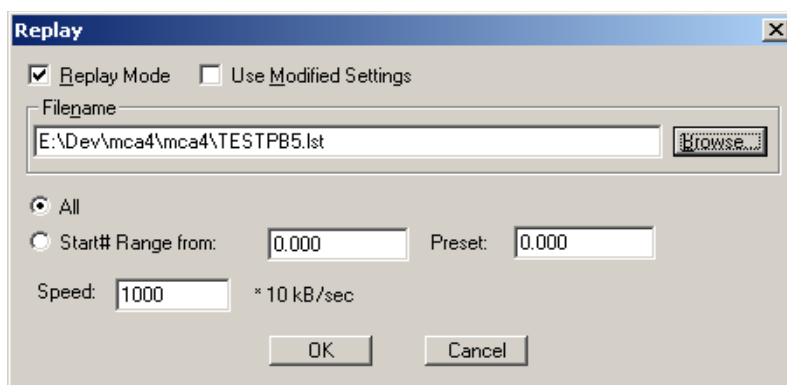


Fig. 4.4: Replay Settings dialog

Enable **Replay Mode** using the checkbox and specify a **File name** of a list file (extension .LST) or search one by pressing **Browse...** With the radio buttons it is possible either to choose the complete list file by selecting **All** or a selected **Start# Range**. Specify the sweep range by editing the respective edit fields **from:** and **Preset:** . The **Replay Speed** can be specified in units of 100 kB per sec. To **Use Modified Settings** enable the corresponding checkbox; otherwise the original settings are used. To start Replay press then Start in the Action menu or the corresponding MPANT toolbar icon.

The menu item **File – A**bout... opens the About MCA4 window where some information of the System can be found. Particularly the serial number is important for support purposes. This serial number is unique for each MCA4 system. A click on the **Installed...** button opens a dialog showing all installed features.

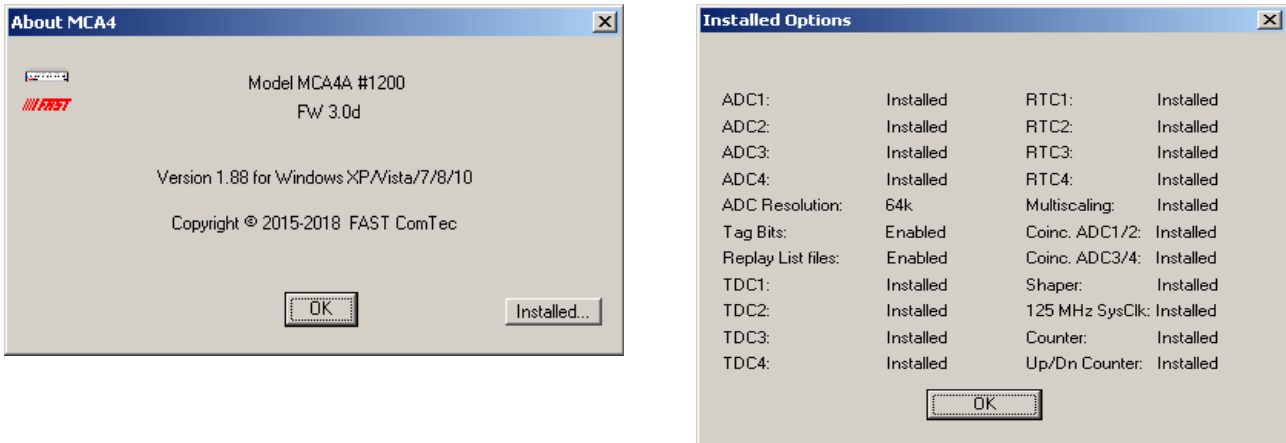


Fig. 4.5: About MCA4 dialog box and Installed items

The MPANT menu item in the file menu starts the MPANT program if it is not running.

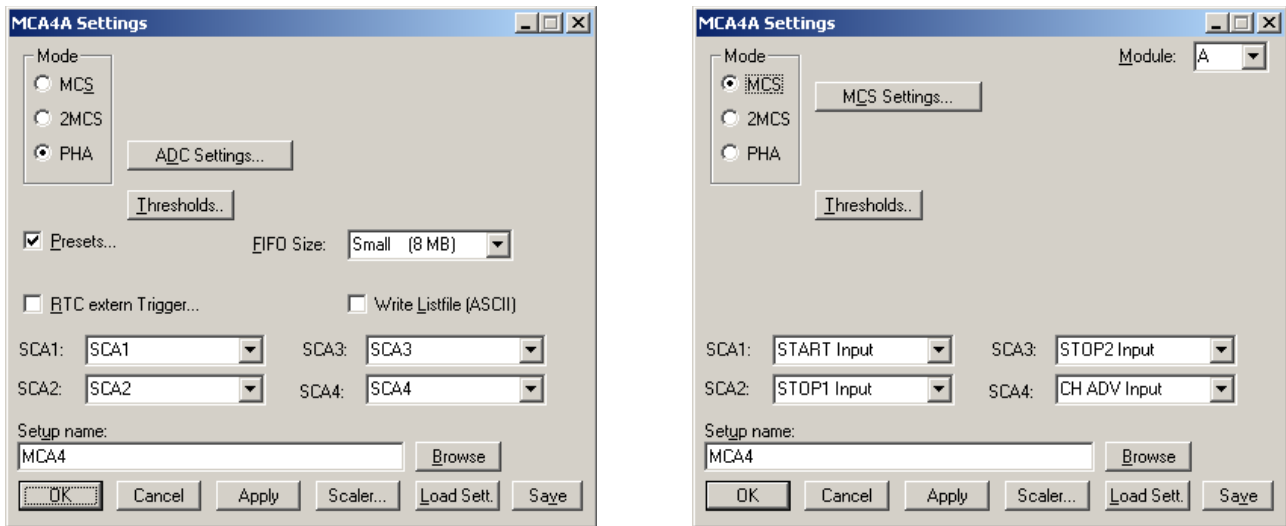


Fig. 4.6: Settings overview, left: PHA mode, right: MCS mode, more than one MCA4 module.

4.1.4 Settings Overview

The Hardware... item in the Settings menu opens the MCA4 Settings dialog box. It enables to switch between PHA mode (left figure) and MCS mode (right). If more than one MCA module are installed, the module (A, B,...) can be selected in the upper right corner as shown right in Fig. 4.6.

Via the **SCA** - combo boxes you can specify the signal you want to observe at the feature connector. See chapter 3.2.4. You can observe for each channel the SCA, ADC BUSY (over-the-threshold), GATE Input, MCS Input or SVA event signals.

In PHA mode a List file can be written by checking the corresponding checkbox **Write Listfile**. The data format can be ASCII or binary depending on the format settings in the data operations dialog. The FIFO size can be selected between Full (128 Mb), Medium (32 Mb) and Small (8 Mb). For real data acquisition the full size is recommended. For adjusting parameters when using Scope mode, a small FIFO size is useful as changes can be observed immediately whereas with a full FIFO size they are delayed by several seconds.

OK takes all settings and closes the dialog. **Apply** sets all settings without closing the dialog. **Cancel** cancels all changes. Pressing **Save** writes all settings into the file **MCA4A.SET** and **MCA4_###.set** resp. the entered setup name. The file name is built by attaching "A.set" and "_###.set" to the shown name, where "###" is the serial number of the MCA4. **Browse** lets you search for settings files, **Load Sett.** loads one. When loading a settings file, the file with the serial number fitting to the device in the file name is loaded if one is found, otherwise or when the software is running in demo mode the file MCA4A.SET for the first MCA4 module.

4.1.5 MCS Settings

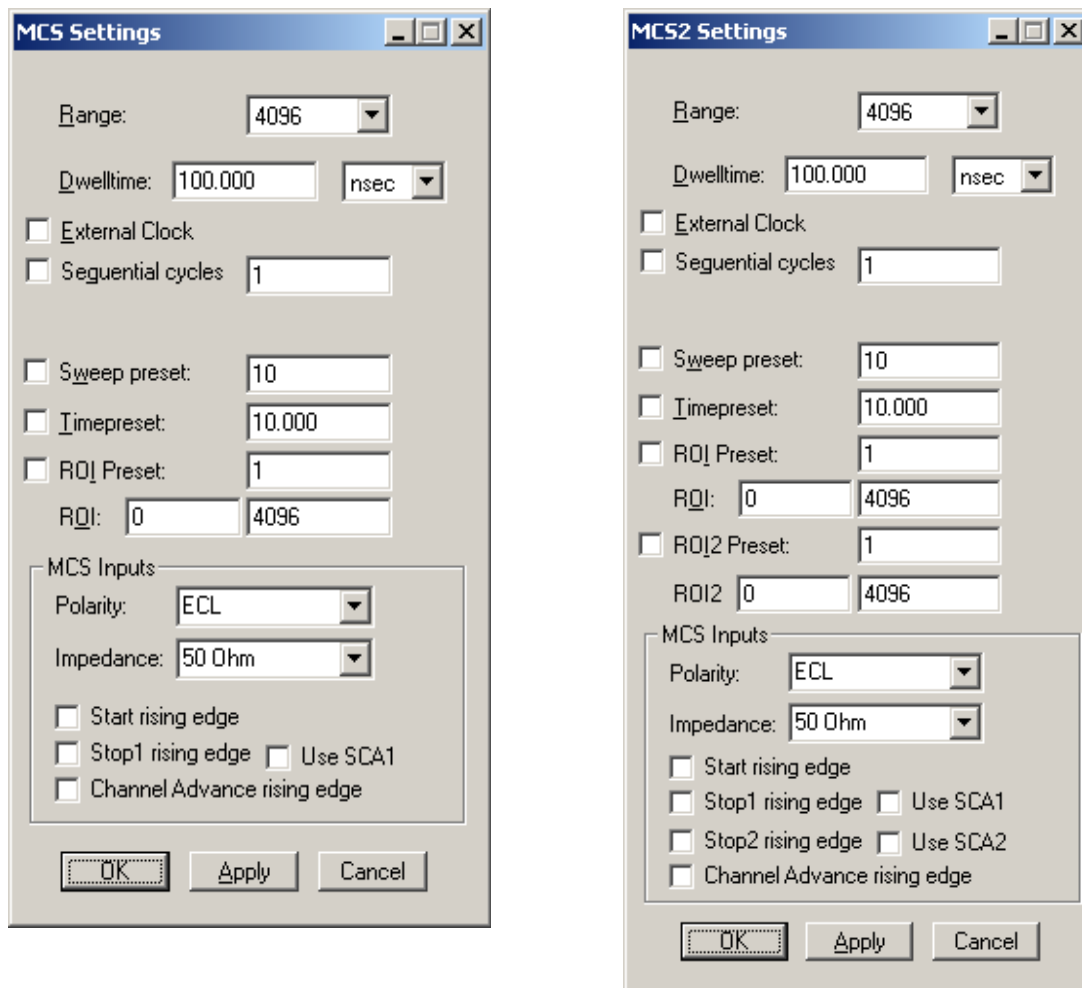


Fig. 4.7: MCS Settings for one (left) or two stop inputs (right)

Clicking the **MCS Settings** button opens the MCS Settings or MCS2 Settings dialog as shown in Fig. 4.7. The **Range** defines the number of time bins, the **Dwelltime** can be set in units of 10 ns from 30 ns (MCS) or 60 ns (MCS2). The measurement can be stopped automatically after a specified number of sweeps by checking **Sweep preset**. If the checkbox **ROI Preset** is marked, the measurement will be stopped after acquiring more events than specified in the interval between the boundaries defined in the corresponding edit fields. The events are counted only if they are within the ROI limits, i.e. \geq the lower limit and $<$ the upper limit. Another possibility is to acquire data for a given time via the **Time Preset**.

A series of measurements can be acquired into separate memory parts by checking **Sequential cycles** and specifying the number of cycles. Each single measurement should be terminated preferably by a sweep preset condition. The complete run stops after performing the specified number of cycles. The spectra are then shown in a two dimensional view row by row in the MPANT program. As the sequential mode is implemented in the firmware, there is only a very small delay of 30 ns between the time spectra when using a sweep preset.

For some measurements it is desired to synchronize the sequential cycles to an external signal from a sample changer or target mover. This can also be done by using the Digital I/O port together with the "waitpin" command. Just define any of the function toolbar icons in the MPANT program by a command like **oncycle waitpin 3000** and click on this toolbar icon. Now in sequential mode when a sweep preset is reached after each individual cycle, the acquisition stops until the voltage level at Dig I/O 7 gets low. If the voltage is still high after 3000 msec, a message box will pop up and warn that a time out happened. See section 4.3 about more details and on the syntax of control language commands.

For the MCS Inputs the **Polarity** can be selected between TTL and ECL (NIM) signals. For each of the MCS inputs you can select between the rising and the falling edge. Usually for TTL signals, rising edge should be active and for (negative) fast NIM signals “rising edge” should be deselected as the leading edge has the better timing. Since version 2.56 instead of the Stop1 or Stop2 inputs also SCA1 or SCA2 can be chosen by crossing the **Use SCA1 / Use SCA2** checkboxes, this may be especially interesting for Mossbauer measurements to reduce background.

OK takes the settings and closes the dialog. **Apply** sets the parameters in the hardware without closing the dialog. **Cancel** closes the dialog and cancels all changes.

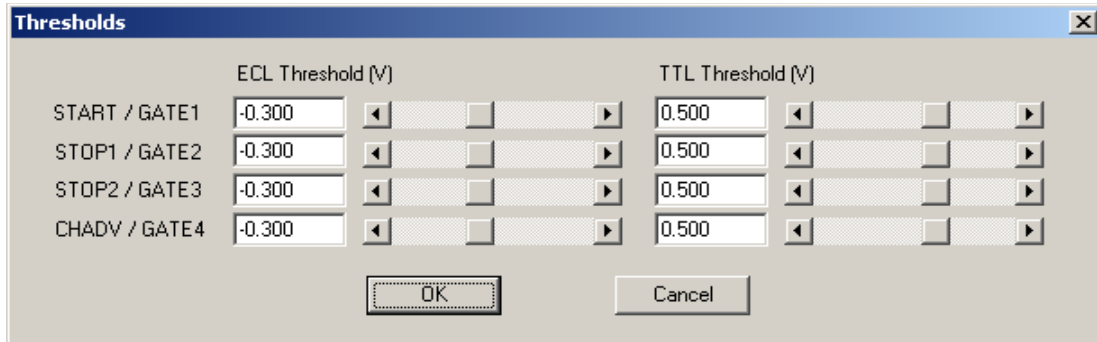


Fig. 4.8: Thresholds Dialog

The input thresholds can be set individually for each input. Click the “Thresholds...” button to open the Thresholds dialog.

4.1.6 Presets

Clicking the Presets checkbox in the Settings dialog opens the ADC Presets dialog as shown in Fig. 4.10. Here all preset parameters for PHA mode and SVA mode can be set. Select the ADC in the drop-down listbox at top right. If the checkbox ROIpreset is marked, the measurement will be stopped after acquiring more events than specified in the corresponding edit field. The Counts are events counted by the hardware that trigger a pulse height analysis, whereas the ROI events are counted in the software if the corresponding data in the spectrum are within the ROI² boundaries, i.e. \geq the lower limit ROI min and $<$ the upper limit ROI max that can be defined in the ADC Settings dialog. Another possibility is to acquire data for a given run time via the RTimepreset or a given live time via the LTimepreset. A measurement will be stopped if the corresponding checkbox is marked.

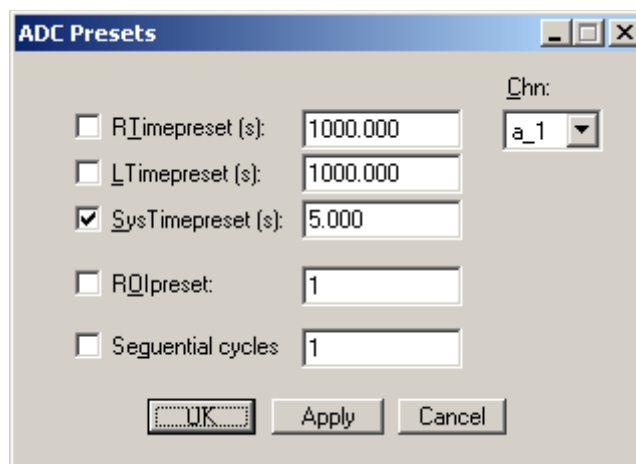


Fig. 4.9: ADC Presets Dialog

2 ROI: Region Of Interest

4.1.7 ADC Settings

Clicking the ADC Settings button in the Settings overview the ADC Settings dialog as shown in Fig. 4.10. Here all parameters like range, etc. for PHA mode and SVA mode can be set. Select the ADC in the drop-down listbox at top right. The Enable checkbox must be checked if the corresponding channel is active. The ADC resolution and length of the spectrum is set by the Range listbox between 2 and 65536 in binary steps. For the Gate input there is a choice between Don't Use, COINC, ANTI and REJ IN.

COINC means that the voltage level at the Gate input must be high to trigger an event, ANTI means that the level must be low. A pull-up resistor in COINC mode and a pull-down resistor in ANTI mode lets the input work if no gate input is connected, see . The gating is evaluated at the time of the threshold crossing. Set All sets the range and Gate mode for all ADC's.

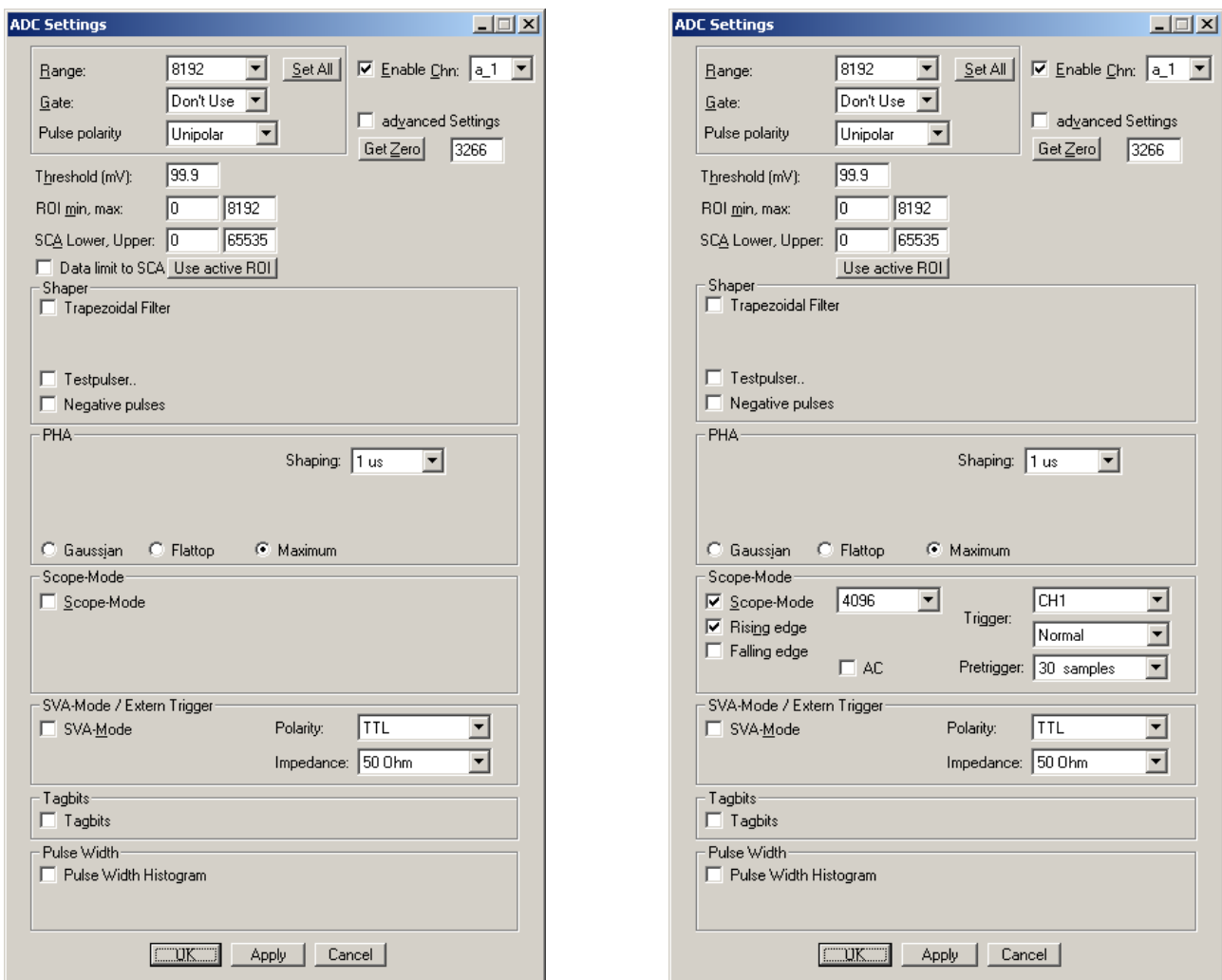


Fig. 4.10: ADC Settings dialog, right: Scope mode

The SCA³ boundaries can be set with the SCA Lower, Upper edit fields. Here the units are in LSB's of the 16-bit ADC, i.e. full range is from 0 to 65535. By pressing the button **"Use active ROI"** the boundaries are recalculated from the ROI boundaries. If an event is recorded inside the SCA boundaries, a pulse signal can be output at the corresponding monitor output at the feature connector. A checkbox **"Data limit to SCA"** limits data transfer to values inside the SCA, i.e. use the SCA as lower and upper level discriminator.

The Threshold for triggering an event can be set in mV, an event will be triggered if the voltage level crosses the zero plus threshold level. The Get Zero button calibrates the offset voltage level. It is recommended to connect no signal when doing this offset calibration. The offsets will be saved into the default settings file

3 SCA: Single-channel-analyzer

MCA4A.SET by clicking “Save” in the hardware settings dialog, and automatically reloaded when the software is started next time. See chapter 2.5.1. As the offset calibration can cause a shift in the spectrum, be careful not to do a new “Get Zero” within a series of measurements. The pulse height spectrum shown in the MPANT program is an histogram of the ADC data, the offset is subtracted. Therefore the maximum pulse height value that can be recorded when using a 64k range is 64k minus the zero value. The maximum spectrum range is therefore somewhat less than the binary range set, so the spectrum length shown in the MPANT program is correspondingly reduced and the spectrum will be new painted when a new “Get Zero” calibration is done.

The MCA4A is able either to record the waveform and send the complete waveform data (Scope Mode) or to do a pulse height analysis and just send the result. The maximum count rate is of course much lower in Scope mode, therefore this mode is used for optimizing the electronics like adjusting a pole zero in the amplifier and also selecting the parameters for the PHA (pulse-height analysis). See chapters 2.5.2 and 2.5.3.

Select Scope-Mode to inspect your signal as shown in the right part of Fig. 4.10. New elements appear now in the dialog for editing the waveform length and the trigger parameters. The waveform length can be selected between 4k, 8k, 16k and 32k. The unit is 8 ns, so the waveform length is between 32 μ s and 262 μ s. For the Triggering you can set the Source and Trigger Mode. For the first check use the channel corresponding to the ADC as trigger source, i.e. CH1 for ADC a_1, and Normal as trigger mode. Now close all dialog boxes with OK and do a short acquisition by clicking the Start toolbar icon in the MPANT program. You should get a view of the waveform of your pulses. If you now go back into the ADC Settings dialog, a new button “guess” will appear beside the Shaping listbox in the PHA box. If you click that button, a recommended value for the shaping time will be selected to get a first choice for the PHA parameters.

The trigger source in Scope-Mode can be selected between Software, CH1, CH2, CH3, CH4 and Extern (START). “Software” means that an event will be triggered by the software as often as possible. “Extern Trigger” means a signal connected at the START / GATE1 input. When using the Extern trigger new edit elements appear in the “SVA-Mode / Extern Trigger” box for the Polarity, Impedance and Edge of the START input discriminator.

The Trigger Mode can be Normal, Auto and Single Shot. “Auto” means that a software trigger will be done if no trigger event is detected. “Single Shot” means that the acquisition is stopped after a single trigger event. To do this, an event count preset of 1 is set. The Pretrigger can be selected between none, 30 samples, 100 samples and 250 samples.

Check SVA-Mode if you want to do a sampled voltage analysis. In this mode an external trigger signal samples a voltage conversion. The same trigger parameters like polarity, impedance, edge can be set as for the external trigger. For ADC1 the START / GATE1 input is used, for ADC2 the STOP1 / GATE2 input and so on. SVA-mode disables automatically Scope-Mode and PHA, in the PHA box “Maximum” will be selected.

The 8 Digital I/O ports (Feature connector) can be used as Tag bits if this checkbox is selected and the respective licence is enabled. A new edit field for the number of tag bits (select between 1 and 8) appears in the Tagbits box and a drop-down select box for One Hot, Binary and Don't show. “Don't show” means that the tag bits are only in the list data but not shown.

“One hot” means that a two dimensional spectrum is acquired with the tag bits as y-coordinate. See Fig. 4.11. For eight tag bits there will be 10 rows of the data array. The first row contains the spectrum if no tag bit was detected. The second row is the spectrum with tag bit 0 set only, the next row is the spectrum with tag bit 1 set and so on, last row contains the spectrum when more than one tag bit is set.

If Sequential mode is selected in the Presets dialog, there is an option “next cycle with tag 0 toggle”. This enables to acquire very short measurements synchronized to an external signal without any deadtime, for example to acquire PIXE images. See in the readme.txt at Server version 1.62 about a DLL that can be called at every cycle switching to do fast evaluations. An example named hvcycle.dll is part of the shipped software, the source code is part of the DLL software interface.

“Binary” means the y-coordinate is the value if the tag bits are seen as a binary number. In this case for 8 tag bits the y-dimension of the 2D spectrum is 256. Tag bits can be used in PHA mode, Scope mode and SVA mode. As the digital I/O port is used, the other functions using the digital I/O port like external controlled Start / Stop of an acquisition are not available together with Tag bits.

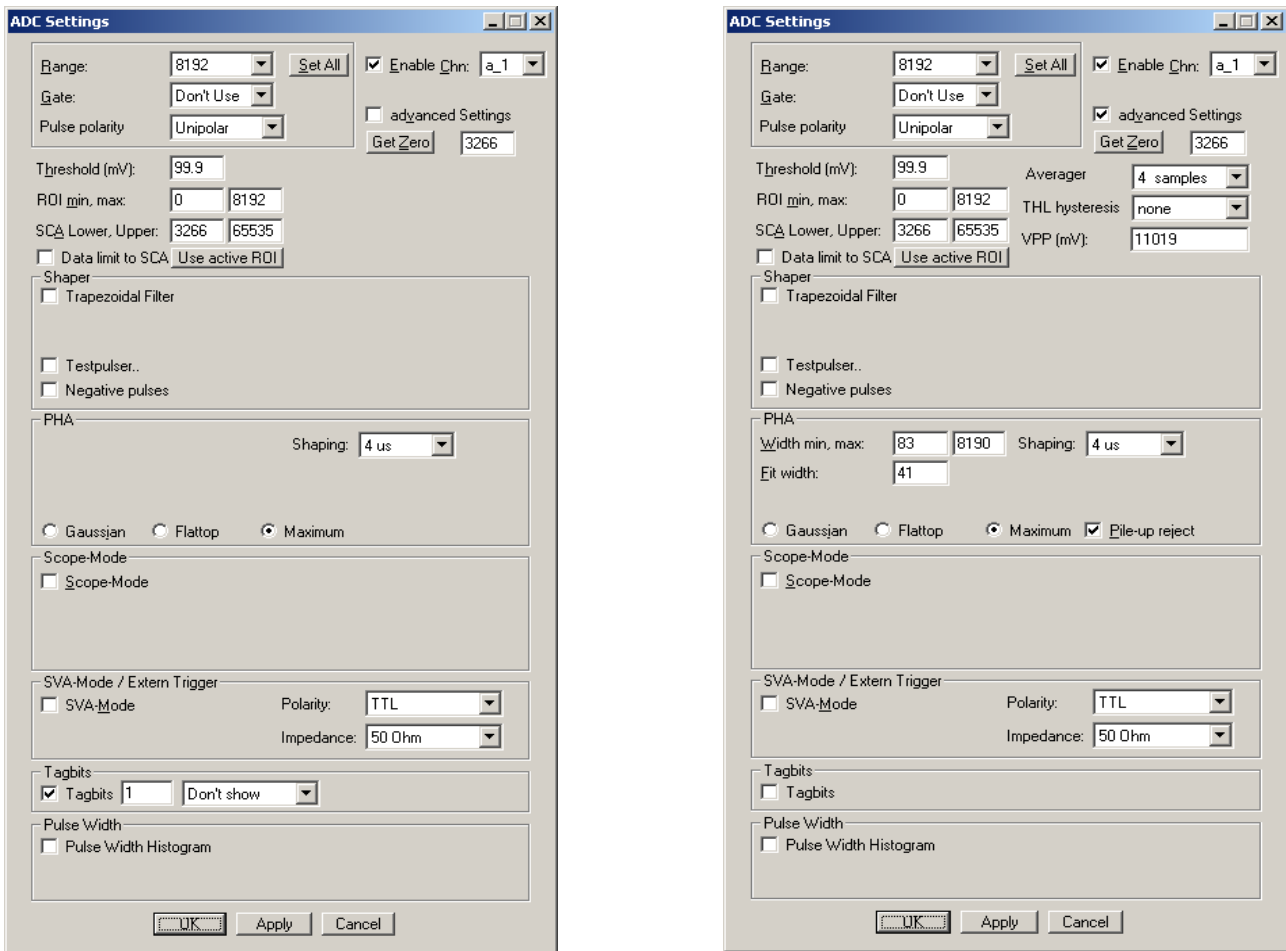






Fig. 4.11: Left: Setting Tag bits, right: advanced Settings

If you check advanced Settings more options will appear for the PHA parameters. VPP (mV) is the voltage corresponding to the full range of the ADC from 0 to 65535 LSB's. You can change it if you have a calibrated value for example by connecting a DC voltage to the ADC input that is known very well for example by using a good digital volt meter. The THL hysteresis can be selected between “none”, “15 LSB”, “30 LSB” and “50 LSB”. The value means that the voltage level at the input must go below the threshold by the respective value before a new pulse will be triggered. The threshold level hysteresis prevents triggering again on spurious signals when crossing the threshold after a pulse. The recommended value depends on the noise level, a value of 30 LSB is recommended, especially for signals that have a slow falling time and when using the “Maximum” method.

The Fit width is the most important parameter that influences the quality of the pulse-height analysis of Gaussian shaped pulses, see Fig. 3.9. This Fit width and the limits Width min, max for the fit range are set when choosing a Shaping time from the list box. You may play with the Fit width parameter and see if you can improve the resolution of your pulses by setting a value corresponding to a shaping time in between of those in the list. If your pulses are shaped with a rather long shaping time, it could be necessary to increase the max Width of the Fit range to get a successful PHA result.

For Flattop pulses it is necessary to specify a ROI at the top of the waveform for calculating the mean value as PHA result, see Fig. 4.12. To get such a zoomed picture of the pulse at its top region, first switch off the Auto scale and Minimum auto scale in the Scale options dialog of the MPANT program (Fig. 4.13). Then set Y-Slice for the Region Shape by clicking the corresponding toolbar icon  and while keeping the right mouse button pressed, drag a zooming rectangle over the interesting region at the top of the waveform, then click the “Zoom in” icon . To define the ROI, switch back to X-Slice for the Region Shape using the corresponding icon . Then drag the ROI boundaries from left to right keeping the right mouse button pressed and click the “Create ROI” icon . To take over the ROI parameters into the flattop PHA, click the from ROI button shown in Fig. 4.14. There is a choice of anchoring the ROI either at the left edge of the pulse, the right edge, or the left ROI boundary at a constant fraction between right and left edge, by selecting CFT in the listbox. The edges are determined by the time when the signal crosses the threshold. It depends

on the electronics creating the flattop pulses, which method of anchoring the ROI gives the best result for a varying pulse height.

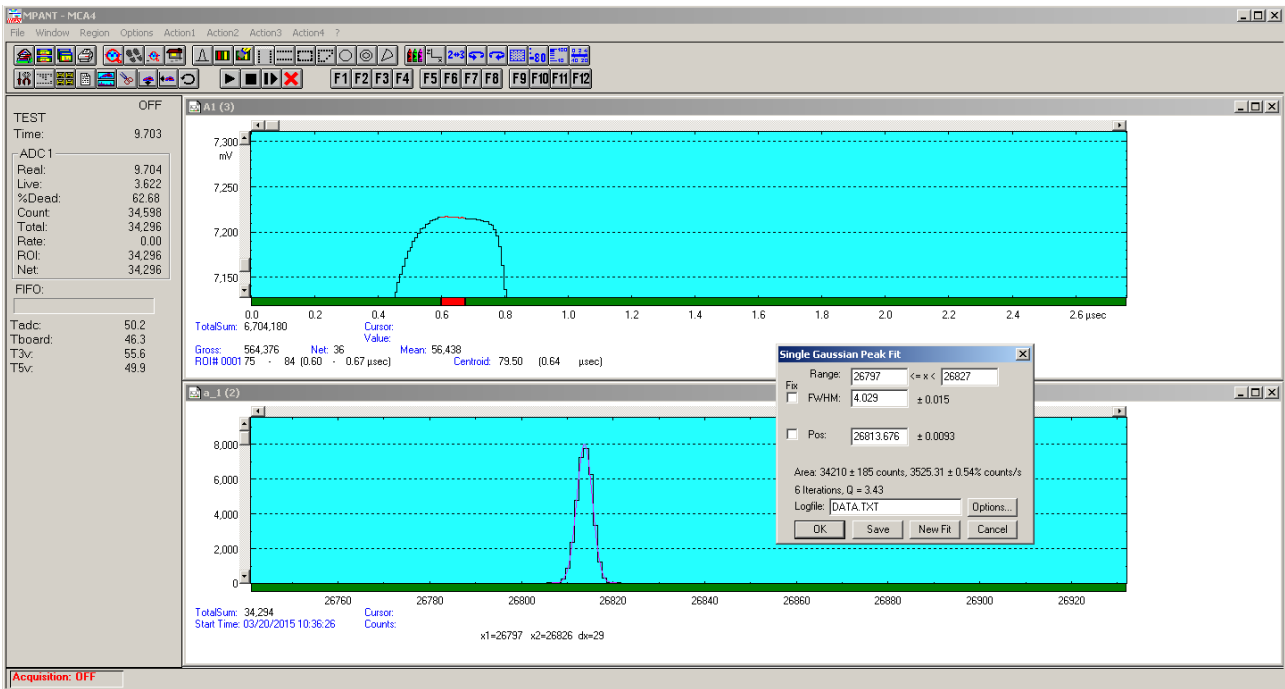


Fig. 4.12: Setting a ROI in the waveform for the pulse-height analysis of flattop pulses

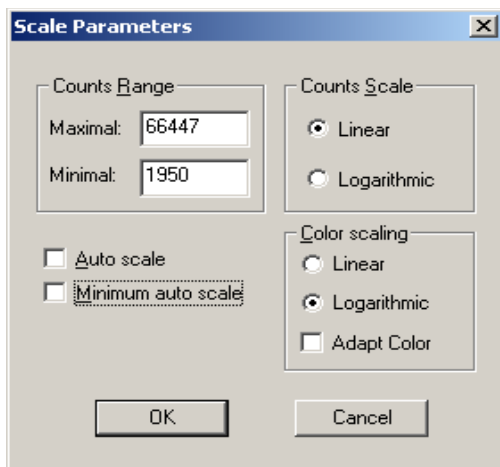


Fig. 4.13: How to switch off Auto Scaling

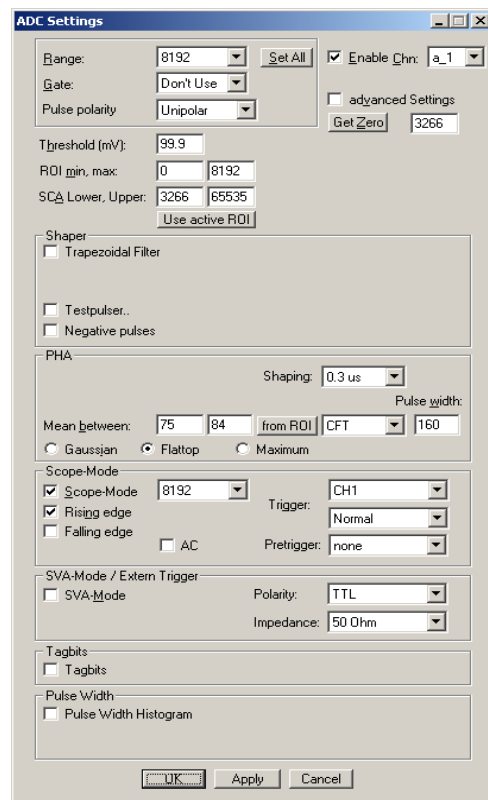


Fig. 4.14: Setting parameters for flattop PHA

4.1.8 Pulse-width histogramming

The pulse width of the analog input pulse at the threshold level can be measured in units of 64 ns, the pulse width data is inserted into the highest 12 bits of the RTC data, the length of the timestamp data is then only 32 bits instead of 44 bits. A histogram of the pulse width distribution can be acquired, and a two dimensional spectrum showing the pulse height versus pulse width can be created. It is possible to reject events if the pulse width is larger than a specified limit. See in Fig. 4.15 how to create a pulse width histogram and set a limit for the maximum pulse width. Note that if "Reject" is enabled, the data is still in the list file and the event is rejected only in the histograms. So when evaluating a list file, the pulse width limit still can be changed.

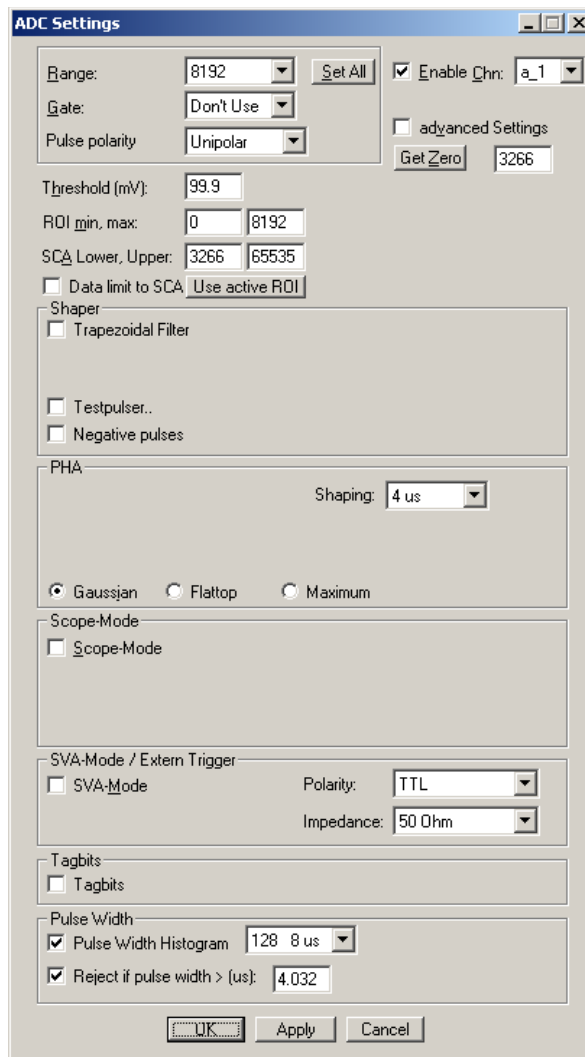


Fig. 4.15: Pulse Width Histogramming

4.1.9 Scaler Setting

The “Scaler Settings...” button in the settings overview opens the respective dialog box to define the use of the 8 125 MHz Scaler/Counters. This is an option, see chapter 3.5 . A utility program mca4scal.exe is started automatically and shows the scaler values, optionally as a ratemeter. The scalers can be individually named, the values are saved into the data file. For an exact synchronizing of start and stop of the scalers with the acquisition it is recommended to enable **Count gated with GO** and **Remove GO when Preset reached**. Furthermore, the three checkboxes concerning the GO-line in the system definition dialog should be crossed, see chapter 4.1.10 .

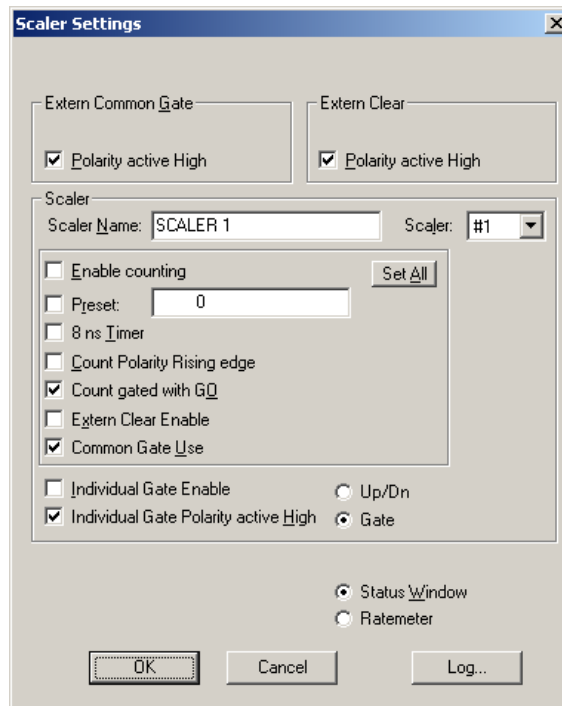


Fig. 4.16 Scaler Settings dialog box

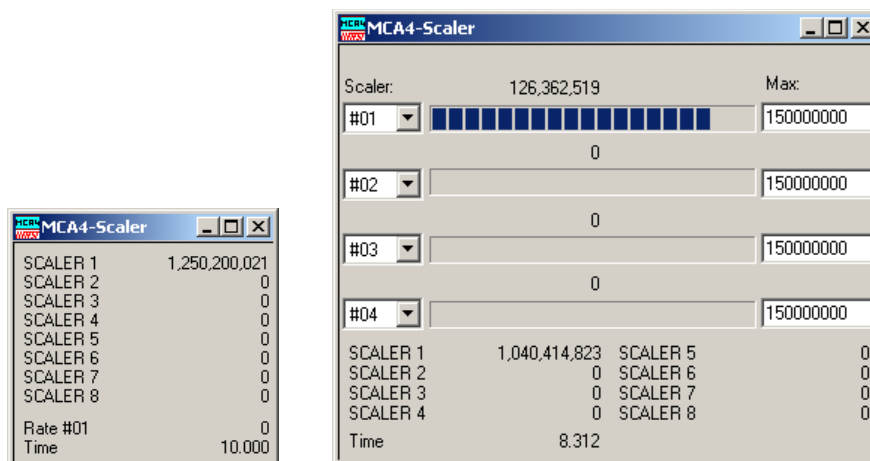


Fig. 4.17 Left: MCA4-Scaler Status utility, right: as ratemeter

4.1.10 System Definition

The "System..." item in the settings menu opens the System Definition dialog box. Several MCA's can be combined to form up to 4 separate systems that can be started, stopped and erased together. The use of the Digital Input / Output and the GO-Line can be defined: It can be used either to show the ON or OFF status of the systems if the checkbox Status Dig 0..3 is marked. At the respective pins +5 Volt are output if an acquisition is running and 0 V if not. Alternatively, it can be used for example with a sample changer by checking "Value inc. at Stop". Here, the 8-bit value entered in the edit field (a number between 0 and 255) is output at the Dig I/O port. This value will always be incremented by 1 if the MCA4 is stopped. The Invert checkbox allows inverting the logical level. See also the commands pulse and waitpin how to handshake a sample changer or synchronize sequential cycles with a target mover, see the control language in 4.3. The Radio buttons Push-Pull and Open Drain define the output mode of the Dig I/O ports.

It is also possible to use the digital input 4 as an external trigger for starting the system (DESY control line). If the corresponding checkbox Start with Dig 4 is marked, a start command for the respective system will not immediately start the system. After the start command, the digital input will be permanently checked for its logical level. If the level changes from high to low, the data for the system is cleared and it will then be started. It will stop if the level returns to high (or vice versa if Invert is marked) and can again be restarted with the next level change. A stop command for the system will finish the digital input checking. By checking Clear before Start the spectra is cleared before the start. As an alternative to controlling Start/stop via digital input 4 also the GO-line can be used by checking Start with GO. This is better as the GO-line controls the hardware directly.

The Use of the GO-Line is controlled via the 3 checkboxes Watch, High when ON, and Low at Stop. The GO line gates directly the hardware. "Low at Stop" means that the GO line is immediately pulled down by the hardware when a sweep preset is reached, or as fast as possible by the software at any other stop. If Tag bits are used, the digital I/O port is used as an input for the tag bits exclusively, see chapter 4.1.7.

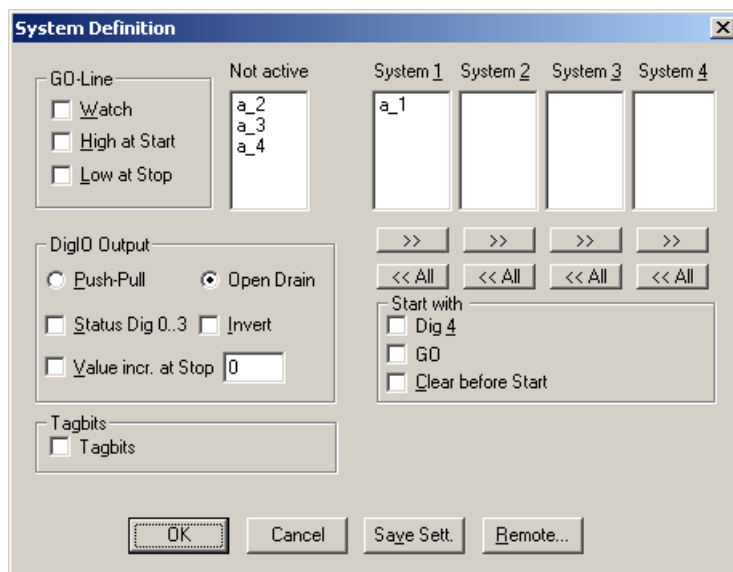
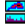


Fig. 4.18: System Definition dialog box

4.1.11 Calculated Spectra dialog

Dependent on installed options like RTC and Coincidence Mode the MCA4 software allows to create several calculated spectra from different channels. By clicking "Spectra..." in the Settings menu or the corresponding MPANT toolbar icon  the Dualparameter and Calculated Spectra dialog box of the MCA4 server is called. Here additional dual parameter and calculated spectra can be added, modified and deleted from the list. Since firmware 5.16 it is possible to create a 2D view of the pulse height versus the pulse width if pulse width histogramming is enabled, by a "Add Shape" button.

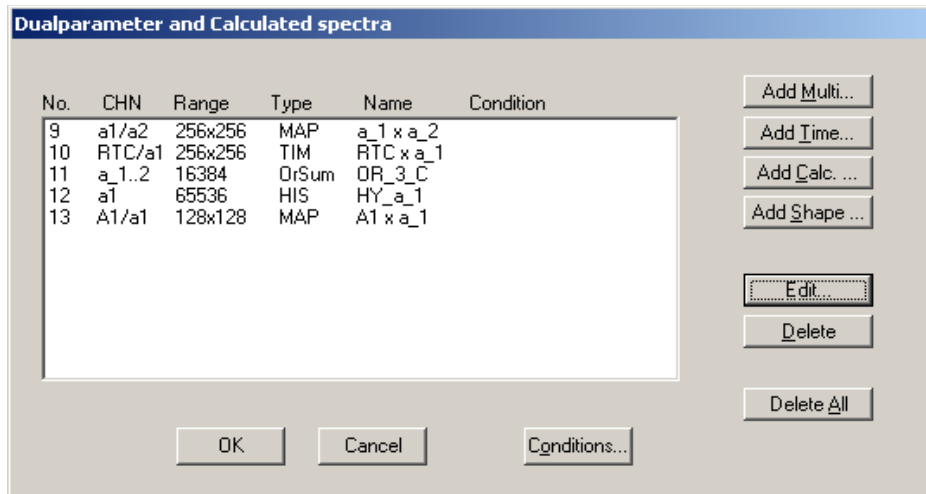


Fig. 4.19: Dual parameter and Calculated spectra dialog box

If the 8 ns RTC for at least two channels and Coincidence mode options are installed, there is a "Add Multi" button. If pressed, the **Multi Display Settings** dialog opens. It allows to define dual parameter spectra showing coincident events from two ADC inputs. The coincidence is recognized by the software if events from different ADC have a RTC value that is inside a specified coincidence time interval.

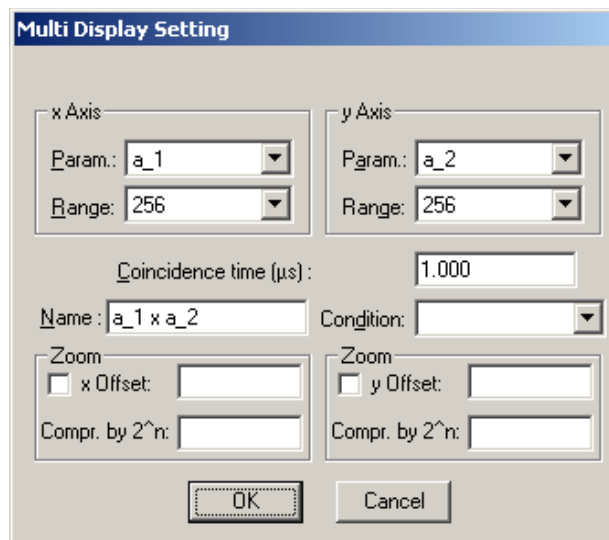


Fig. 4.20: Multi Display Setting

For the x Axis and y Axis a Parameter and Range must be chosen from the respective drop down list. You can edit a Name of the new multi spectra, or you can left the 'Name' edit field untouched to choose an automatically generated name like a_1 x a_2. The length of the Coincidence time window can be defined in units of μs . If Conditions are defined, one can be selected as a necessary Condition for incrementing a channel in this dual-parameter spectra. See more about Conditions later in 4.1.12. For zoomed Maps check the x Offset and y Offset checkbox and enter a value for the respective Offsets and Compression. The Offset

point marks the lower left corner of the zoomed map, and the Compr. by 2^n means the power of two by which the spectra are compressed. A value of zero means full resolution.

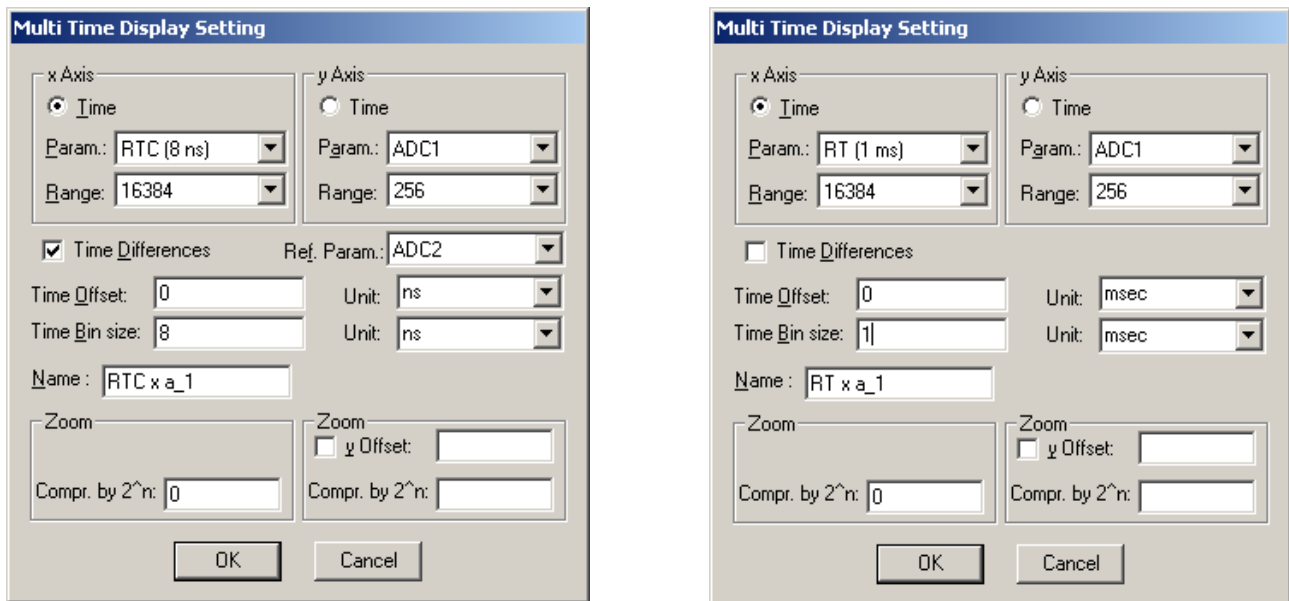


Fig. 4.21: Multi Time Display Setting, left with installed RTC option and right without RTC option.

Press Add Time to define a new Time spectrum in the Multi Time Display Settings dialog. For the x Axis and y Axis a parameter and range must be chosen from the respective drop down list. The time axis can be either the x Axis or the y Axis, select the respective radio button. The time parameter can be either a 1 msec Real time or the 8 nsec RTC clock if the respective option is available. Define the Time Offset and Time Bin size in the corresponding edit fields and Unit selection boxes. A Name can be entered for the new time spectra, or you can left the 'Name' edit field untouched to choose an automatically generated name like RTC x a_1. The checkbox Time Differences enables recording time spectra relative to an event in any selected ADC. For zoomed Maps check the x Offset or y Offset checkbox and enter a value for the respective Offsets and Compression. The Offset point marks the lower left corner of the zoomed map, and the Compr. by 2^n means the power of two by which the spectra are compressed. A value of zero means full resolution. For a one-dimensional time spectrum just use 1 for the ADC range.

Press Add Calc from the Map and Calculated Spectra dialog to define a new calculated single spectrum in the Calculated Spectrum Setting dialog. Depending on installed RTC and Coincidence options, this dialog comes up either as shown left or right in Fig. 4.22.

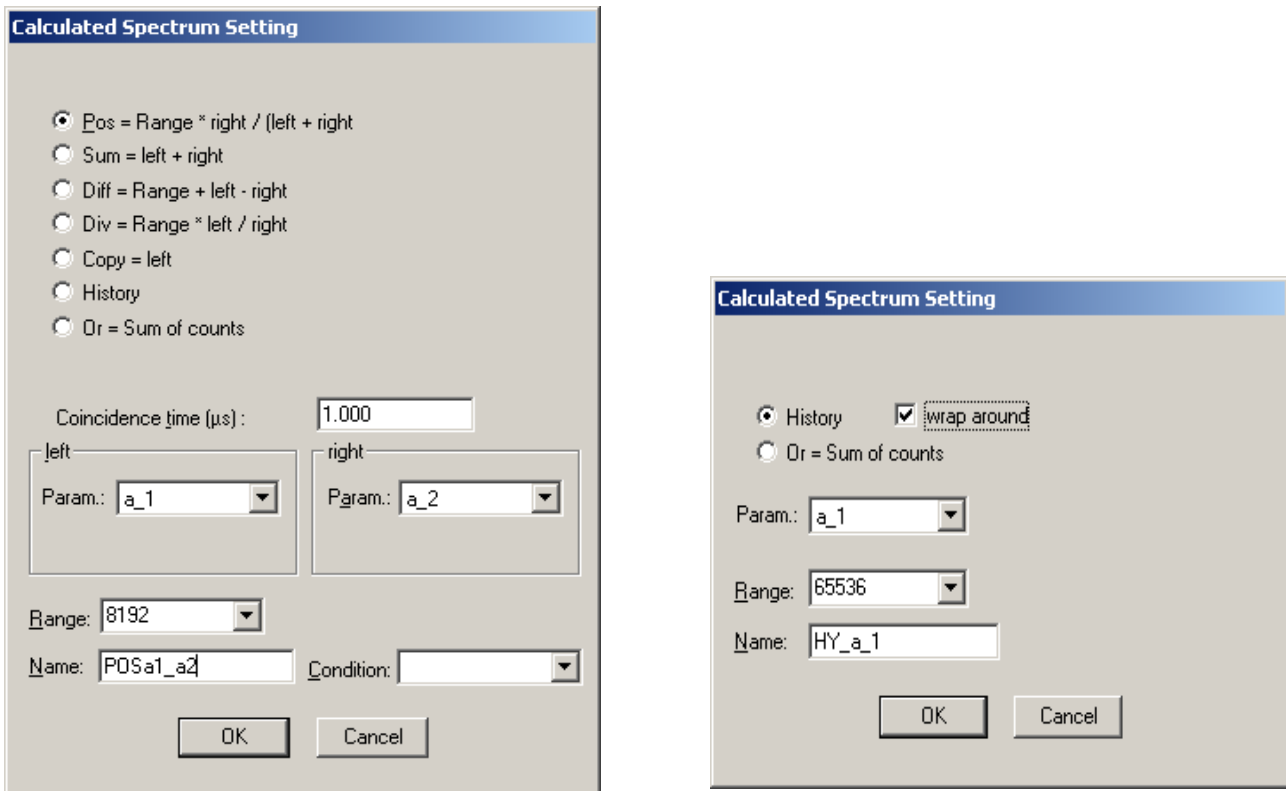


Fig. 4.22: Calculated Spectrum Setting, with (left) and without (right) RTC and Coincidence options

You have a choice between several formulas to combine two parameters: $\text{Sum} = \text{left} + \text{right}$ makes a sum spectrum, and $\text{Diff} = \text{Range} + \text{left} - \text{right}$ can be used to subtract two spectra. The Range and Name can be defined in the edit fields or default values will be taken. Any such calculated spectrum can be used as a parameter for a new dual parameter spectrum. If Conditions are defined, one can be selected as a necessary Condition for incrementing a channel in this spectra. For ADC inputs in addition the $\text{Pos} = \text{Range} * \text{right} / (\text{left} + \text{right})$ formula is available, that is often used for position dependent detectors read out from both sides, and the Division formula $\text{Div} = \text{Range} * \text{left} / \text{right}$ can be used to divide two spectra. The Range and Name can be defined in the edit fields or default values will be taken. Any such calculated spectrum can be used as a parameter for a multi parameter spectrum. If Conditions are defined, one can be selected as a necessary Condition for incrementing a channel in this spectra. The COPY type of spectra allows to set conditions for copies of primary ADC spectra. History allows to display a history plot of an ADC as a new spectrum, i.e. you can see the data of the ADC as they come in versus the event number. The length of the history plot can be defined free, just edit the Range number. Wrap around lets the history plot start new from the left overwriting the old curve, when the right end is reached.

Spectra containing the sum of counts from several ADC's can be used to create a "superdetector" from multiple smaller detectors. The sum spectrum is calculated from the complete single spectra. Select the "Or = Sum of counts" radio button.

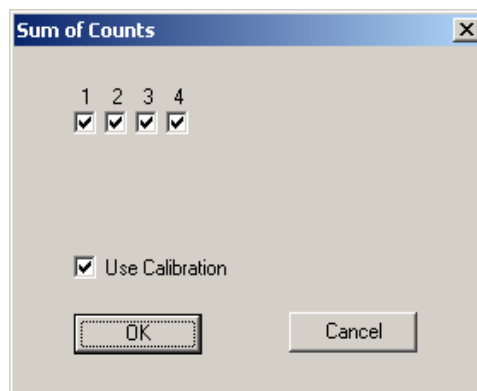


Fig. 4.23: Sum of Counts dialog

In the "Sum of counts" dialog you can select which ADCs will contribute to the superdetector. It is recommended to enable "Use Calibration" and to perform a careful calibration of all selected ADC's:

Make a short acquisition using a pulser or a calibration source so you get a peak in each used ADC. In each ADC zoom into the spectrum, then by keeping the right mouse button pressed drag from left to right over the peak to mark a fit region and click on the "Fit" toolbar icon to make a Gaussian peak fit. Note the position value of the peak in the first ADC. For the first ADC you can just enable "Use Calibration" in the calibration dialog and use the standard parameters $p_0=0$, $p_1=1$. In the next ADC enter two calibration points: for the first point enter 0 for the Channel and 0 for the value and click "Add>>". For the second calibration point use the result of the fit: Click "Fit" to enter the Channel value of the fitted peak position. For the "Value" enter the fit result of the peak position in the first ADC. Click "Calibrate" and enable "Use Calibration", then click "OK". Proceed this way until all used ADCs are calibrated. You will immediately get the effect: Instead of a peak multiplet you will see a single sharp peak in the sum spectrum.

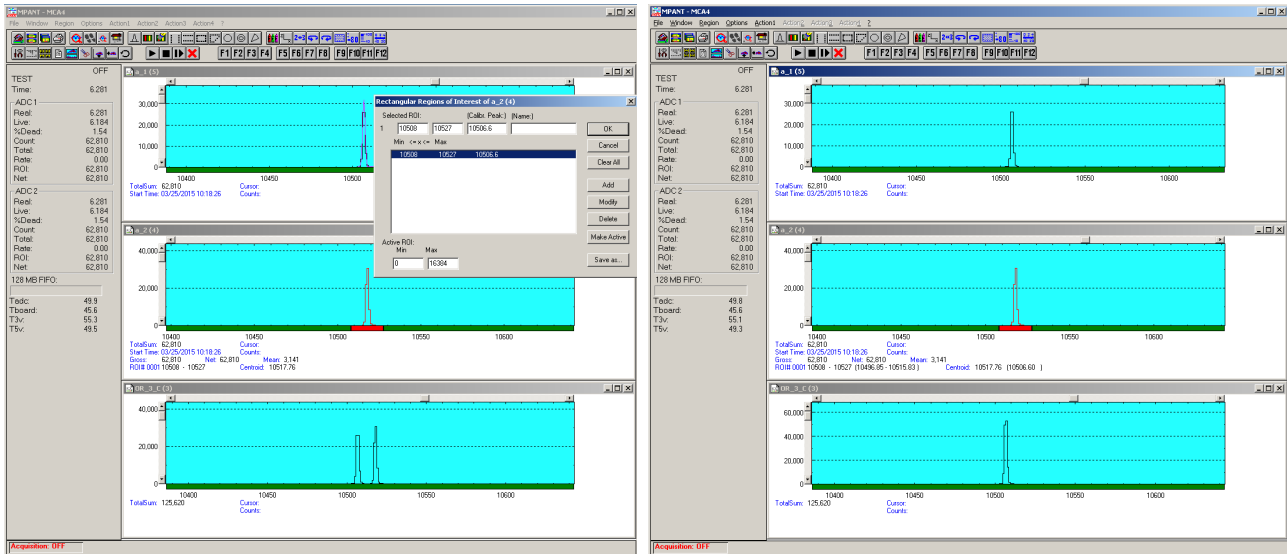
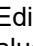


Fig. 4.24: Sum of counts spectrum before and after calibration

Another good method to enter the calibration is to draw in the second spectrum a ROI over the peak by keeping the right mouse button pressed and clicking the "Create ROI" icon , then select "Region..Edit" item from the MPANT menu, select the ROI in the list, enter the fit result from the first spectrum for the value of the Calibration Peak and click "Modify" and "OK". Then choose Region..Auto Calib in the MPANT menu.

The automatic generated spectra name is for example "OR_3_C". The 3 is here an hexadecimal bit pattern indicating that ADC1 (bit 0) and ADC2 (bit 1) contribute to the OR sum. The C at the end of the name means that "Use Calibration" is enabled.

4.1.12 Conditions

If the RTC and Coincidence options are enabled, press Conditions... from the "Dualparameter and Calculated Spectra" dialog to define or edit ROI Conditions in the Conditions dialog.

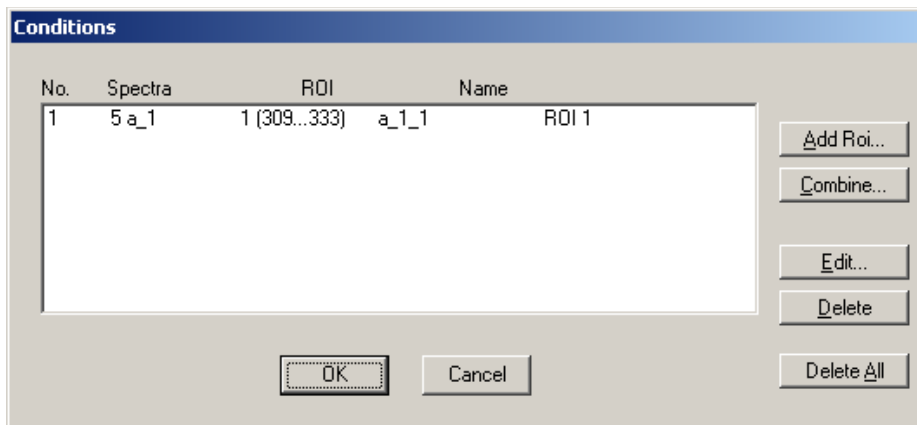


Fig. 4.25: Conditions

To define a new condition, press **Add Roi...** to open the ROI Condition dialog.

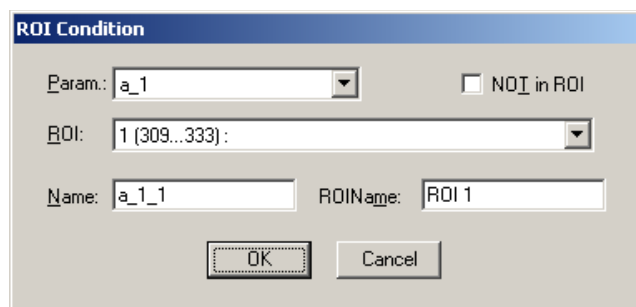


Fig. 4.26: ROI Conditions dialog

Here it is possible to define a Condition as an event inside or outside of an ROI in any parameter spectrum i.e. a calculated parameter or a "copy" of a primary spectrum. There are two drop-down list boxes for the **Parameter** and **ROI**. Of course the ROI must be defined before by using the MPANT program, it is not possible here to add any new ROI. A checkbox labeled **NOT in ROI** allows to define a condition as an event outside the ROI. Using the edit field **Name** it is possible to define a name for the Condition. A default name built by the parameter name and ROI number will be set automatically.

Press **Combine...** to open the **Combine Conditions** dialog.

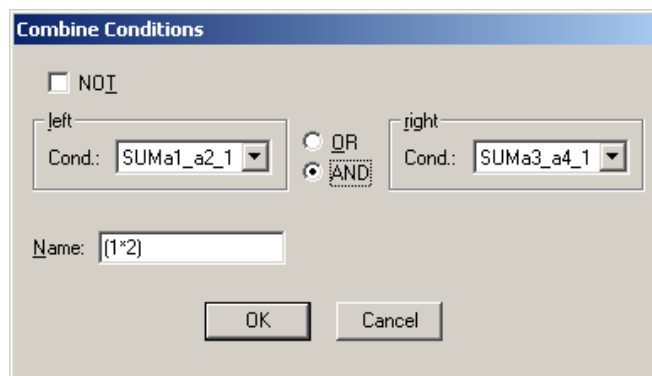


Fig. 4.27: Combine Conditions dialog

Here a Condition can be defined as a combination using the Boolean operators **NOT**, **OR** or **AND** of already defined conditions. The OR will be symbolized in the automatically generated name by a plus sign "+", the AND by an asterisk "*".

If tag bits are enabled in the user licence, it is also possible to define tag bit conditions. The conditions dialog contains then a button **Add Tag...**

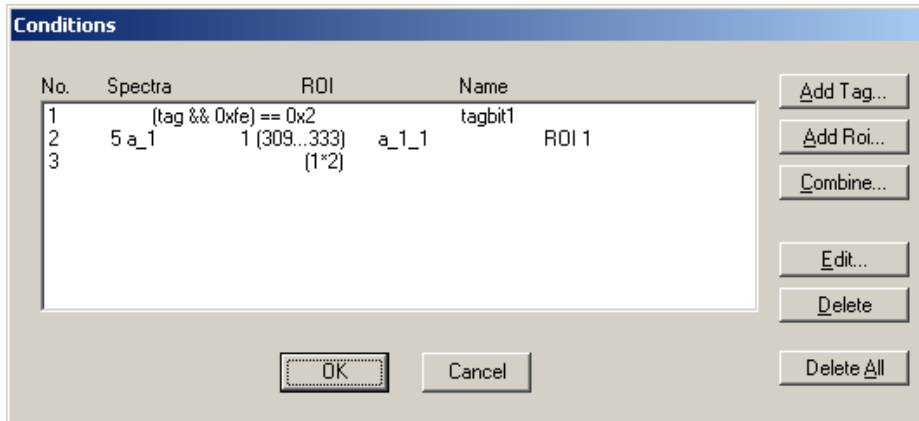


Fig. 4.28: Conditions dialog with tag bits enabled

To define a new tag condition, press **Add Tag...** to open the Tag Condition dialog.

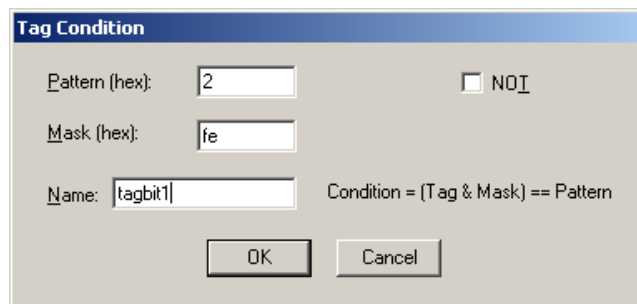


Fig. 4.29: Tag Conditions dialog

Here a Condition can be defined that is true for a special tag bit pattern. Two parameters **Pattern** and **Mask** can be edited in hex format with two hexadecimal digits. The Mask defines which tag bits are used, and the Pattern defines for which tag bit pattern the condition is true (or false if **NOT** is checked). The shown example with Pattern=2 and Mask=fe is a condition that is true if tag bit 1 is ON and tag bits 2 to 7 are all OFF, tag bit 0 does not care.

4.2 File formats

The .mpa format is used to save all spectra in a single file. It starts with an ASCII header containing the settings and then the spectra follow one after the other, each preceded with a header line like

```
[DATA0,16384 ]
```

(This means the first single spectrum with a length of 16384 channels.)

```
[DATA1,16384 ]
```

(This means the second single spectra...)

```
[CDAT0,16384]
```

(This is the first dual parameter or calculated spectra with a length of 16384 channels. The enumeration starts again at zero after all "physical" spectra.)

The format of the data itself can be ASCII, Binary, GANAAS, EMSA or CSV (see below). Individual spectra can also be saved into single files. If "separate header" is checked in the data operations dialog, such spectra data is written into two separate files, one with extension .mp containing configuration data and one containing pure spectra data with an extension indicating the chosen format. The .mp file contains the settings in ASCII format using the control language described in section 4.3

Spectra data files with extension .asc contain in each line one decimal number in ASCII containing the corresponding count value in the histogram.

Binary data files with extension .dat are written with 4 bytes per data value, as usual in the Intel world in reverse order i.e. the least significant byte comes first.

GANAAS data files with extension .spe contain 10 data in ASCII format after a simple header. This format is used by the IAEA gamma spectra evaluation program GANAAS.

EMSA data files are in ASCII with extension .emsa and follow the NIST EMSA spectrum file format.

Another ASCII file format is the x y format with extension .csv. It can be read for example with Excel and contains the channel number and content as two decimal numbers in ASCII per line separated by a TAB character.

List files have the extension .lst and start with a header containing the usual report and configuration data in ASCII as in the .mpa or .mp files. The header ends with a line containing [DATA].

Then follows the data, depending on the format chosen for the data file either in ASCII or binary. In ASCII format one data word is written in hex format per line. In binary format each event is written with 8 bytes, as usual in the Intel world in the reverse order, i.e. the least significant byte comes first. A single event is a 64 bit number (8 bytes), or in scope mode a multiple of such 8 byte numbers.

The lowest two bits 0 and 1 are the channel number 0..3 for ADC1..ADC4.

Bit 2 is a flag for pileup detected.

Bit 3 indicates scope mode. If it is on, more data words follow.

Bit 4 – 47 is a 44 bit event time. If the RTC option is enabled, it is in units of 8 ns, otherwise in units of 1 msec. With tag bits enabled, the 8 most significant event time bits are used for tagging and represent the 8 dig I/O input bits. If pulse width measurement is enabled, the most significant 12 bits are used for the pulse width and represent the pulse width at the threshold level in units of 64 ns.

Bit 48 – 63 are the 16 bit ADC value or in scope mode the length of the waveform data minus one in units of 16 bit words. This size value can be 4095, 8191, 16383 or 32767. The waveform data follow then next.

Here is an example of list file data:

```
[DATA]
adf20000000266f0
ad4f0000000266f2
ae0c0000000266f3
ae820000000266f1
ae850000000573a1
ad500000000573b2
ae0e0000000573a3
```

The software adds in some time intervals and at the end of a list file a small data block containing the counter values.

The format follow the scheme of scope mode data, it can be distinguished by the waveform length that is only 12 16-bit words.

The lowest two bits 0 and 1 are the channel number 0..3 for ADC1..ADC4.

Bit 2 indicates that more data words follow. (This is different from the MCA4, here it was bit 3.)

Bit 4 – 47 are zero.

Bit 48 – 63 are the length of the counter data minus one in units of 16 bit words. This size value is 11 or in hex format b for 12 16-bit words or 3 64-bit counter data. The counter data follow then next, real time, live time and event counter. The counter data are 64 bit values

Here is an example of a counter data block:

```
000b00000000000004
0000000007735943
00000000075aacad
0000000000002710
000b00000000000005
0000000007735943
00000000075aaed3
0000000000002710
000b00000000000006
0000000007735943
00000000075ab725
0000000000002710
000b00000000000007
0000000007735943
00000000075aaf24
0000000000002710
```

4.3 Control Language

A sequence of commands that are stored in a file with extension .CTL can be executed by the MCA4 server program or MPANT with the "Load " or "Run" command. Also the configuration files MCA4A.SET or the header files with extension .MP contain such commands to set the parameters. Each command starts at the beginning of a new line with a typical keyword, the case is ignored. Any other characters in a line may contain a value or a comment.

Following methods are available to execute commands:

- Define one of programmable toolbar icons F1..F12 in the MPANT program with the command and click the icon, see Options..Function keys in the MPANT menu.
- Load the command file using the Load command in the file menu.
- Enable remote mode in the server and send commands via the serial connection. A special DLL is necessary which is part of the optional available MPANT EXTERNAL CONTROL software.
- Open a DDE connection and send the commands via DDE as described in section 4.4. The application name for opening the DDE connection with the standard MCA4 server program MCA4.EXE is MCA4, the topic is MCA4-. Implemented are the DDE Execute to perform any command, and the DDE Request with items RANGE and DATA.
- Send the commands over a TCP/IP net using a remote shell and the optional available MPANT EXTERNAL CONTROL software. It is necessary to have a TCP/IP Winsock installed and that the remote shell daemon program MPA4NET is running. See the readme file on the installation disk.
- Send the commands via the DLL interface from LabVIEW, a Visual Basic program or any other application (software including the complete source code of the DLL and examples optional available).
- From your own Windows application, register a Windows message and then send the command as can be seen in the DLL source code.

The file MCA4A.SET contains a complete list of commands for setting parameters. An example is:

```
[MCA4A] 104           ; Device and serial number
samplemode=10f       ; bit 0..3 chn1..4 0: scope-mode, 1: phamode
                     ; bit 4..7 chn1..4 pile-up reject
                     ; bit 8: Software-PHA
enable_dither=1      ; enable dithering, improves the differential linearity of the ADC chip.
pr_ena=0             ; Preset enable flags, MCS mode:
                     ; bit 0: runtime preset enabled
                     ; bit 1: realtime preset enabled
                     ; bit 2: sweep preset enabled
                     ; bit 3: ROI preset enabled
                     ; bit 5: ROI2 preset enabled
enb_ch=1             ; bit0..3: enable ADC chn1..chn4
thr1=ad20c62         ; LOWORD: threshold 1 HIWORD: zero in 64 bit ADC LSB's.
thr2=b370cc7         ; threshold 2 and zero
thr3=9dc0b6c         ; threshold 3 and zero
thr4=9dc0b6c         ; threshold 4 and zero
width1=531ffe        ; LOWORD: maxwidth1 samples per pulse, HIWORD: minwidth1
width2=531ffe
width3=531ffe
width4=531ffe
digio=0              ; LOWORD: Use of Dig I/O, GO Line:
                     ; bit 0: status dig 0..3
                     ; bit 1: Output digval and increment digval after stop
                     ; bit 2: Invert polarity
                     ; bit 3: Push-Pull output
                     ; bit 4: Start System 1 with Dig 4
                     ; bit 5: Start System 1 with GO
                     ; bit 8: GOWATCH
                     ; bit 9: GO High at Start
```

Software Description

```

; bit 10: GO Low at Stop
; bit 11: Clear before triggered start
digval=0      ; 0..255 DIG I/O Output value
trigger=6    ; bit 1..4 control_trigger for scope mode.
             ; bit 0: TRG_EDGE_DOWN
             ; bit 1: TRG_EDGE_UP
             ; bit 2..4: TRG_SRC, 0=none (softw.), 1..4=CH1..4, 5=extern
             ; bit 5: 0=normal, 1=auto
             ; bit 6: single shot
mcsmode=0    ; bit 16: mcs mode
             ; bit 17: 2mcs mode
             ; bit 18: sequential mode
             ; bit 19: ext. clock
cycles=1     ; cycles for sequential mode
mcspol=10    ; bit 0: start with rising edge
             ; bit 1: stop1 rising edge
             ; bit 2: stop2 rising edge
             ; bit 3: ch adv rising edge
             ; bit 4: FAST NIM
             ; bit 5: 4k7 MCS Inputs
svapol=f     ; for external trigger and SVA:
             ; bit 0: trigger / ch1 with rising edge
             ; bit 1: ch2 rising edge
             ; bit 2: ch3 rising edge
             ; bit 3: ch4 rising edge
             ; bit 4: FAST NIM
             ; bit 5: 1=4k7, 0=50 Ohm
dwellunit=0  ; choice of units in dialog 0=ns, 1=usec, 2=msec, 3=sec
dwelltime=10 ; in units of 10 ns
swpreset=10  ; sweep preset value
timepreset=10.000 ; time preset
unipolar=1   ; unipolar or bipolar (0)
event_cnt_mode=1 ; count only stored events
fifosize=0   ; 0=128 MB, 1=32 MB, 2=8 MB
sca1=ffff   ; LOWORD: upperlevel SCA 1, HIWORD: lowerlevel SCA 1
sca2=ffff
sca3=ffff
sca4=ffff
phafit=0    ; bit0..1 chn1: 0=max, 1= fit, 2 flattop, 3 reserved
             ; bit2..3 chn2
             ; bit4..5 chn3
             ; bit6..7 chn4
             ; LOWORD: fitwidth gaussian, HIWORD, pulsewidth flattop
fitwidth1=29
fitwidth2=29
fitwidth3=29
fitwidth4=29
memsegs=aa  ; Bit0..1: SEGMENTS_CH1, 0..3 = 1x32k/2x16k/4x8k/8x4k segments
             ; Bit2..3: SEGMENTS_CH2
             ; Bit4..5: SEGMENTS_CH3
             ; Bit6..7: SEGMENTS_CH4
ctm8=125    ; coincidence time in units of 8 ns
mwidth1=0   ; flattop: range for Mean calc LOWORD upper, HIWORD lower
mwidth2=0   ; LOWORD: Bit0..13 meanrange, Bit(15,14) Method 0=left, 1=cft, 2=right
mwidth3=0   ; HIWORD: roleft
mwidth4=0
gate=0      ; bit 0..3: GATE_ADC1..4, 1=Gate must be active for pulse processing
             ; bit 4..7: GATE_POL1..4, 1=active high
             ; bit 8..11: REJECT_ADC1..4, 1=use GATE as reject, e.g. pile-up detected,
             ; cor. GATE_ADCx must be 0
             ; bit 12..15: SVA mode
ctlsca=0    ; bit 0..3: SCA_CH1, 0=SCA, 1=Over-Threshold, 2=GATE Input, 3=MCS Input
             ; bit 4..7: SCA_CH2,
```

Software Description

```

; bit 8..11: SCA_CH3,
; bit 12..15 SCA_CH4
autoinc=0      ; 1=Enable Auto increment of filename
savedata=0    ; bit 0: auto save after stop
              ; bit 1: write listfile
              ; bit 2: listfile only, no evaluation
              ; bit 5: drop zero events

mpafmt=asc    ; format used in mpa datafiles
              ; (dat=binary, asc=ASCII, csv=CSV, spe=GANAAAS, emsa=EMSA)

sephead=0    ; seperate Header
fmt=dat      ; data format used in separated spectra (extension MP)
              ; (dat=binary, asc=ASCII, csv=CSV, spe=GANAAAS, emsa=EMSA)

mpacaluse=0  ; use calibration for shifted summing .mpa files according to calibration
smoothpts=5  ; Number of points to average for a smooth operation
wndwidth=160 ; Sets width of server window
wndheight=24 ; Sets height of server window
sysdef=0     ; Defines which module belongs to a system, relevant only for more
              ; modules

[MCS1]       ; MCS1 parameters
active=1     ; active
range=4096   ; spectrum length
[MCS2]       ; MCS2 parameters
active=0     ; active
range=4096   ; spectrum length (same as MCS1 range)
[CHN1]       ; Scopemode Channel 1 parameters
active=0     ; active
range=4096   ; waveform length
[CHN2]
active=0
range=4096
[CHN3]
active=0
range=4096
[CHN4]
active=0
range=4096
[ADC1]       ; ADC1 parameters
active=1     ; active
range=16384  ; spectrum length
prena=0     ; bit 0: realtime preset enabled
              ; bit 1: livetime preset enabled
              ; bit 2: Event preset enabled
              ; bit 3: ROI preset enabled
              ; lower ROI limit for ROI preset
roimin=0    ; upper limit: roimin <= channel < roimax
roimax=16384
evpreset=1  ; ROI preset or event count preset
rtpreset=1000.000 ; real time preset (sec)
ltpreset=1000.000 ; live time preset (sec)
bitshift=2  ; bit shift between ADC LSB's and spectrum resolution
caloff=0.000000 ; calibration parameter: offset
calfact=1.000000 ; calibration parameter: factor
calfact2=0
calfact3=0
calunit=keV ; calibration unit
caluse=0    ; bit 0=1: use calibration, higher bits: calibration formula
[ADC2]     ; the following section concerns parameters of ADC2
...

```

The following commands perform actions and therefore usually are not included in the MCA4A.SET file:

Software Description

start	; Clears the data and starts a new acquisition of system 1. ; Further execution of the .CTL file is suspended until any ; acquisition stops due to a preset.
halt	; Stops acquisition of system 1 if one is running.
cont	; Continues acquisition of system 1. If a time preset ; is already reached, the time preset is prolonged ; by the value which was valid when the "start" command ; was executed. Further execution of the .CTL file ; is suspended (see start).
erase	; Clears all spectra of system 1.
savecnf	; Writes the settings into MCS6A.SET
savempa	; Saves all configuration and spectra data. ; An existing file is overwritten.
loadmpa	; Loads mpa data; the filename ; must be specified before with a command mpaname=...
addmpa	; Adds mpa data to actual spectra; the filename ; must be specified before with a command mpaname=...
submpa	; Subtracts mpa data from actual spectra; the filename ; must be specified before with a command mpaname=...
MC_A	; Sets actual input channel to MC_A (STOP1) for the rest of ; the control file.
MC_B	; Sets actual multichannel analyzer to MC_B (STOP2)
savedat	; Saves data of actual channel as separated ; spectrum (extension .MP) An existing file ; is overwritten.
loaddat	; Loads data of actual channel, the filename ; must be specified before with a command datname=...
adddat	; Adds data into actual spectra; the filename ; must be specified before with a command datname=...
subdat	; Subtracts data from actual spectra channel; the filename ; must be specified before with a command datname=...
smooth	; Smooths the data in actual spectra
erase	; Clears the data of system 1.
exit	; Exits the server (and MPANT) programs
alert Message	; Displays a Messagebox containing Message and an OK ; button that must be pressed before execution can continue.
waitinfo 5000 Message	; Displays a Messagebox containing Message, an OK ; and an END button. After the specified time (5000 msec) ; the Messagebox vanishes and execution continues. OK ; continues immediately, END escapes execution.
beep *	; Makes a beep. The character "*" may be replaced with '?', '!' or ; left empty. The corresponding sound is defined in the WIN.INI ; file in the [sounds] section.
delay 4000	; Waits specified time (4000 msec = 4 sec).
pulse 100	; Outputs a pulse of 100 ms duration at dig 3.
waitpin 4000	; Waits 4000 ms for going the level at dig 7 going low. ; After a timeout a Messagebox warns and waits for pressing OK ; Can be used for connecting a sample changer.
waitpin3 4000	; waits for going the level at pin 3 going low. All Dig I/O ports 0..7 ; can be used this way.

Software Description

loadcnf configfile	; Loads a configuration stored in a config file like mca4a.set.
run controlfile	; Runs a sequence of commands stored in controlfile. This ; command cannot be nested, i.e. it is not possible to execute ; a run command from the controlfile called with a run command, ; only if it was called with a load command.
onstart command	; The command is executed always after a start action when the ; acquisition is already running. The command can be any valid ; command, also 'run controlfile' is possible.
onstart off	; Switches off the 'onstart' feature. Also a manual Stop command ; switches it off.
onstop command	; The command is executed always after a stop caused by a ; preset reached or trigger. This can be used to program measure ; cycles. For example the command 'onstop start' makes a ; loop of this kind.
onstop off	; Switches off the 'onstop' feature. Also a manual Stop command ; switches it off.
oncycle command	; The command is executed in sequential mode after each cycle. ; can be used together with 'waitpin' command to synchronize with a ; beam- or target mover.
oncycle off	; Switches off the 'oncycle' feature.
lastrun=5	; Defines the file count for the last run in a measure cycle. After a ; file with this count or greater was saved with autoinc on, instead ; of the 'onstop command' the 'onlast command' is executed.
numruns=5	; Defines the file count for the last run in a measure cycle. The ; last count is the present one plus the numruns number. After a ; file with this count was saved with autoinc on, instead of the ; 'onstop command' the 'onlast command' is executed.
onlast command	; The command is executed after a stop caused by a preset ; reached or trigger instead of the 'onstop command', when the ; last file count is reached with autoinc on. This can be used to ; finish programmed measure cycles.
onlast off	; Switches off the 'onlast' feature. Also a manual Stop command ; switches it off.
exec program	; Executes a Windows program or .PIF file. ; Example: exec notepad test.ctf ; opens the notepad editor and loads test.ctf.
fitrois	; Makes a single peak Gaussian fit for all ROIs in the active ; Display of MPANT and dumps the result into a logfile. This is ; performed by the MPANT program and therefore can be ; made only if this application is running.
fitrois MC_A	; Similar to the fitroi command, but using the argument allows to ; specify which spectrum should be evaluated independently of ; which child window is activated in MPANT. MC_A means STOP1
autocal	; Makes a single peak Gaussian fit for all ROIs in the active ; Display of MPANT for which a peak value was entered in the ; MPANT Region Edit dialog and uses the results for a ; calibration. This is performed by the MPANT program and ; therefore can be made only if this application is running.
autocal SP_1	; Similar to the autocal command, but using the argument allows ; to specify which spectrum should be evaluated independently of ; which child window is activated in MPANT
deleteallrois	; deletes all ROIs in the active MPANT window.

Software Description

deleteallrois MC_# ; allows to delete all ROIs of a specified spectra (1 for # means
; A1..).

The following commands make sense only when using the serial line or TCP/IP control or DLL control interface:

MC_A? ; Sends the status of MC_A via the serial port and make it actual.

...

MC_H? ; Sends the status of MC_H via the serial port and
; make it actual.

? ; Send the status of the actual channel

sendfile filename ; Sends the ASCII file with name 'filename' via the serial line.

The execution of a control file can be finished from the Server or MPANT with any Halt command.

4.4 Controlling the MCA4 Windows Server via DDE

The MCA4 program can be a server for DDE (Dynamic Data Exchange). Many Windows software packages can use the DDE standard protocols to communicate with other Windows programs, for example GRAMS, FAMOS or LabVIEW. In the following the DDE capabilities of the MCA4 program are described together with a demo VI ("Virtual Instrument") for LabVIEW. It is not recommended to use the DDE protocol for LabVIEW, as also a DLL interface is available that is much faster. The following should be seen as a general description of the DDE conversation capabilities of the MCA4 server program.

4.4.1 Open Conversation

application: MCA4

topic: MCA4-

Any application that wants to be a client of a DDE server, must open the conversation first by specifying an application and a topic name. The application name is MCA4 and the topic is MCA4-.

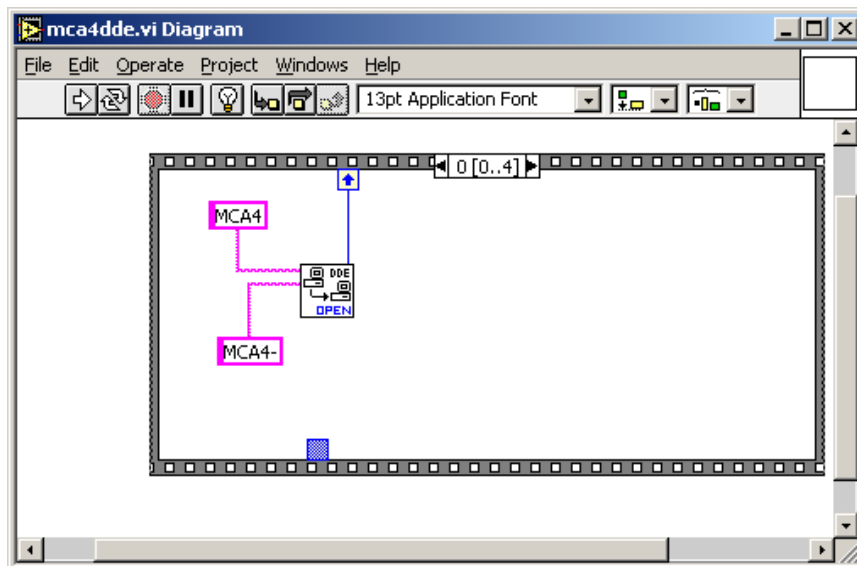


Fig. 4.30: Opening the DDE conversation with the MCA4 server in LabVIEW

4.4.2 DDE Execute

The DDE Execute command can be used to perform any action of the MCA4 program. Any of the Control command lines described in section 4.3 can be used. For example a sequence of control commands saved in a file TEST.CTL can be executed by specifying the command

```
RUN TEST.CTL
```

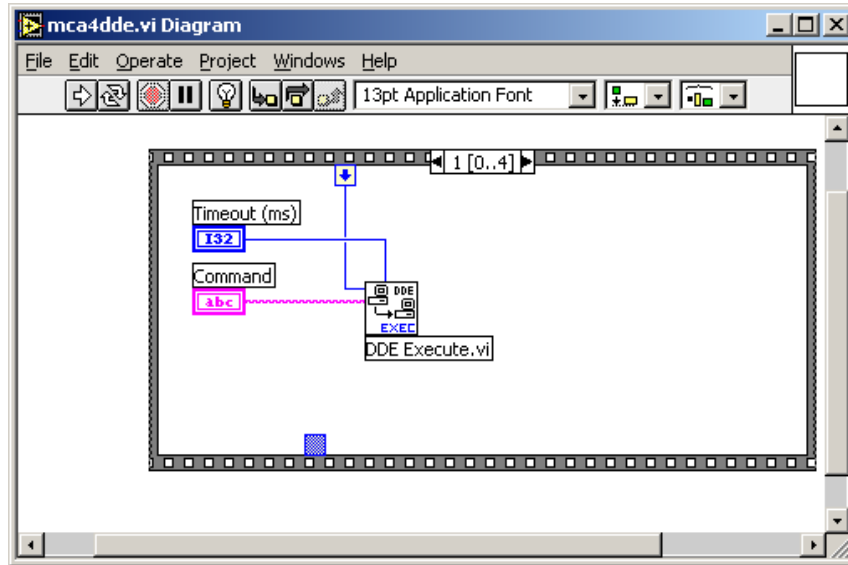


Fig. 4.31: Executing a MCA4 command from a LabVIEW application

The MCA4 program then executes the command and, after finishing, it sends an Acknowledge message to the DDE client. This can be used to synchronize the actions in both applications. DDE Request

The DDE Request is a message exchange to obtain the value of a specified item. Only two items are defined for DDE request up to now: RANGE and DATA. The value is obtained as an ASCII string, i.e. it must be converted by the client to get the numbers. All other parameters concerning the setup can be obtained by the client application by reading and evaluating the configuration file.

RANGE

The RANGE item can be used to obtain the total number of data.

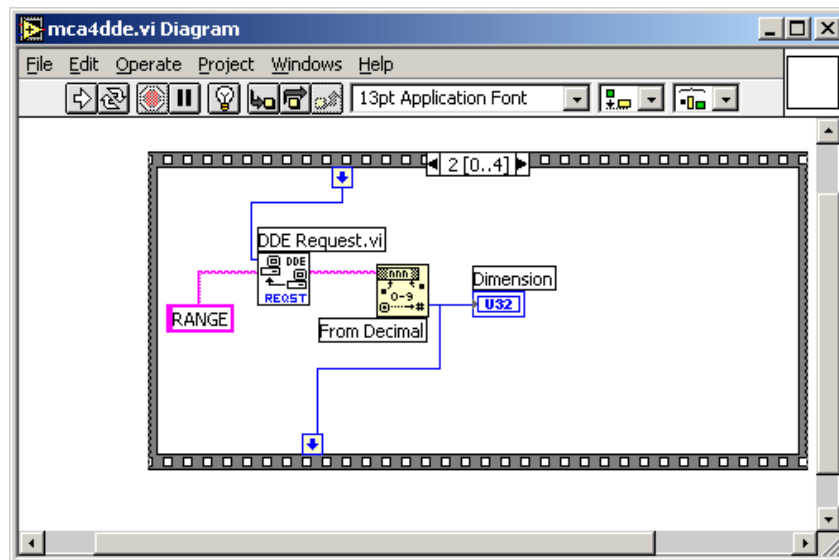


Fig. 4.32: Getting the total number of spectrum data with LabVIEW

DATA

With the DATA item the data is obtained. The value of this item is a multiline string that contains in each line a decimal number as an ASCII string.

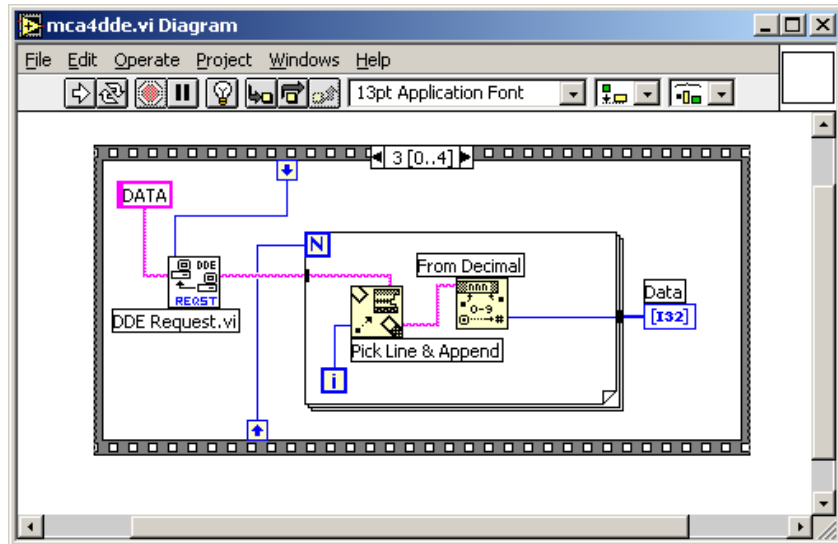


Fig. 4.33: Getting the data with LabVIEW

4.4.3 Close Conversation

After finishing the DDE communication with the server program, it must be closed.

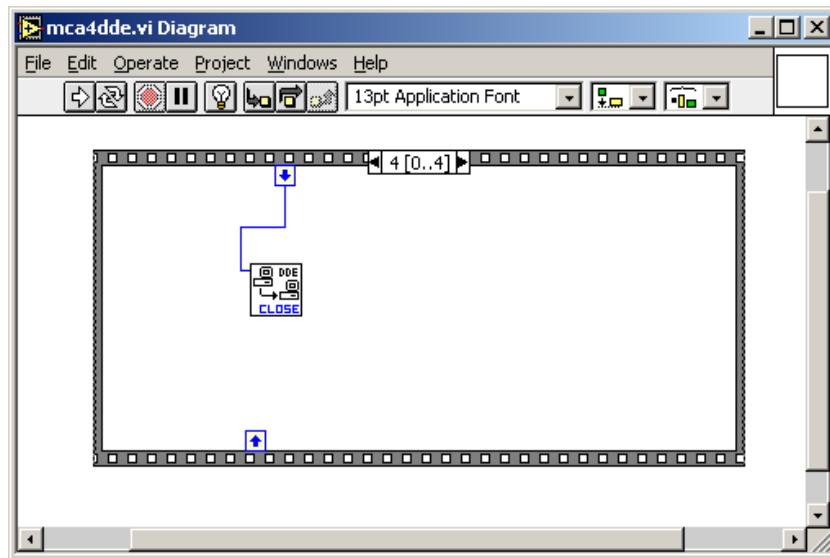


Fig. 4.34: Closing the DDE communication in LabVIEW

The following figure shows the “Panel” of the described VI for LabVIEW.

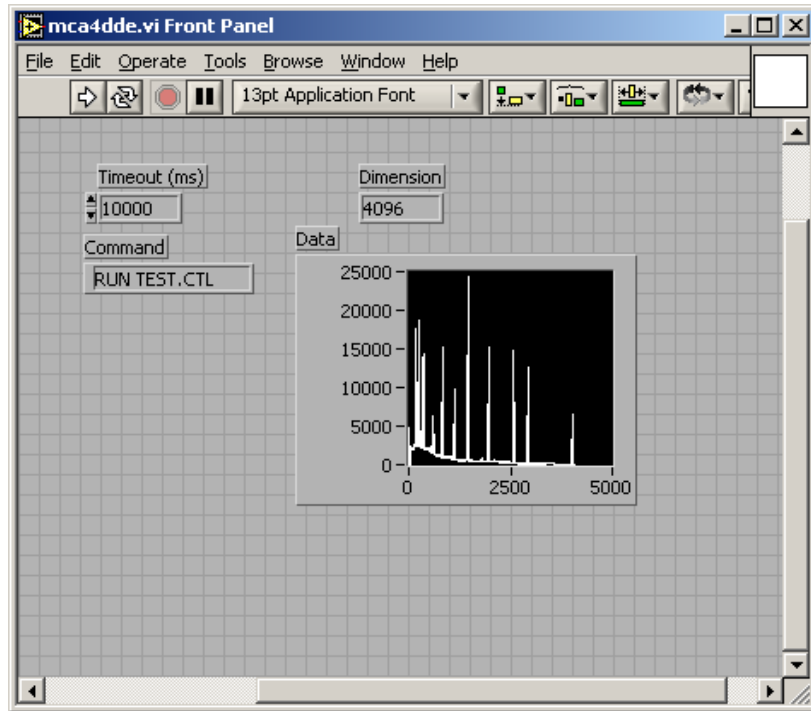


Fig. 4.35: Control Panel of the demo VI for LabVIEW

4.5 Controlling the MCA4 Windows Server via DLL

The MCA4 server program provides access to all functions, parameters and data via a DLL ("dynamic link library"). So the server can be completely controlled by the MPANT software that provides all necessary graphic displays.

NOTE:

The complete documented source code of the DLL including a library with fundamental VI's and an example VI for LabVIEW and an example Visual Basic and C program is available as an option.

The DLL DMCA4.DLL is an interface providing functions to communicate with the server program. Most of these functions send messages to the server as you do it when operating the server program by sending Windows messages via mouse clicks. Please do not expect any functions in this DLL for controlling the hardware directly. All software using the DLL interface requires that the server program is running. The DLL was mainly developed as an interface between the server program and MPANT, not as a nice developing tool for customers. But by looking at the programming examples and some hints in the DLL documentation it should be easy to develop own programs that are able to control the server like MPANT does.

The server program has a built-in command interpreter. The syntax of these commands is described in chapter 4.3, and in the MPANT on-line help (look for: "How to use the command language.."). It is recommended to send commands like "range=16384" to the server via the RunCmd DLL function, i.e. RunCmd(0, "range=16384"); if you want to set parameters like a spectra length. The alternative method is to store all settings parameters into the DLL by calling the DLL function StoreSettingData(setting, 0); and then calling NewSetting(0); to send a message to the server to read new settings from the DLL. This method will not work for changing a spectra length to avoid the problem of any undefined memory pointers. The range parameter should be changed by the server program only (or by sending a "range=.." command). The recommended usage of the DLL is reading parameters like Status, Settings, Strings, Cnt numbers, ROI boundaries using the corresponding DLL functions, but for any actions or setting any parameters the command interpreter should be used.

It is important that the DLL is loaded first by the Server program and that it is loaded from the same path by all programs using it. Otherwise it does not work to access the shared memory. The dmca4.dll is installed by the installation program into the standard Windows folder for 32-bit or 64-bit DLLs, respectively. Please make sure that there is nowhere else any file dmca4.dll.

5 MPANT Software

The window of the MPANT program is shown here. It enables the full control of the MCA4 via the server program to perform measurements and save data, and shows the data on-line in several windows.

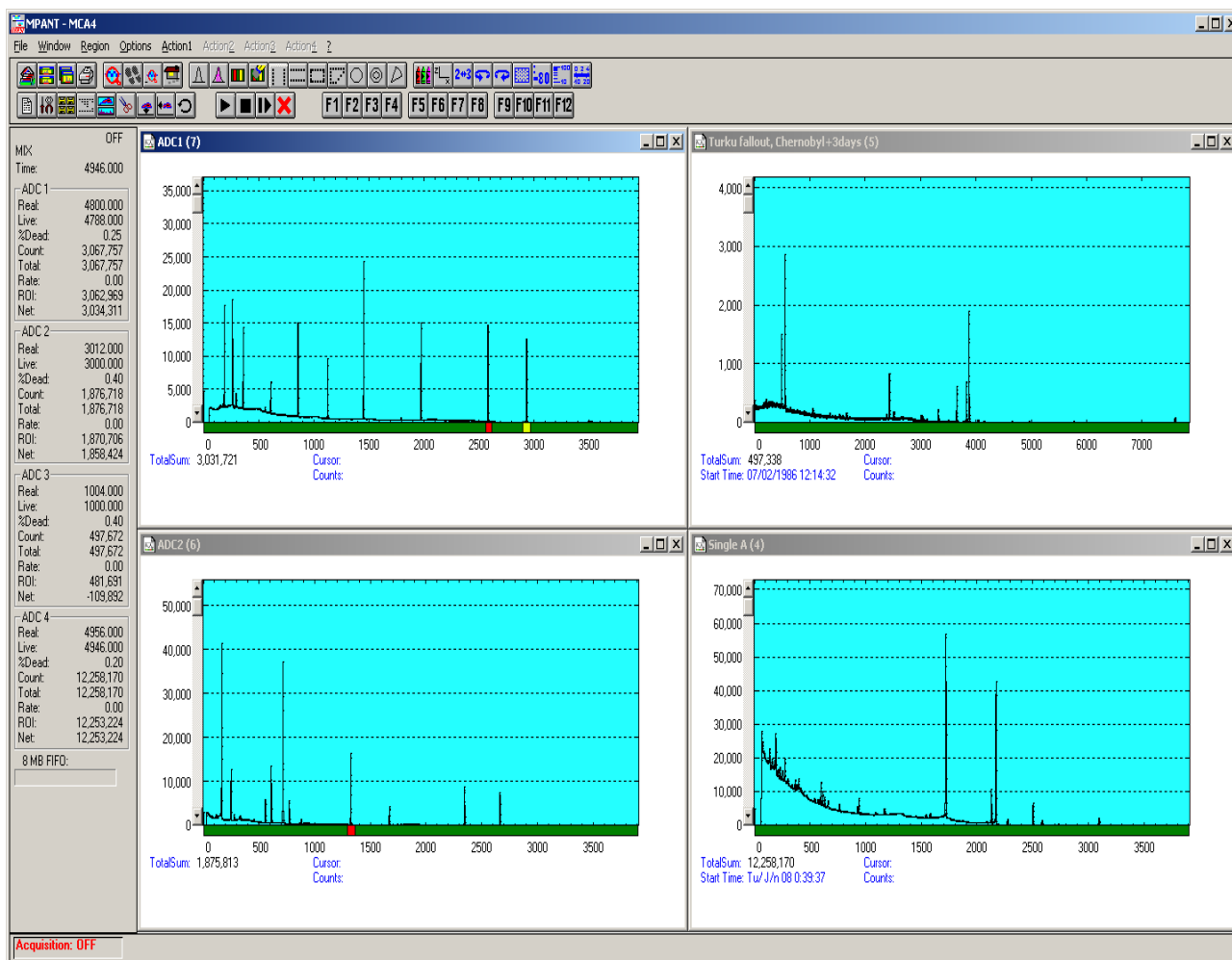


Fig. 5.1: MPANT main window

The server program MCA4.EXE automatically starts MPANT. If you try to start MPANT before the server is started, a message box warns that you should start the server first. If you want to use the server without the MPANT display program, you may rename the file "mpant.exe" into "mym pant.exe". It will then not be found by the server and not started. MPANT can be closed at any time without affecting the data acquisition and can be restarted out of the server menu or by a shortcut on the desktop.

MPANT has viewing capabilities for single- and two-dimensional spectra. For ROI's the ROIsum, Net sum and the mean counts per channel (i.e. ROIsum divided by number of channels) and the centroid is displayed (<x> for single spectra, <x> and <y> for two dimensional). A single spectrum can be converted into a two-dimensional one by specifying the x dimension in the display option dialog. It is possible to drag a rectangle and zoom into this rectangle. Rectangular ROI's can be set and the ROIsum and Net ROIsum is displayed. The Net Sum is calculated the same way like in the single view, by subtracting a linear interpolated background from the both outmost channels in x-direction. This Net sums are then summed up in y-direction. The ROI editing dialog is changed into a Rectangular Editing dialog for MAP and ISO displays. The Cursor can be moved in x and y direction using the mouse and the arrow keys, in ISO display only using the arrow keys.

A status window at the left side gives all information about the status of the MCA4. Below the different counting rates a progress bar shows the filling of the FIFO.

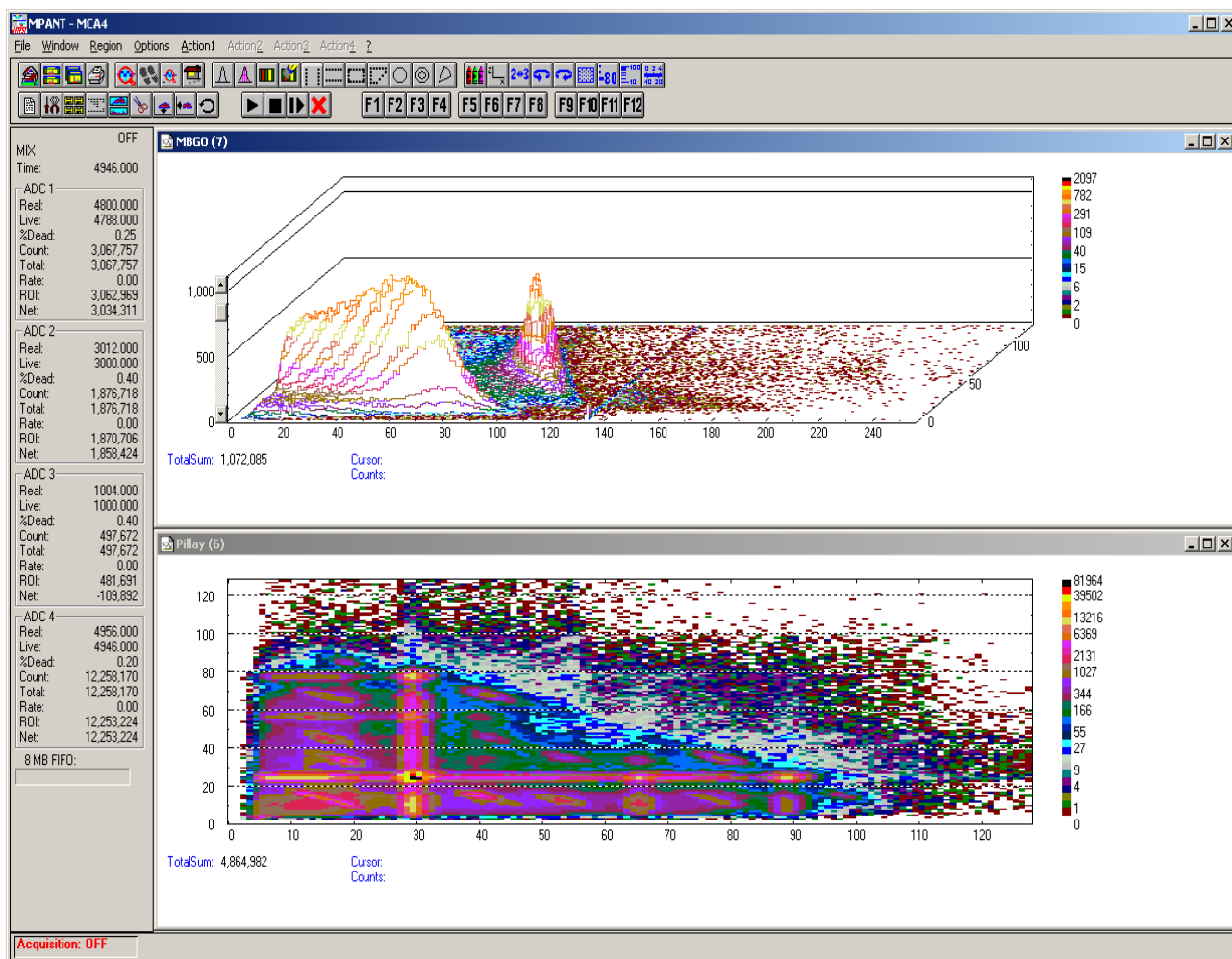


Fig. 5.2: MPANT Map and Isometric display

A toolbar provides fast access to many used functions in the menu. A status bar at the bottom gives help about the meaning of the toolbar icons. A cursor appears when clicking the left mouse button inside the graphics area. The cursor can be moved using the arrow keys. To get rid of the cursor, make a double click with the right mouse button outside the graphics area. To define a region, press the right mouse button, and while keeping the button pressed, drag a rectangle. In zoomed state a scroll bar appears that allows scrolling through the spectrum.

In the following the several menu functions are described together with the corresponding toolbar icons.

5.1 File Menu

Load..., Add..., Save Display As..., Save MPA As...

These menu items provide the usual functions for loading and saving data, either a complete data set when loading or saving a .MPA data file, or separated spectra when loading or adding a .MP data file. The 'Load' function creates then a new spectrum for compare, whereas the 'Add' functions adds the data to the spectrum selected by the active window. For saving a selected spectrum into a seperated file, click on the respective display window and select "Save Display As.." from the MPANT File menu. If you want to have the data file seperated from the header containing the configuration, first open the data operations dialog using the toolbar icon showing the 4 yellow diskettes, select "Selected Spectra" and check "Seperate Header" as shown in Fig. 4.3. When saving data, you have the choice between several formats like binary (.DAT), ASCII (.ASC), GANAAS (.SPE), EMSA (.emsa) and CSV (.CSV) format by selecting it from "Save type As...". To read a CSV file with Excel, load it as a text file and configure TAB characters as separators.

When you load data, select a header file (extension .MP). This file contains the information about the size and format of the data file, which is then automatically read. With "Add" the data is added to the present data.

The data read from a file is shifted according to the calibration, if it is available and 'Use Calib' is checked in the MCA4 Data Operations dialog, see Fig. 4.3. This dialog can also be used to load spectrum data from a file into an existing display window.

New Display:



With the Open New menu item or the corresponding icon a new Display window can be created and shown as the active window. In the “Open New Display” dialog box the spectrum for the new display can be selected.

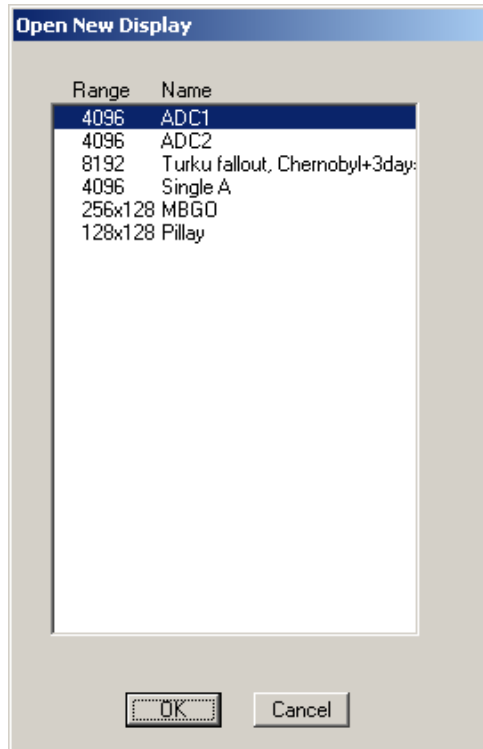


Fig. 5.3: File New Display dialog

Open All

By selecting the Open All menu item, all available Displays are shown. The windows of the last opened Display become active.

Close All

By selecting the Close All menu item, all Display windows are closed.

Compare...

The Compare... menu item opens a dialog that enables to show two spectra with different colors in the same window. The active display window can be checked as a primary window and there is a list of other spectrum windows that can be inserted as secondary windows.

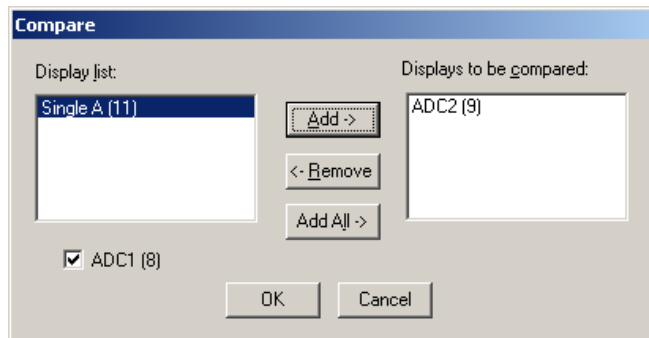


Fig. 5.4: Compare dialog

Using the features of the software it is possible to create a shifted spectrum that has a characteristic peak at the same place as the primary spectrum, see 4.1.11 . This way a view as shown in Fig. 5.5 can be created.

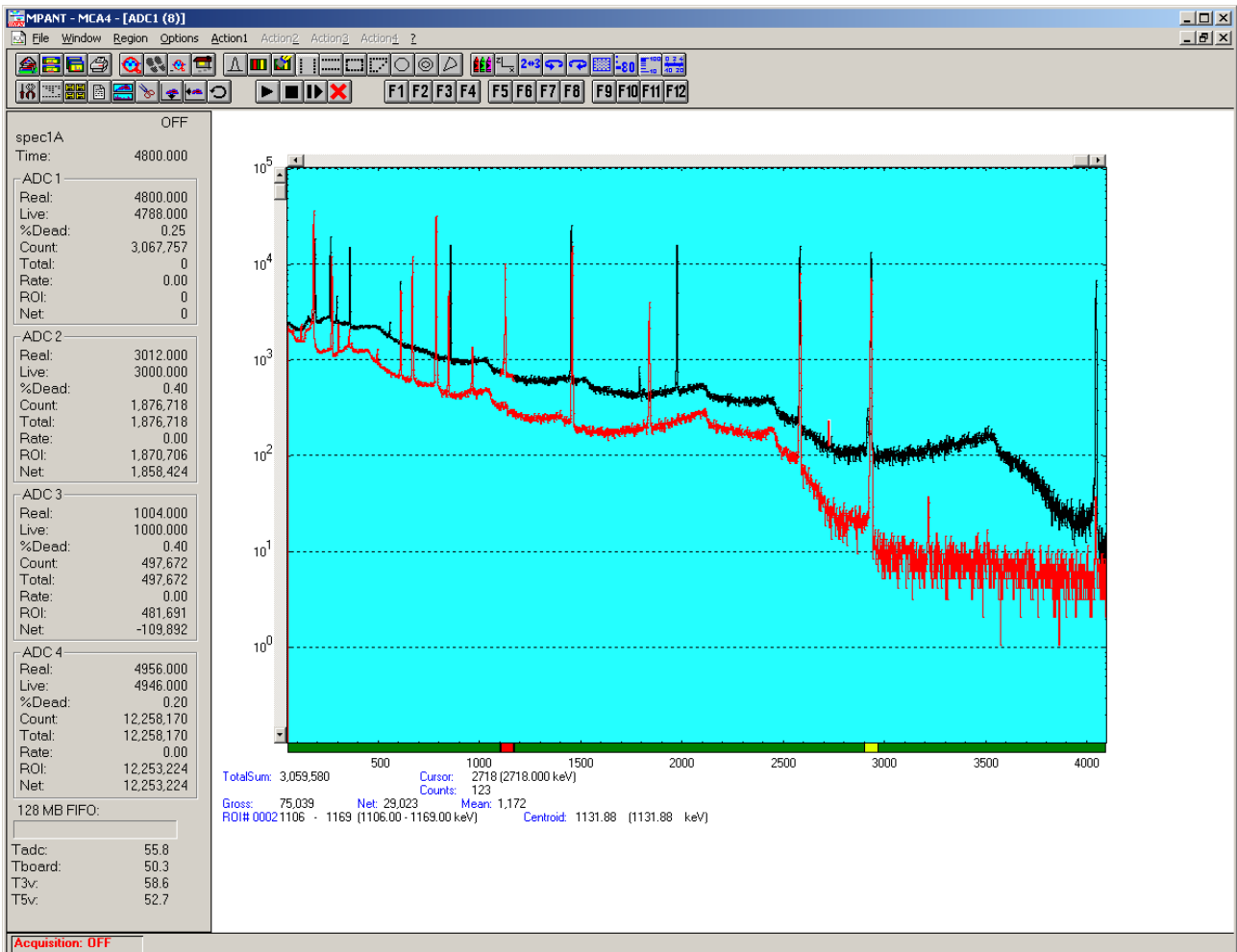


Fig. 5.5: Compare two spectra

Print...



The Print menu item opens the print dialog. It allows to arrange several pictures on a page into zones. The number of zones in vertical and horizontal direction can be specified. The Color can be black/white, RGB (colored) or Gray scale. RGB is recommended also for black laser printers. Some info lines containing date, file name and title can be added. For each page a temporary file PRINT1.WMF, PRINT2.WMF... will be created. This file is in Windows Meta file format and can be exported into some other Windows applications.

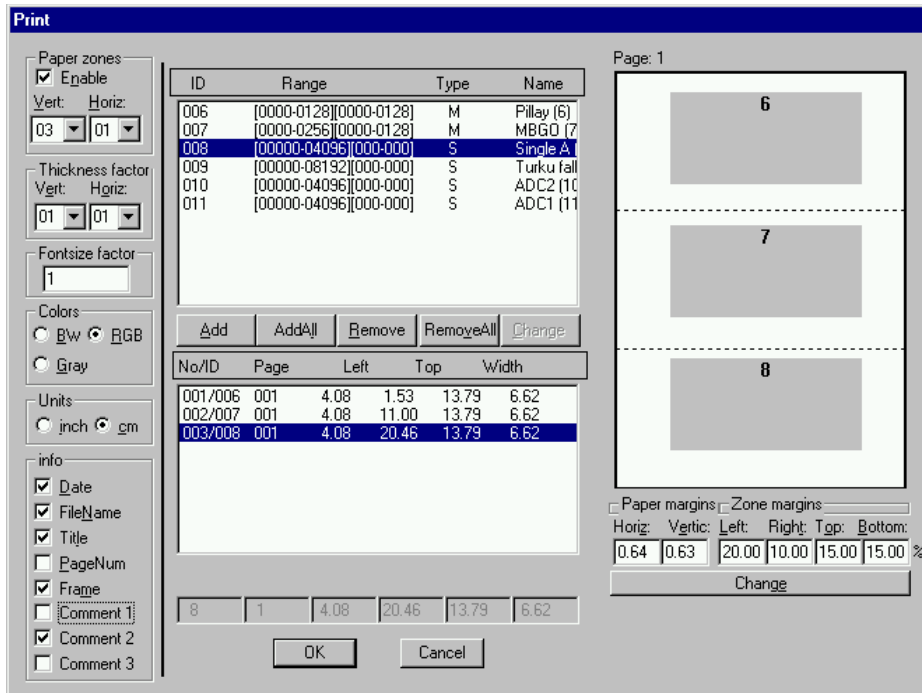


Fig. 5.6: Print dialog

Setup Printer...

The Setup Printer menu item enables to configure the printer.

Exit

The Exit menu item exits the MPANT program. The server and any measurement data remain active, MPANT can be easy restarted out of the file menu of the server program.

5.2 Window Menu

The Window menu enables arranging the Display windows.

Tile



With the Tile menu item or by clicking the corresponding icon, all opened and displayed MPANT Display windows are arranged over the full MPANT client area trying to allocate the same size for each window.

Cascade



The Cascade menu item or respective icon arranges all windows in a cascade display.

Open All


By selecting the Open All menu item, all Display windows are opened.

Window list

At the end of the Window menu, all created Display windows are listed with their names, the current active window is checked. By selecting any of the names, this window becomes the active window and is displayed in front of all the others.

5.3 Region Menu

The Region menu contains commands for Regions and ROI's (Regions of Interest). A Region can be defined by marking it in a display, with the mouse using the right mouse button and dragging a rectangle over the area one is interested in. A ROI, i.e. an already defined region in a single spectrum can be shown zoomed by double-clicking with the left mouse button on the corresponding colored area in the bar at the bottom of the spectra display. A single mouse click with the left button on the corresponding colored area makes this to the selected ROI and lets the counts contained in this ROI be displayed in the information lines of the respective window. Advance the selected ROI to the next or previous one using the "+" or "-" key from the numeric keypad.

Polygonal and curved ROI's in MAP displays are implemented in addition to the rectangular ROI's. The boundaries can be polygonal, circular, annular or pie. Make a choice of the shape using the Menu Region...Shape or the corresponding toolbar icon. Define the ROI using the right mouse button, the TAB and the END key (please read the help in the status bar) and save the ROI using the toolbar icon "Create new ROI" . Advance the selected ROI to the next or previous one using the "+" or "-" key from the numeric keypad. The ROI boundary is displayed by lines or circles and the ROI area can in addition be displayed in inverted color by enabling the checkbox "Invert ROI Pixel" in the MAP Display options dialog. The sum and net sum of counts within the selected ROI is displayed. The net sum is calculated by subtracting a background integral by summing line integrals $S(y)$ for each y coordinate using the mean value of the counts at the lowest and highest x coordinate of the ROI.

Zoom



The Zoom item or respective icon enlarges a Region to the maximum Spectrum Display size.

Back



The Back menu item or clicking the corresponding icon restores the last zoom view. Each time a Back command is clicked the view is stepped back one step.

Zoom Out



The Zoom Out menu item or clicking the corresponding icon enlarges the actual zoom view by a factor 2, if possible. You can also use a mouse wheel for zooming in at out.

Home



Clicking the Home menu item or the corresponding icon restores a Display to the basic configuration showing the complete spectrum.

Shape

Selecting the Shape menu item opens a sub menu with the items Rectangle, X-Slice Y-Slice and Rectangle, Polygon, Circle, Annular and Pie to choose the ROI shape.

Rectangle

Sets the region shape to a rectangle with arbitrary dimensions. To enter the rectangular region, press the right mouse button, drag a rectangle, and release the button to define the region.

X-Slice

Sets the Region shape to the rectangle with maximal height.

Y-Slice

Sets the Region shape to the rectangle with maximal width.

Polygon

Sets the Region shape to polygonal. To enter a polygonal Region using the mouse, press the right button for each point, the TAB key for a next polygon and the END key to close it.

Circle

Sets the Region shape to circular. To enter a circular Region using the mouse, click the right button at the center, move the mouse and press the END key to finish.

Annular

Sets the Region shape to annular. To enter an annular Region using the mouse, click the right button at the center and first circle and press the END key at the second circle.

Pie

Sets the Region shape to pie. To enter a pie Region using the mouse, click the right button at the center and first radial end-point, press the END key at the second radial line end-point.

Create

The Create menu item creates a new ROI from the current marked Region.

Delete

By selecting the Delete menu item or the respective icon, the current active ROI is deleted and the previously defined ROI is activated.

Edit...

With the Edit item, a dialog box is opened which enables editing the ROI list, i.e. create a new or delete, change and activate an existing ROI. Also the peak values for an automatic calibration can be entered here. A ROI can be edited and added to the list. It can also be made to the "Active ROI", that is the special ROI that is used by the server program to calculate the events within this ROI and look for an event preset. The ROI list can be cleared and it can be written into a file with extension .CTL, which can be directly loaded into the server to restore the ROI list.

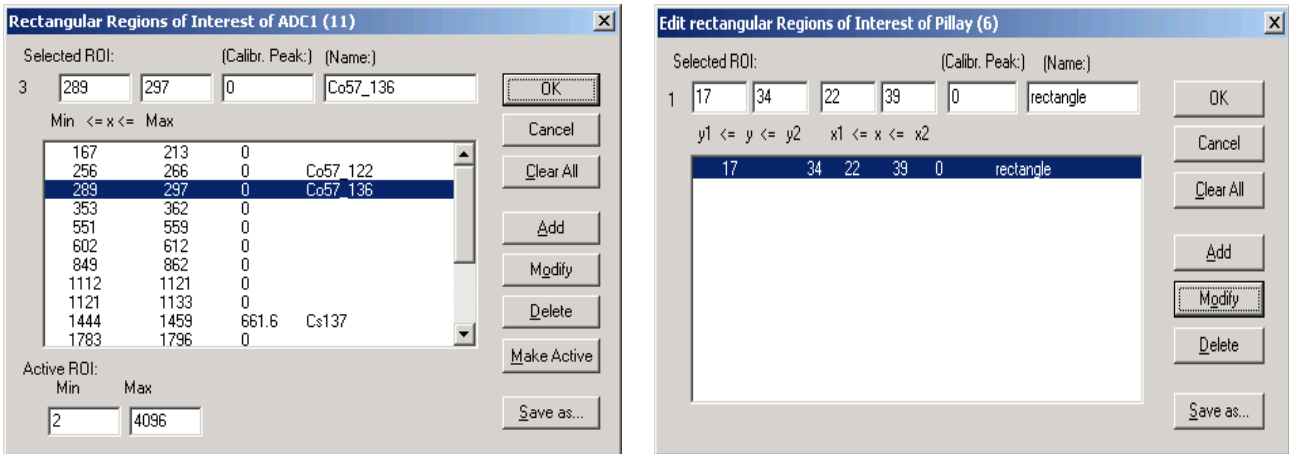


Fig. 5.7: Slice and rectangular ROI Editing dialog

Clicking on it in the ROI list can change the selected ROI. In the MPANT spectrum display the total and net sum of the selected ROI is displayed.

ROI names are implemented. The name can be entered in the ROI editing dialog. Press "Modify" to insert a new name from the edit field of the selected ROI into the list.

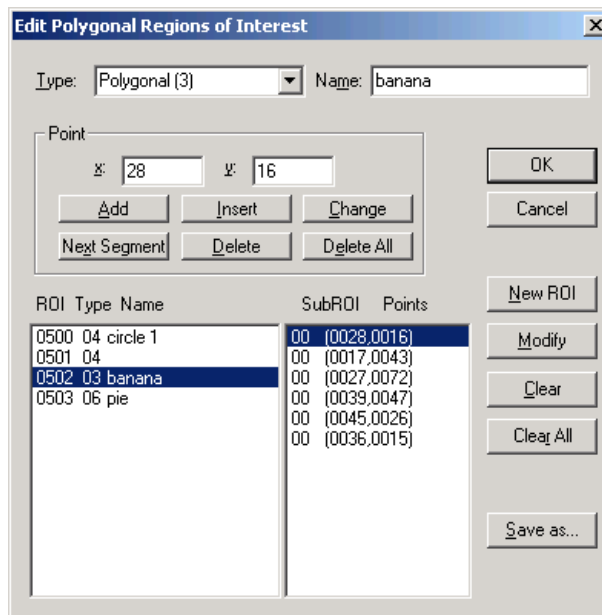


Fig. 5.8: Polygonal ROI Editing dialog

For polygonal, circular, annular and pie ROI's there is a special dialog to select and edit all the parameters. The ROI definition parameters can be saved into a .CTL file and reloaded using the File menu like any control file. Note that the selected ROI can be advanced in any activated spectra display by pressing the "+" and "-" key from the numeric keypad.

X-projection, Y-Projection



Projections of ROI's in dual parameter spectra to the x- and y-axes are implemented. There are menu items "X-Projection" and "Y-Projection" in the Region menu. They create new single spectra displays showing the respective projection of the selected ROI in a dual parameter spectrum in the active window.

Fit...



By selecting the Fit... menu item or the respective icon, a single Gaussian peak fit with linear background is performed for the currently marked region. The fitted curve is displayed and a dialog box shows the results:

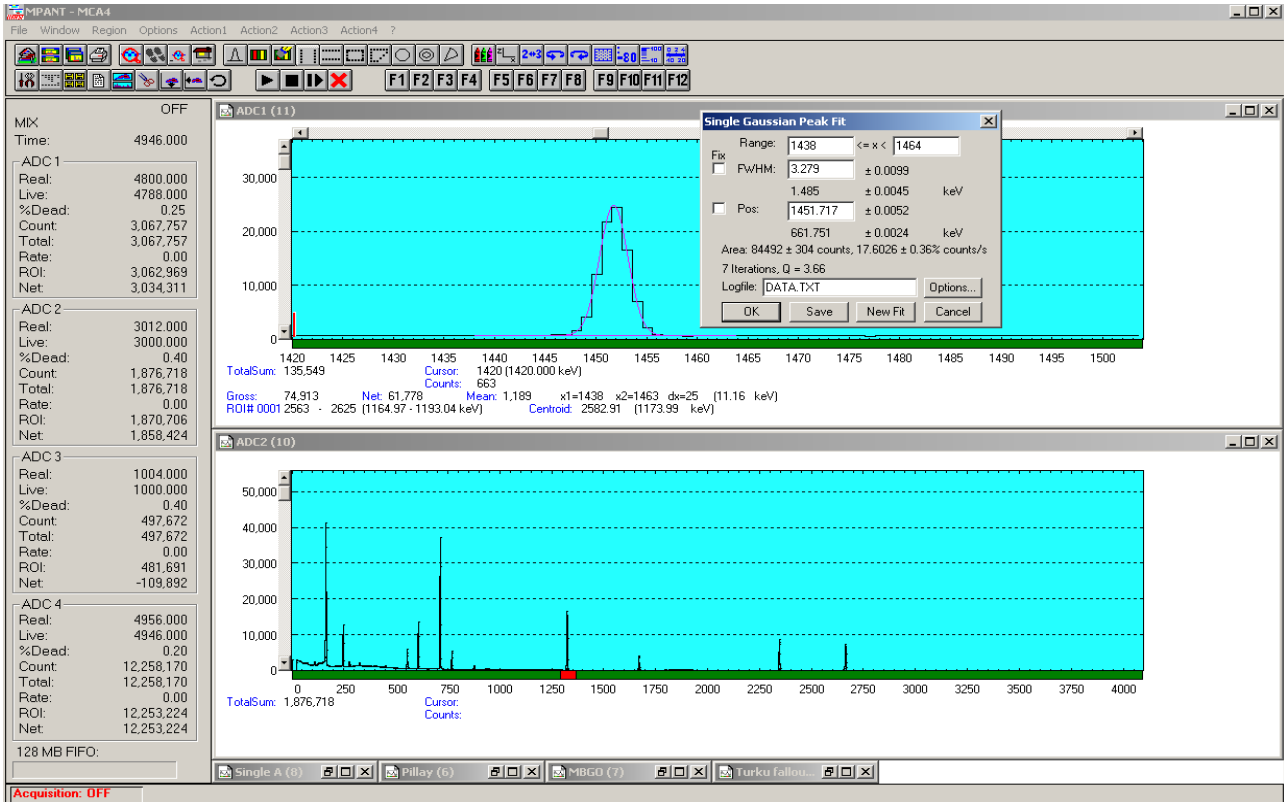


Fig. 5.9: Single Gaussian Peak Fit

The full width at half maximum FWHM and Position of the Gaussian can be changed and a New Fit can be performed, they even can be fixed to the entered value by marking the respective check box. The Position and FWHM are displayed in channels and also in calibrated units, if a calibration is available. The area of the Gaussian is also shown. For all values also the standard deviations are given. The value of Q is the normalized χ^2 . To take into account the systematic error of the line shape, you may multiply the errors with the square root of Q. Click on Save to append a line containing the results to a Log file with the specified name. OK closes the dialog and lets the fitted function in the display also if it is refreshed, whereas after Cancel the curve no longer will be shown in a refreshed display. Options... opens a new dialog box to define the information in the log file.

Peak Evaluation...



By selecting the Eval... menu item or the respective icon, a general peak evaluation with linear background is performed for the currently marked region. The calculated background curve is displayed and a dialog box shows the results:

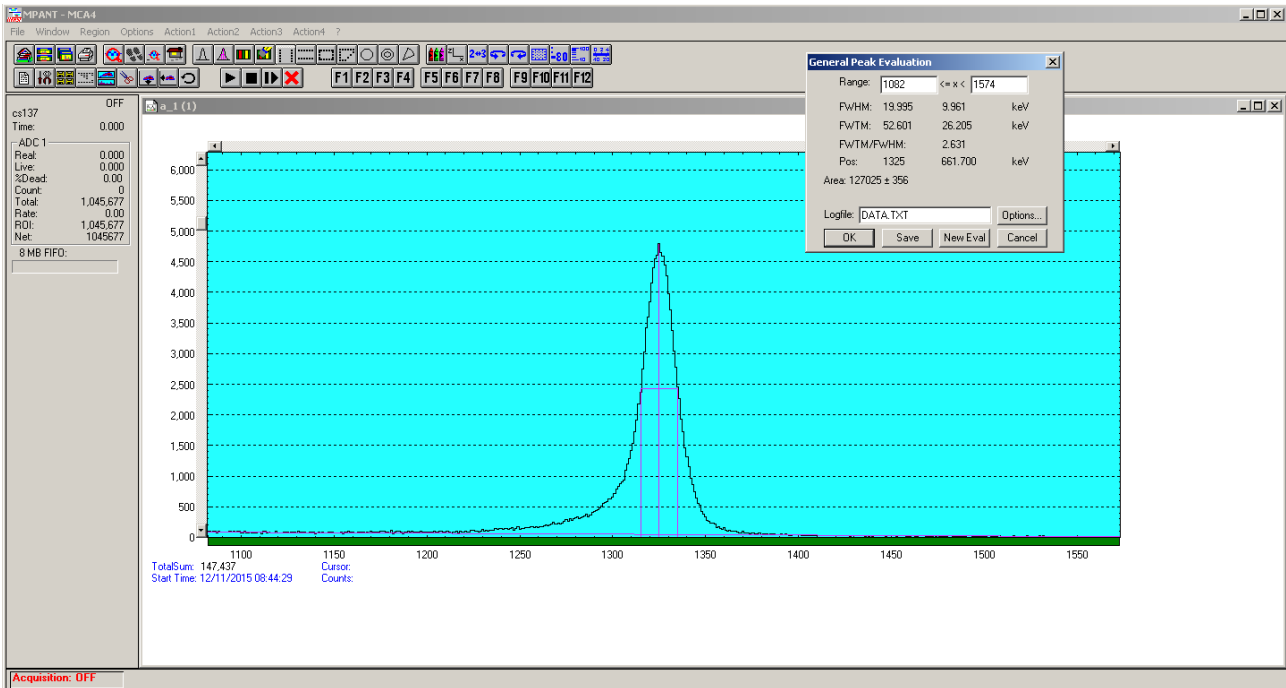


Fig. 5.10: Evaluation of a non-Gaussian peak

The Position and FWHM are displayed in channels and also in calibrated units, if a calibration is available. The area of the peak is shown, as well as the FWTM and FWTM / FWHM ratio. Click on Save to append a line containing the results to a Log file with the specified name. OK closes the dialog. Options... opens a new dialog box to define the information in the log file:

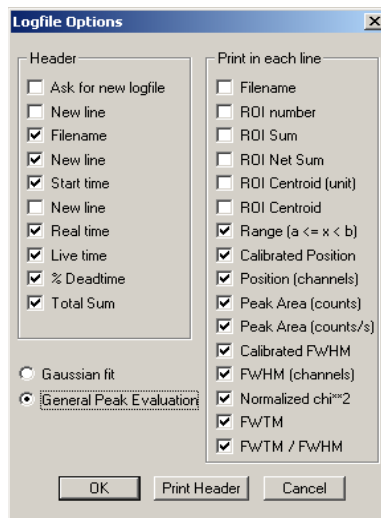



Fig. 5.11: Log file Options

Two radio buttons enable to select between Gaussian fit and general peak evaluation. The several quantities are written in standard text format with Tabs as separators and a Newline character at the end of each line, so the file can be read with standard calculation programs like EXCEL as a text file. Click on Print Header to write a header line.

Fit ROIs

With the Fit ROIs item, for all ROI's a Single Gaussian Peak Fit is performed and the results are dumped into the log file.

Auto Calib

Makes a Gauss fit for all ROI's in the active Display for which a calibration value was entered in the ROI editing dialog, and performs a calibration using the fit results. You have then only to check "use calibration" in the calibration dialog opened with  to use it.

5.4 Options Menu

The Options Menu contains commands for changing display properties like scale, colors etc., hardware settings, calibration and comments.

Colors...



The Colors menu item or respective icon opens the Colors dialog box.

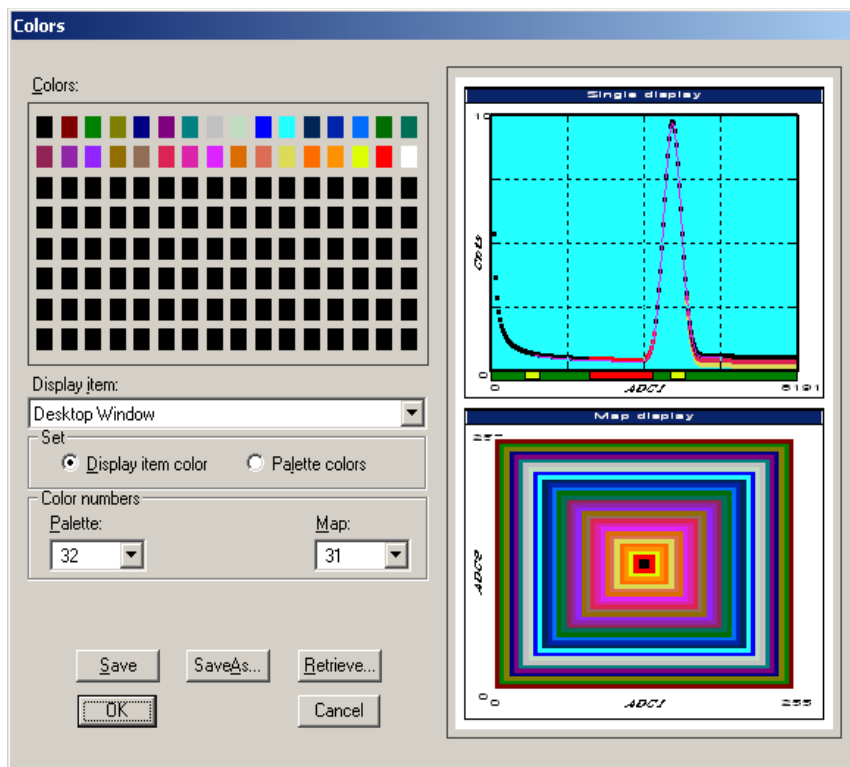


Fig. 5.12: Colors dialog

It changes the palette or Display element color depending on which mode is chosen. The current color and palette setup may be saved or a new one can be loaded.

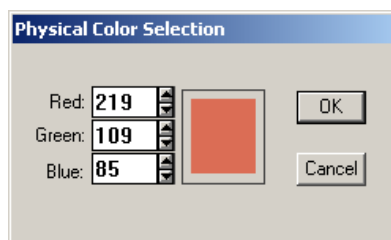


Fig. 5.13: Color Palette dialog

To change on of the colors, select "Palette colors" and click on one of the colors. In the Color Palette dialog box the RGB values can be edited or for a 256 color video driver one of the Physical palette colors can be chosen.

Display...



The Display menu item or the corresponding icon opens for single spectra the Single view dialog box.

Here the graphic display mode of single spectra can be chosen. The '**Type**' combo box gives a choice between dot, histogram, spline and line. The '**Symbol**' combo box gives a choice between None, Circle, Triangle down, Triangle up, Cross, Snow-flake and Diamond. The symbols can be filled by checking **Fill**, error bars can be displayed by checking **Error Bar**.

'Dot' means that each spectra point is shown as a small rectangle or the specified symbol, the size can be adjusted with the **size** combo box. 'Histogram' is the usual display with horizontal and vertical lines, 'spline l' means linear interpolation between the points, and 'line' means vertical lines from the ground to each spectra point.

If the displayed spectra range contains more channels as pixel columns are available in the video graphic display, usually only the maximum value of the channels falling into that pixel columns is displayed. But it can also explicitly specified by marking the check boxes "**Max Pixel**", "**Mean Pixel**", "**Min Pixel**" or "**All**" which value will be displayed. It is also possible to display all three possible Pixel values in different colors that can be chosen in the colors dialog. For the "Mean Pixel" a Threshold value can be entered; channel contents below this value then aren't taken into account for the mean value calculation. As plotting all values costs a lot of computing time for very large spectra, it is replaces automatically by "Max" and "Min" pixel for spectra with more than 1 million counts.

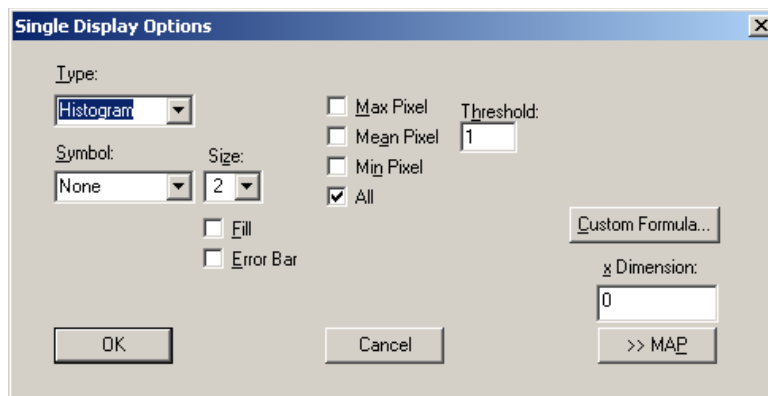


Fig. 5.14: Single Display Options dialog

The button labeled **Custom Formula...** provides spectra that are calculated in a special DLL named FMPA3.DLL. This DLL allows to calculate spectra and to display it with the MPANT program. The source code of this DLL is part of the optional DLL software interface, so it is possible to program customer defined functions. Select a single spectrum and open the Display options dialog. Press then the button labeled "Custom formula...". A dialog box "Custom-transformed spectra" is opened. Here the formula for the calculation can be selected and the parameters can be edited. A new display window showing the calculated spectra can be created. Error Bars can be activated in the display options dialog, for the calculation of the error bars also the DLL is used. Formula #6 calculates a Fourier transform of a wave form acquired in scope mode, see Fig. 2.14.

It is possible to change to a two dimensional view of the spectrum by specifying the x Dimension and clicking the button ">> MAP" from the Single Display Options dialog.

For MAP displays the Display Options dialog is changed and enables a choice between four Graphic types: **bitmap dot**, **vector dot**, **bitmap contour** and **vector contour**. Bitmap Dot is recommended as a standard, because it makes a good and fast display. The check box **Invert Roi Pixel** changes the color inside ROI's to make them better visible. Vector Contour paints colored contour lines. To calculate the lines takes a lot of time and causes the mouse pointer changing to an hourglass. But it gives very impressive colored pictures suited especially for presentation and when looking carefully at spectra details.

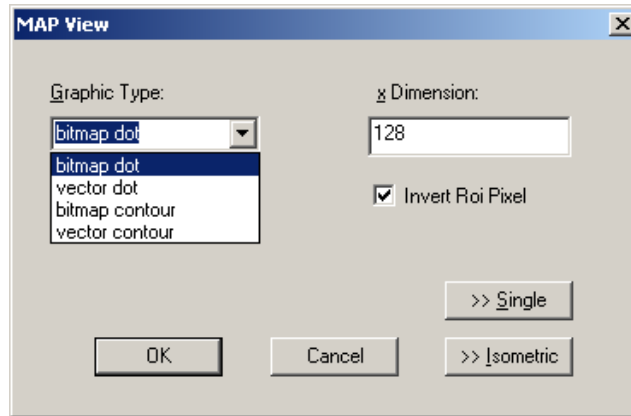



Fig. 5.15: Map Display Options dialog

From the MAP View dialog it is possible to change to Single view by clicking ">> Single" or change to Isometric View by clicking ">> Isometric". This can also be done directly via the respective toolbar icon .

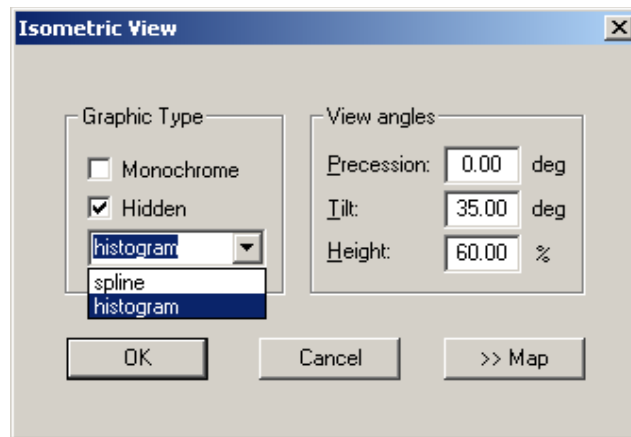




Fig. 5.16: Isometric Display Options dialog

In isometric mode several single spectra are drawn behind each other. The **Precession** angle around the vertical axis can be chosen in multiples of 90 degrees. This rotation can also be done directly using toolbar icons  and . The **Tilt** angle is between the x and y-axis and can be chosen between 15 and 89 degrees. The **Height** specifies the percentage of the z-axis length respective to the whole drawing, it can be entered between 0 and 99. With **hidden** it can be specified whether the hidden parts are not drawn. If "Monochrome" is checked, the spectra are painted monochrome, otherwise in color.

Axis...



By the Axis... menu item or the respective icon, the Axis Parameters dialog box is opened.

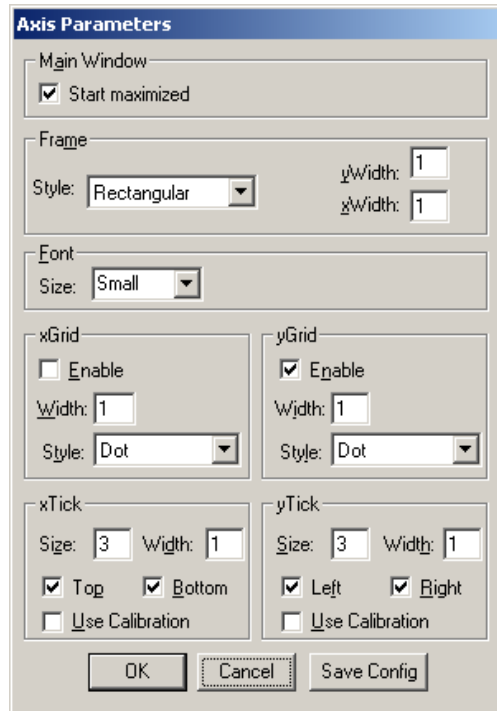


Fig. 5.17: Axis Parameters dialog

It provides many choices for the axis of a display. The frame can be rectangular or L-shape, the frame thickness can be adjusted (xWidth, yWidth). The font size can be chosen between Small and Large. A grid for x and y can be enabled, the style can be chosen between Solid, Dash, DashDot and DashDotDot. Ticks on each of the four frame borders can be enabled, the tick length and thickness can be chosen. The style of the axis labeling depends on enabled ticks at the bottom respective left side: If no ticks are enabled there, only the lowest and highest values are displayed at the axis, otherwise the ticks are labeled. 'Use calibration' changes the axis ticks and labels into calibrated values. A check box "Start maximized" allows to define if the MPANT main window should start in maximized size. It is necessary to save the MPANT options to make these settings permanent.

Scaling...



The Scaling menu item or the corresponding icon opens the Scale Parameters dialog box.

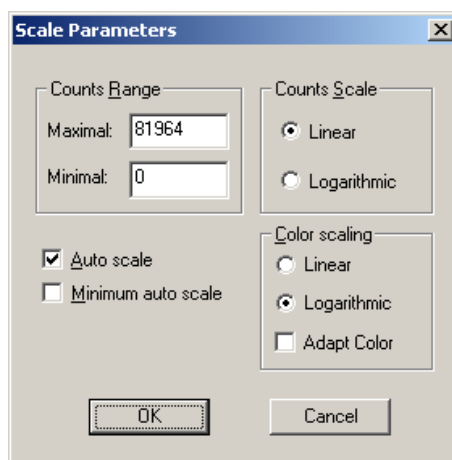


Fig. 5.18: Scale Parameters dialog

It enables changing the ranges and attributes of a Spectrum axis. By setting the Auto scaling mode, the MPANT will automatically recalculate the y-axis maximum value for the visible Spectrum region only. To keep

the same height of the visible region for a longer time, deselect the Auto scaling mode. Then with the scroll bar thumb one can quickly change the visible region scale, otherwise the scale will be changed automatically. The Minimum auto scale mode helps to display weak structures on a large background.

Lin / Log scale



For a Lin scale all data intervals have the same size. With Log scale the intervals will be small for small y values and large for large y values. All options have effect only on the active Display.

Calibration...



Using the Calibration menu item or the corresponding icon opens the Calibration dialog box.

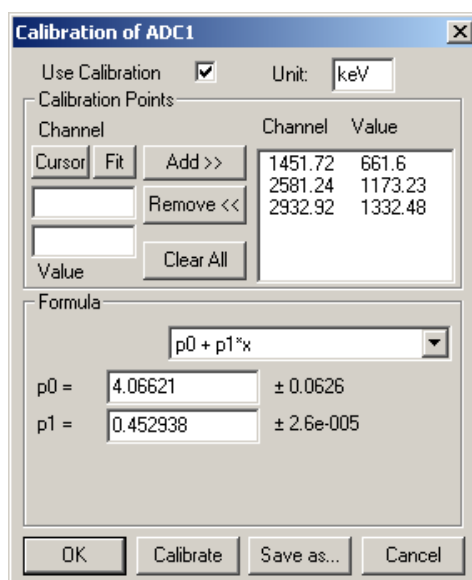


Fig. 5.19: Calibration dialog

Choose between several calibration formulas. Enter some cursor positions and the corresponding values. The actual cursor position can be entered by pressing 'Cursor' or the last fitted peak position by pressing 'Fit'. Click on Add to insert the calibration position into the list, then on Calibrate. The obtained coefficients can be inspected together with the statistical error, or they can be changed and entered by hand. If 'use calibration' is on, the calibrated values are displayed together with the channel position of the cursor. For MAP spectra the dialog is slightly different and allows to select by two radio buttons between the calibration of x- and y-axis. For a selected parameter there is only one calibration: If any ADC is already calibrated in a single spectrum and it is used also in a MAP spectrum, the calibration is the same here.

The fastest way to get a good calibration is to create ROI's around peaks with known energy and enter the energy in the ROI editing dialog, see Fig. 5.7. Then use the "Auto Calib" function in the region menu. You will immediately get a calibration and the dialog opens as shown in Fig. 5.19.

Comments...



Up to 13 comment lines with each 60 characters can be entered using the Comments dialog box. The content of these lines is saved in the data header file. The first line automatically contains the time and date when a measurement was started. The second line contains the title of the spectra. Note that editing this line changes the spectra name. This is the only method to change the name of physical ADC spectra. The description of each comment line can be changed by editing the file COMMENT.TXT.

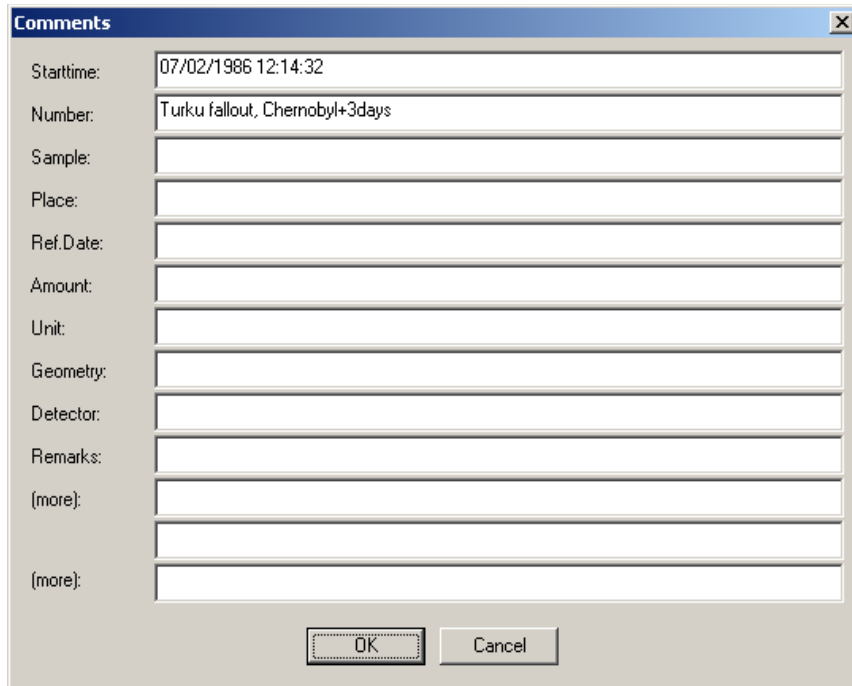


Fig. 5.20: Comments dialog

Hardware...



The Hardware settings dialog box allows to make all the respective settings (ref. 4.1.4).

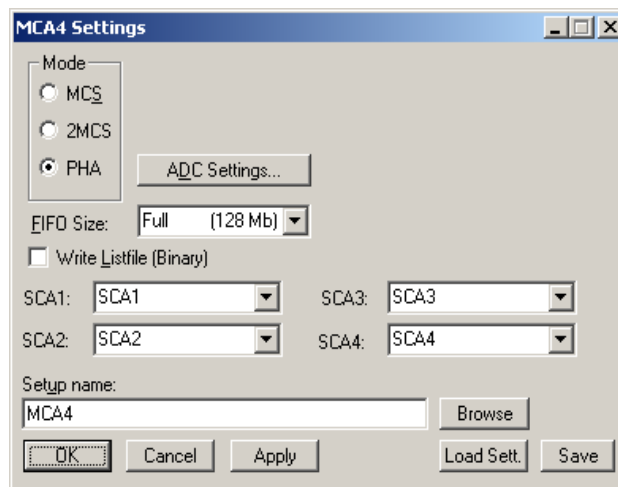


Fig. 5.21: Settings dialog

Data...



The Data Operations dialog allows to edit format settings and to load and save spectra data, ref. 4.1.3 .

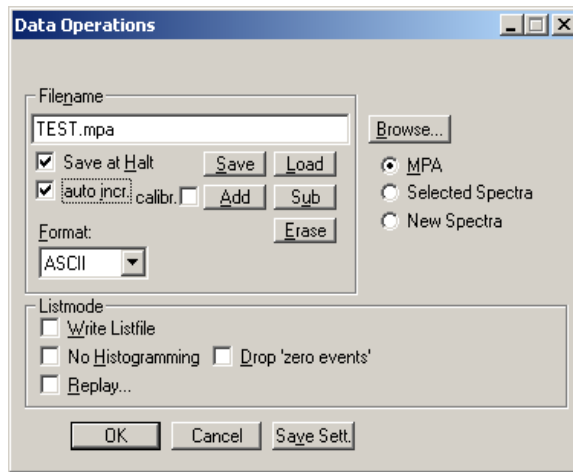


Fig. 5.22: Data Operations dialog

System...



The System Definition dialog box allows to make all the respective settings, see 4.1.10 .

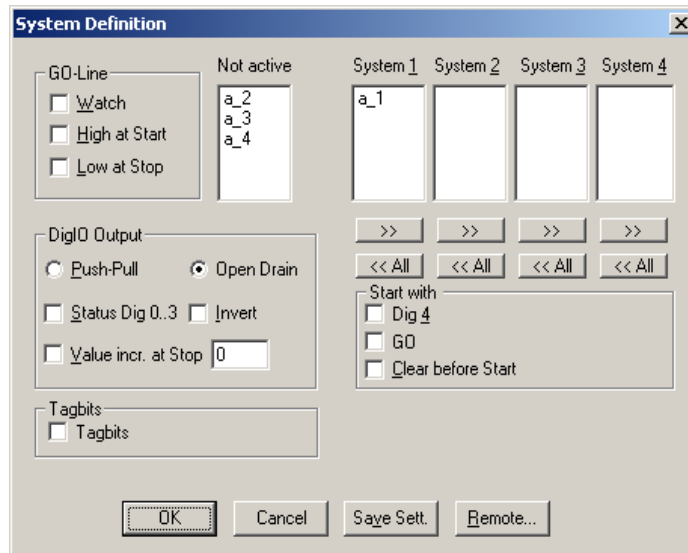


Fig. 5.23: System Definition dialog

Spectra...



The Spectra dialog box allows editing the list of calculated and dual parameter spectra, see chapter 4.1.11

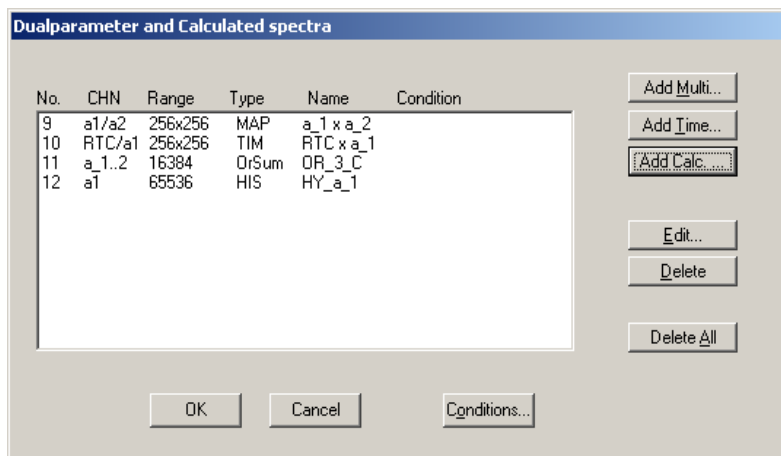


Fig. 5.24: Spectra dialog

Slice...



The Slice option allows to create new single spectra displays, showing a slice in a dual parameter spectrum. Click on a dual parameter spectrum display to make it the active window, then select the Slice menu item or toolbar icon. The Slice dialog box is displayed.

Select "x=const" or "y=const" for the slice direction, and the coordinate. Clicking the "create" button creates the new display window. In the title bar of the new window the name of the dual parameter spectra and the slice coordinate is shown.

The slice position can be changed using the scroll bar in the Slice dialog, or by entering the value in the edit field and pressing the button, which is labeled "Set" after creation of the slice view.

Clicking its close field can close the Slice dialog. Created slice spectra displays remain visible and their coordinates can be changed later using the Slice utility again. The position of the Slice dialog with respect to the MPANT main window can be saved in the MPANT.CFG file. Rectangular ROI's are visible in the slice spectra display and can be created here.



Fig. 5.25: Slice dialog

Replay...



The menu item Options – Replay... opens the Replay settings dialog (ref. 4.1.3).


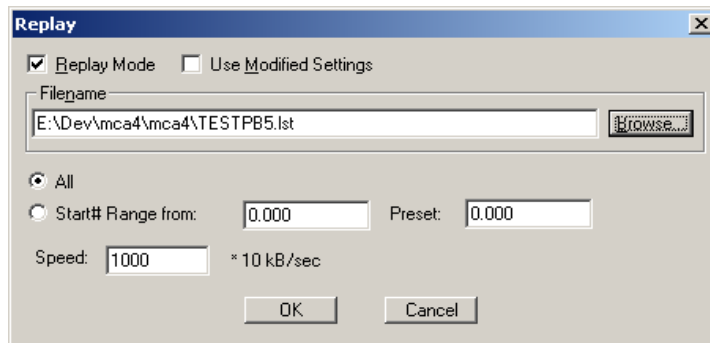


Fig. 5.26: Replay dialog

Tool Bar...

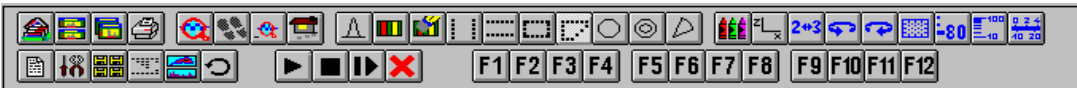

Selecting the Tool Bar Menu item opens the Tool Bar Dialog Box. It enables arranging the icons in the Tool Bar.

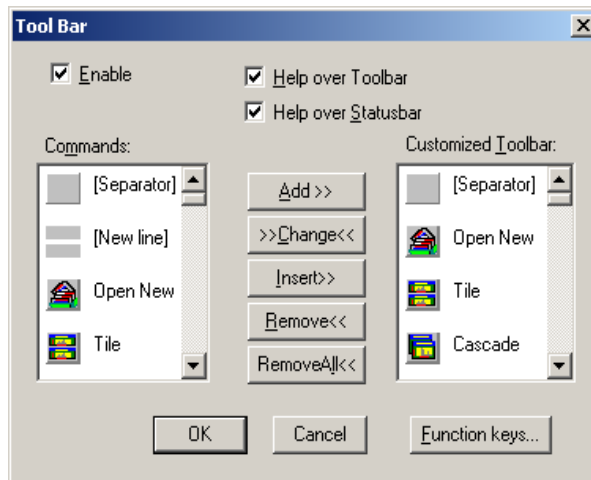


Fig. 5.27: Tool Bar dialog

If it is enabled, an array of icons in the MPANT Menu is shown. Clicking the left mouse button with the cursor positioned on an icon, the user can perform a corresponding MPANT Menu command very quick.

It is also possible to include icons for free programmable function keys F1...F12 into the Toolbar. The function keys can be programmed in the Function keys dialog with any one of the commands described in 4.3 . It can be accessed either by clicking the "Function keys..." button or directly from the Options menu.

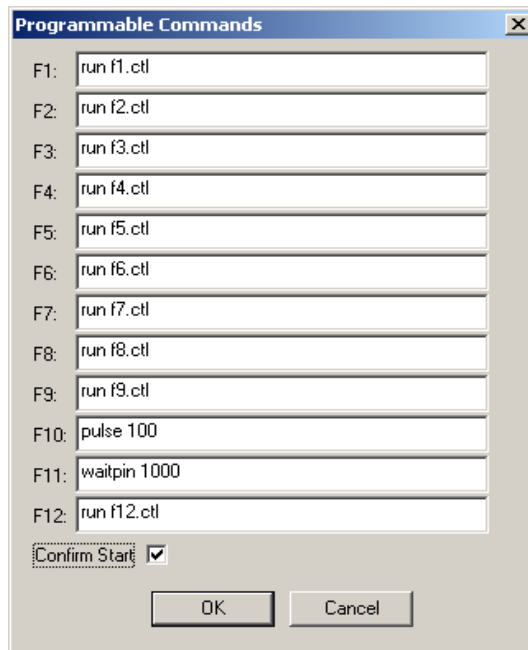


Fig. 5.28: Function keys dialog

The functions can be executed by clicking the corresponding icon in the toolbar or by the corresponding function key on the keyboard simultaneously with the CTRL key. The MPANT window must be the active on the desktop and have the focus, i.e. it must be the top window. Uncheck the “Confirm Start” check box if you don't like the message box asking for a confirmation before starting a new acquisition.

Status bar

With this menu item the Status bar at the bottom of the MPANT main window can be switched on or off. A corresponding check mark shows if it is active or not. The Status bar usually shows if an acquisition is active. When the left mouse button is pressed while the mouse cursor is within a toolbar icon, it displays a short help message what the meaning of the toolbar icon is.

Status window

The status window at the left side of the MPANT main window can be shown in a small or large size or can be hidden.

Save As...

Saves all parameters defined in the Options menu to the MPANT.CFG or a user defined config file.

Open...

Loads a new configuration.

5.5 Action Menus

Start



The Start toolbar button erases the data and starts a new measurement.

Halt



The Halt toolbar button stops a measurement.

Continue



The Continue toolbar button continues a measurement.

Erase



The Erase toolbar button erases the data.

If more than one systems are defined for starting and stopping combined MCAs together, more groups of start / halt / continue / erase icons appear in the toolbar for the respective systems.

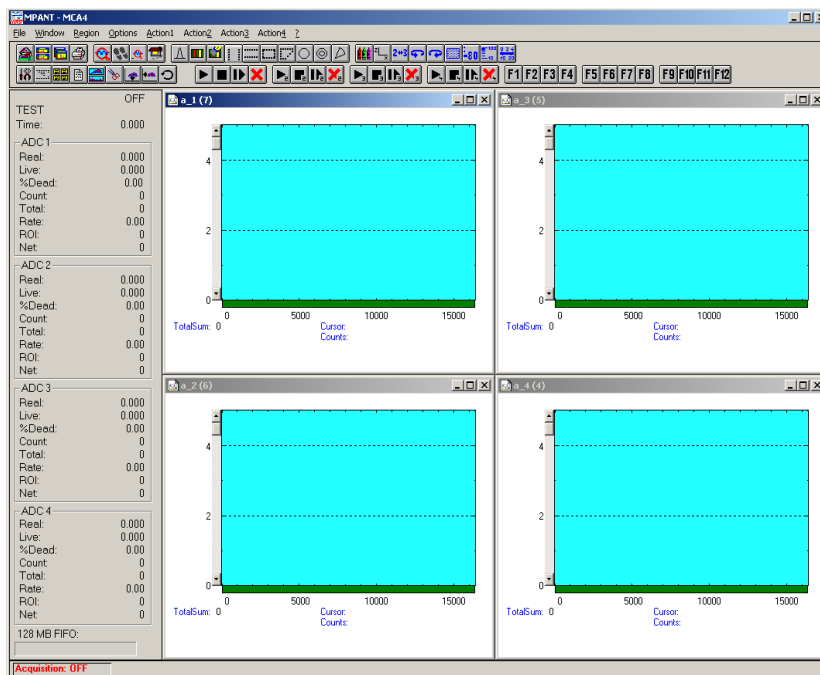


Fig. 5.29: MPANT with four systems enabled

6 Programming and Software Options

The MCA4A can be controlled by user-written programs using the DLL software interface with example programs for Visual Basic, LabVIEW and C that is available as an option. Furthermore, LINUX software is available as an option containing a driver, library and console test program. A Windows software similar to the LINUX package that runs without the server using a stand-alone DLL is also available on demand for customers who own one of the two available library packages.

7 Appendix

7.1 Absolute Maximum Ratings

Input voltage:	any ADC input: -16 ... +17 V
	any GATE / MCS input (short duration only): -1.7 ... +5.5 ⁴ V
	any FEATURE I/O port: -0.5 ... +5.5 ⁵ V
	any Counter input: -0.5 ... +7.0 V
	DC power supply: +14 V
DC input current:	any ADC input: ±17 mA
	any GATE / MCS input: ±70 mA
	any DIGITAL I/O port: ±10 mA
	any Counter input: ±50 mA
Input power:	any GATE / MCS input: max. 0.25 W
	any Counter input (termination resistors): max. 0.5 W

7.2 Recommended Operating Conditions

Power supply: nominal +12 V _{DC} / 15 W
	voltage range: +11 ... +14V _{DC}
Operating temperature range: 0 to +50°C
GO-Line load: min. 1 kΩ to +5.0 V
	or min. 2 kΩ to GND
Digital I/O input voltage: 0 ... +3.3 V

7.3 Specifications

7.3.1 ADC inputs

Input impedance:	10V, standard 1.33 kΩ
	(8V, option) 1.07 kΩ
	(2.5V, option) 1.31 kΩ
	(2V, option) 1.05 kΩ
Input voltage range 10V standard:	unipolar nominal 0 ... +10 V
	bipolar nominal ±5 V
Other ranges settable by jumpers:	
8V range:	unipolar nominal 0 ... +8 V
	bipolar nominal ±4 V
2.5V range:	unipolar nominal 0 ... +2.5 V
	bipolar nominal ±1.25 V
2.0V range:	unipolar nominal 0 ... +2 V
	bipolar nominal ±1 V

7.3.2 GATE / MCS inputs

Input impedance (selectable):	ECL 50 Ω
	TTL (+ 3.3V) 50 Ω / 4.7 kΩ

⁴ respect maximum input power limit

⁵ Maximum DC undershoot below GND must be limited to either -0.5V or 10mA. During transitions, undershoot to -2.0V or overshoot to +7.0V is allowed, provided this lasts less than 10ns and with the current is limited to 200mA

Input voltage range:	positive (TTL, slow NIM) 0 ... +3.3V (5 V compatible) negative (fast NIM) 0 ... -800mV
Input threshold voltage:	TTL settable in steps of 1 mV, rising or falling edge..... -2V..+2V ECL adjustable in steps of 1 mV, rising or falling edge..... -2V..+2V
Pulse width (high, low):	(MCS START and CHADV) min. 1.0 ns (MCS STOP 1 and 2) min. 0.8 ns
Frequency:	(MCS START and CHADV) max. 100 MHz (MCS STOP 1 and 2) max. 400 MHz
Slew rate requirement: min. 25 V/ μ s

7.3.3 SCA outputs

Output voltage (at IC pin):	V_{OH} ($I_{OH} = -4.0mA$) min. 2.4 V V_{OL} ($I_{OL} = 8.0mA$) max. 0.4 V
Output impedance:	(series resistor) 22 Ω
SCA output pulse width: 120 ns

7.3.4 Digital I/O ports

R Pull:	(default pull-up to +3.3V) 1.0 k Ω
R I/O:	(default series resistor) 22 Ω
Input voltage (at IC pin):	V_{IH} min. 2.0 V V_{IL} max. 0.8 V
Output voltage (at IC pin):	V_{OH} ($I_{OH} = -4.0mA$) min. 2.4 V V_{OL} ($I_{OL} = 8.0mA$) max. 0.4 V

7.3.5 Reference clock I/O

Input frequency:	(auto dtected) 5 ... 100 MHz
Output frequency: 10 MHz
Impedance: (dc blocked) 50 Ω
Output HIGH voltage:	(50 Ω load) typ. 2.0 V
Output LOW voltage:	(50 Ω load) typ. 0.4 V
Input amplitude:	(clipped sine or sine wave) typ. 3.3 V_{PP}

7.3.6 Octal Counter inputs

Count frequency: max. 125 MHz
Input impedance: 50 Ω / 0.5 W
Input pulse width:	(high or low) min. 4 ns
Input HIGH voltage: min. 2 V
Input LOW voltage: max. 0.4 V

7.3.7 RTC RESET input

Input impedance: 50 Ω / 0.5 W
Input pulse width:	(high or low) min. 8 ns
Input HIGH voltage: min. 2 V

Input LOW voltage: max. 0.4 V

7.4 Performance

7.4.1 General

DDR2 Memory: 1Gbit = 32M x 32 bit
 Memory modes: (FIFO) Listmode
 Data throughput: (USB transfer rate) MB/s
 ADC data buffer memory: (each input) 32k x 16 bit
 segmentable in 1 / 2 / 4 / 8 equal sized blocks
 Basic operating modes⁶: Pulse Height / Multichannel Analysis PHA
 Oscilloscope SCOPE
 Frequency Domain FFT
 Sampled Voltage Analysis SVA
 Multichannel scaling MCS
 RTC (Real Time Clock) resolution: 1 ms
 (optional) 8 ns
 Reference oscillator: Stability: ±50 ppm
 Aging: ±3 ppm / year

7.4.2 Analog-to-Digital Converter

Conversion time: abs. Maximum PHA 8 ns
 Gaussian fit PHA $(2 \times i1^7 + 96) \times 8$ ns
 Flattop fit $(M^8 + 42) \times 8$ ns
 Live- / realtime resolution: (standard) 1 ms
 (RTC option) 8 ns
 Timer preset range: (standard) $< 2^{48} \times 1$ ms
 (RTC option) $< 2^{48} \times 8$ ns
 Reference voltage: temperature drift: max. 3 ppm/°C
 Resolution: Gaussian fit (32k range) typ. 1.7 ... 2.8 FWHM
 Integral non-linearity⁹: (for 97% of the channels, ref. Fig. 7.1 and Fig. 7.2) $< \pm 3.5$ %
 Differential non-linearity¹⁰: (Gaussian fit and 97% of the channels, ref. Fig. 7.3) $< \pm 1$ %
 GATE to threshold crossing: t_{setup} TBD
 t_{hold} TBD
 TAG bits sampling period: 2.5 μs

⁶ some modes may be used simultaneously

⁷ fit is performed through $(2 \times i1 + 1)$ sample data ($1 < i1 < 64$)

⁸ M = number of samples averaged ($0 < M < 256$)

⁹ $INL(n) = [\text{COUNTS}(n) / \text{Mean}(\text{Range})] - 1$, with $n = 300 \dots 31300$, Range = 31001
 includes non-linearities of the ramp generator, tail pulser and shaping amplifier

¹⁰ $DNL(n) = \text{COUNTS}(n) / [0.5 \times \text{COUNTS}(n-1) + 0.5 \times \text{COUNTS}(n+1)]$, with $n = 300 \dots 31300$

Appendix

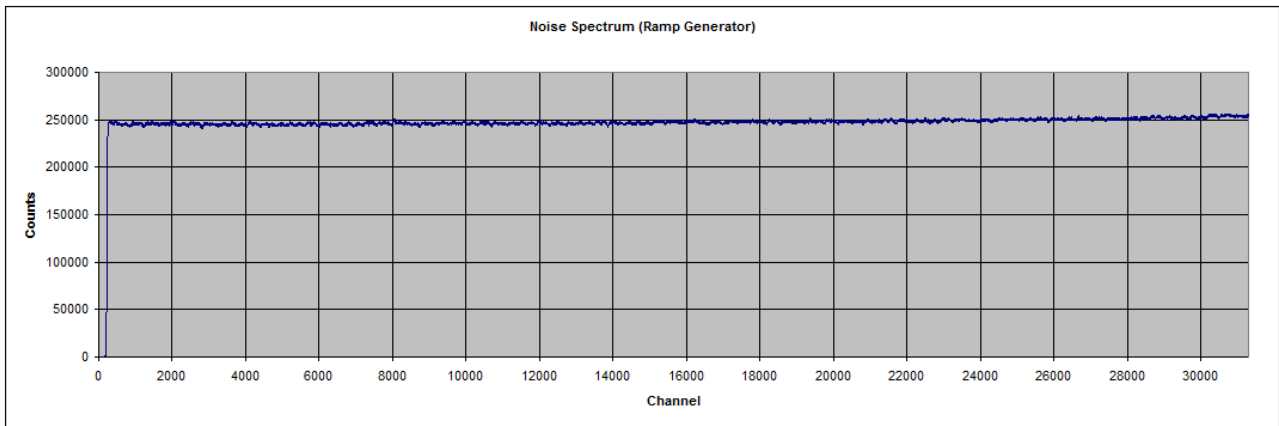


Fig. 7.1: ADC noise spectrum (generated with a ramp generator and Gaussian fit)

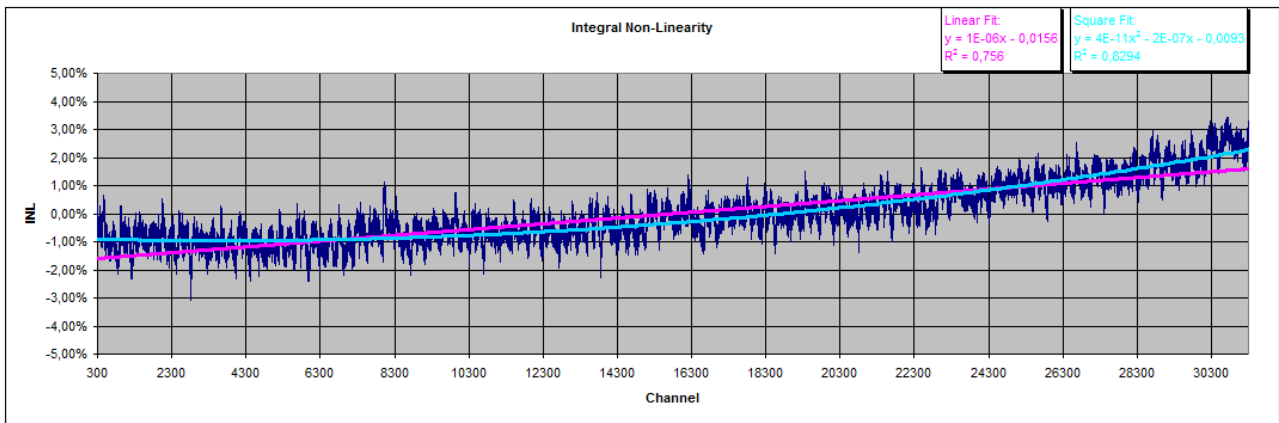


Fig. 7.2: ADC integral linearity plot for 32k ADC range

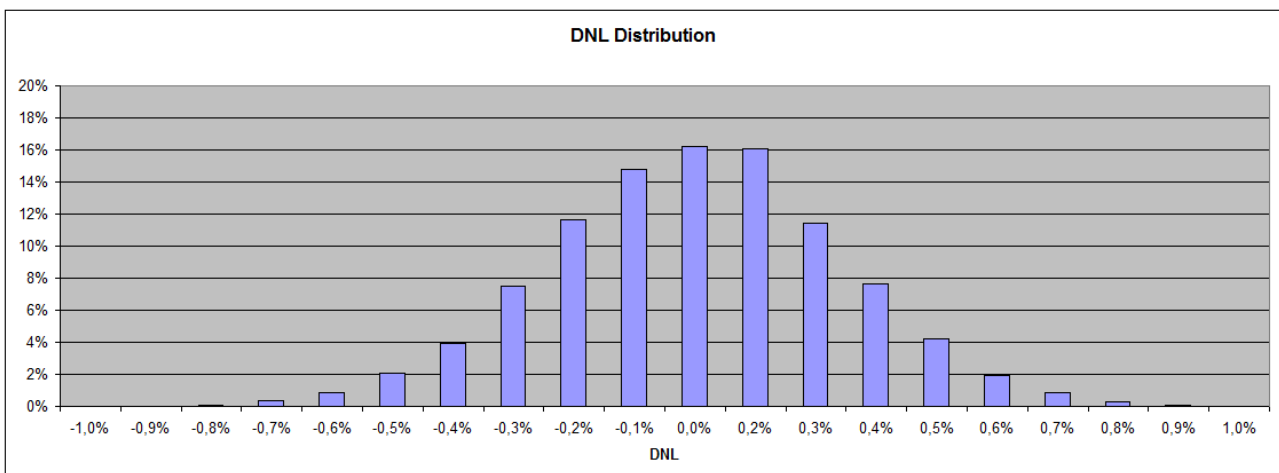
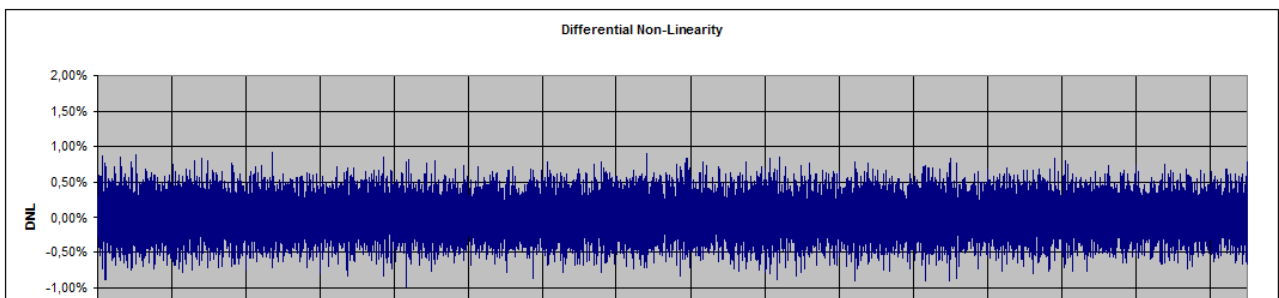


Fig. 7.4: ADC DNL distribution

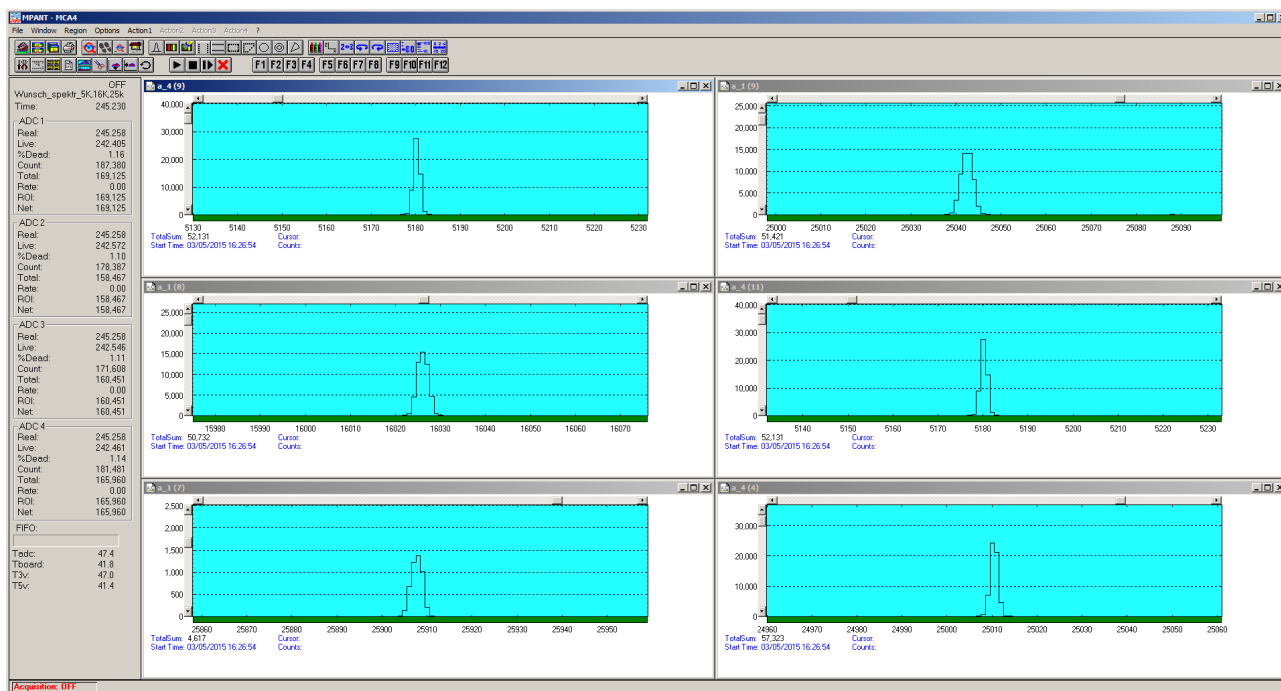


Fig. 7.5: Typical peak widths of Gaussian fitted pulses (32k range, 4 μ s shaping)

7.4.3 Counters

Connector:	15 pin high-density female D-SUB
Count frequency:	max. 125 MHz
Location of Carry outputs:	ADC GATE connector
Input Impedance:	50 Ω / 0.25 W
Input pulse width:	(high or low).....	min. 5 ns
Input HIGH voltage:	min. 2.0 V
Input LOW voltage:	max. 0.8 V
Output impedance:	(series resistor, Carry outputs).....	33 Ω / 0.1 W
Output HIGH voltage:	$I_{OH} = -0.1$ mA.....	min. 3.1 V
Output LOW voltage:	$I_{OL} = 0.1$ mA.....	max. 0.2 V

7.4.4 Multichannel scaler

Dwell time modes:	internal (auto), external
Internal dwell time:	(1 channel MCS)	min. 30 ns
	(2 channel MCS)	min. 60 ns
	max. $(2^{48} - 1) \times 10$ ns = 32 days
	in steps of:	10ns
External channel advance:	edge sensitivity:	rising / falling
	period (internally asured):	min. 100 ns
	uncertainty:	± 5.0 ns
Trigger uncertainty:	± 5.0 ns
Trigger (START) delay:	45.0 ns ± 5.0 ns
Count (STOP) delay:	41.25 ns ± 1.25 ns

Appendix

Dead time:	End-of-sweep:	10 ns
	between time bins:	none
Sweep repetition time:	\geq (Range + 10ns)
Differential non-linearity:	(channels 2 ... 32000, 50ns dwell time, ref. Fig. 7.7)	$< \pm 0.1 \%$
Sweep counter:	48 bit
Sweep counter preset range:	1 ... 2^{48}
Spectrum sequencing:	max. 65536 spectra
End-of-sequence dead time:	(auto sequence mode)	30 ns

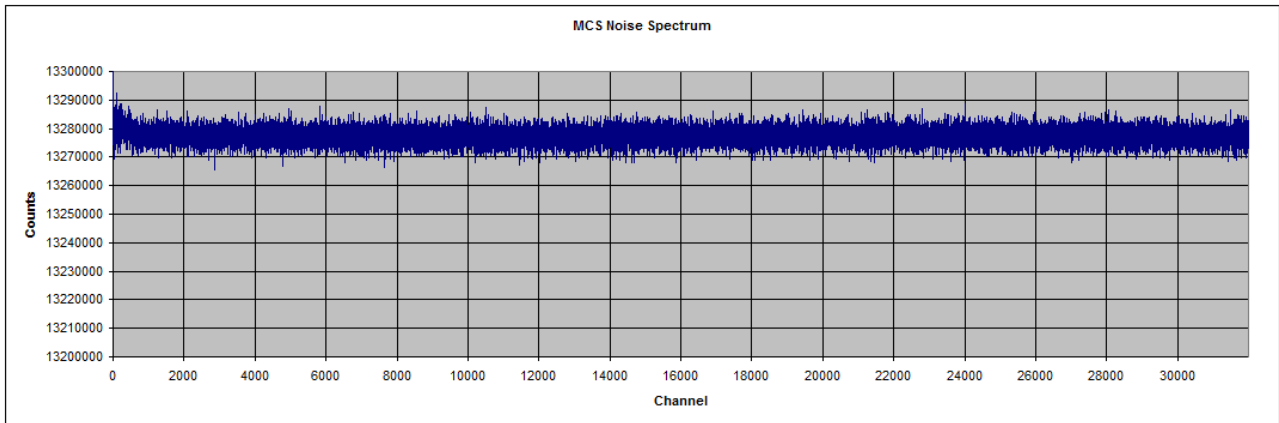


Fig. 7.6: MCS noise spectrum (32k range, 50ns dwelltime)

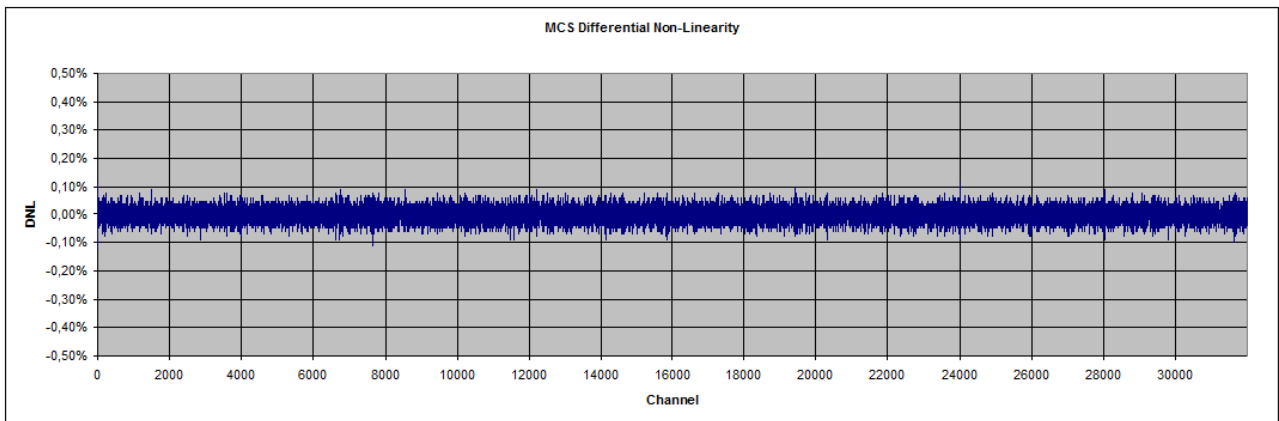


Fig. 7.7: MCS differential non-linearity plot (32k range, 50ns dwelltime)

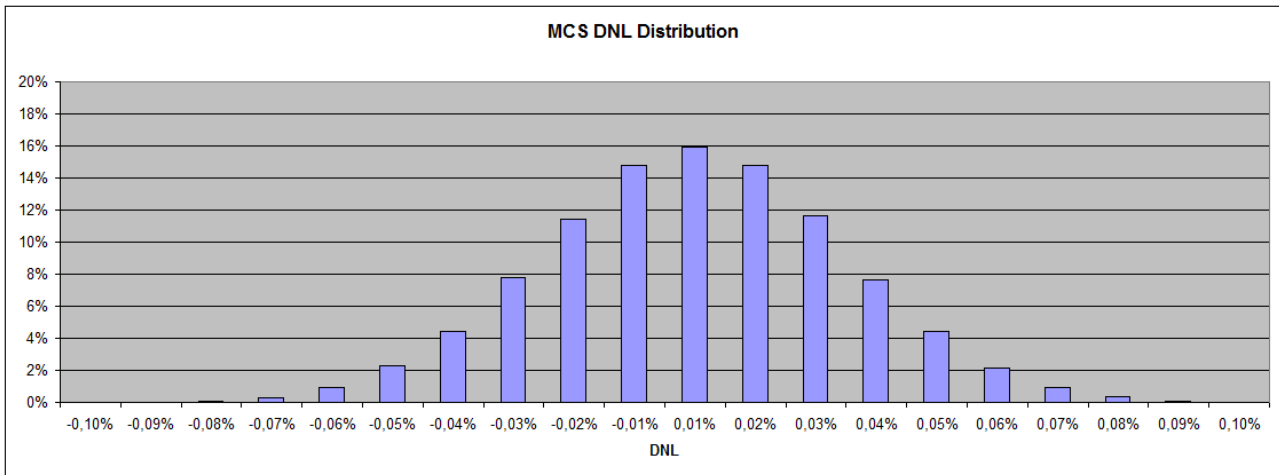


Fig. 7.8: MCS DNL distribution

7.5 Physical

Case material:	aluminum
Size:	260 mm x 48 mm x 275 mm
Weight:	1.7 kg
Shipping case:	Paper box with two foam pack, typ. 420 mm x 320 mm x 160 mm, 3.4kg	

7.6 Accessories

Included:

- USB 3.0 (A/A) cable, 3m , USB 2.0 compatible,
- Power supply 100-240 V_{AC} / 12 V_{DC} , min 2.0 A
- Power supply cord (if necessary)
- D-sub15 HD plug with plastic cap

Optional:

- Oven-disciplined internal 10 MHz reference oscillator
- Counter cable BNC or Lemo

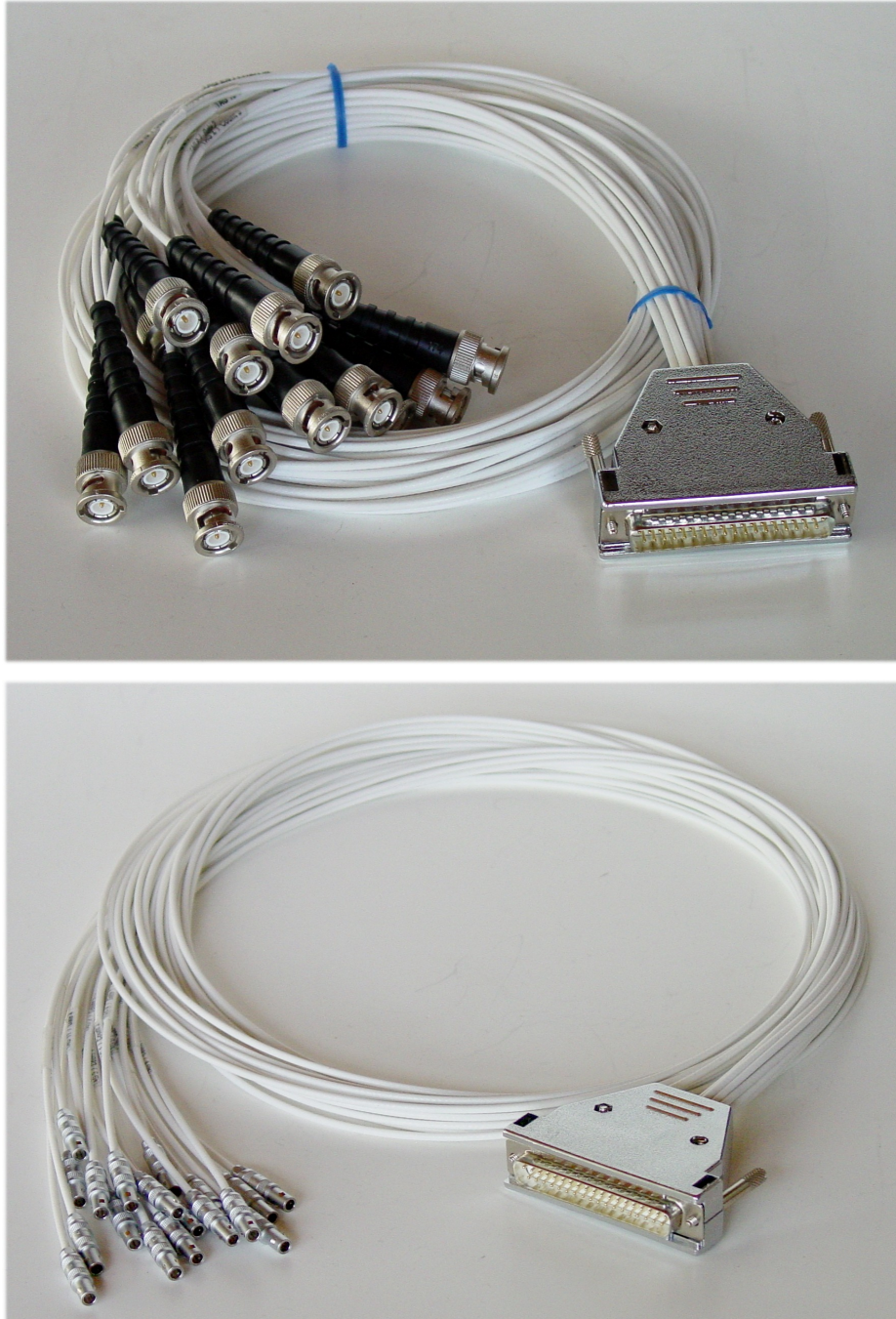


Fig. 7.9: Counter Cable BNC or Lemo

7.7 Trouble Shooting

- The 'POWER ON' LED does not light up:
 - 1) check that power is applied to the external power supply
 - 2) check the 12V supply cable is properly plugged in
 - 3) check that the external power supply works – check it's output voltage
- The MCA4A is not found by the host computer
 - 1) check that the USB cable is properly plugged in
 - 2) Remove the MCA4A power supply, wait 1s and reapply power

7.8 Personal Notes